



MAR 22, 2023

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.bp2l69b91lqe/v1

Protocol Citation: maggie.bowman 2023. Core Receiving and Splitting . **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bp2l69b91lqe/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Oct 06, 2022

Last Modified: Mar 22, 2023


PROTOCOL integer ID:
70972

Core Receiving and Splitting

In 2 collections

maggie.bowman¹

¹Environmental Molecular Sciences Laboratory

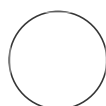

nicholas.sconzo


ABSTRACT

The method covers how soil cores were received and depth stratified.







MATERIALS











- Site Cores (Cores A, B, C1-4)
- Core Extruder (Long stick with flask disk on end by broom in 1446)
- Aluminum Foil
- 70% ethanol
- Kimwipes
- Gloves
- Bench Scraper
- Ruler/Measuring tape in cm
- Sharpies/writing utensils
- 4mm sieves and base
- Processing kit (Box with tubes, bags, tins, and labels)
 - 50 mL conical tubes (for both top and bottom)
 - 1 pH
 - 1 P - extraction
 - 6 N (3 for extraction and 3 for extractant)
 - 9 for microbial biomass (3 MicB, 3 MicB_Back, 3 MicB_Fum; only MicB gets filled with soil)
 - 6 unlabeled 50 ml tubes for P extract (will get labeled following pH measurements)
 - Whirlpak bags
 - DNA (Top and Bottom)
 - Removed rock and roots (Top, Middle, and Bottom)
 - Large aluminum pans (4 total for Top (2), Middle, Bottom; Extras are stored in 1106 if needed)
 - Small aluminum weigh boats
 - GWC – 3 each for top and bottom

- 1 When cores are received first check that there are cores labeled Core A, Core B, and Cores C1-4 included in the sample cooler.

- 2 Place Core A and C4 back into the site bag and store vertically in the  4 °C refrigerator in 1446.
- 3 Spray down the table with a 70% EtOH solution and wipe down using kimwipes.
- 3.1 Keep each core, if processing multiple simultaneously, in its own dedicated area/table.
- 4 Wearing gloves, tear two large sections of foil (~2.5 ft in length) and place them overlapping by half onto the cleaned tabletop.
- 5 Using a permanent marker write the site code at one end of the foil.
- 6 Lay the Core B core on the foil horizontally and remove the caps (top and bottom) noting which section of the core is top on the foil.
- 7 Clean the core extruding tool with 70% EtOH, align the flat end of the extruder with the bottom exposed soil.
- 8 While maintaining a horizontal position on the bench slide the tube up the core extruders handle to release the core onto the foil.

- 9 Take photographs of the soil core with and without a measuring tape making sure to include the site label in the photos.
- 10 Record the cores total length and how well it maintained the shape of the tube in the core processing files.
 - 10.1 [Processing_Template](#) *actual template
 - 10.2 Make sure to save the Processing_Template document with the site codes for example: "Processing_PRS1_PRS2.xlsx"
- 11 In the Processing file record the core receiving and splitting information.
- 12 Using a bench scraper and a measuring tape slice the core at 10 cm from the top and 10 cm from the bottom. Record the depth range for each section on the processing forms.
- 13 Using 70% EtOH wipe out clean (rinsed with water) 4 mm sieve. Stacking the sieve on top of the lower pan scrape each section (top, middle, and bottom) into its own clean and labeled sieve.
 - 13.1 In only the top 10 cm section add cores C1-3 during sieving to homogenize.
 - 13.2 Do not combine the top middle and bottom sections at any point.



- 14** Gently shake and stir (with a clean sterilized gloved hand) the soil in the sieve being careful not to force soil through sieve.
- 14.1** If the soils have a high clay content you may need to use a sterilized paint brush to break up the particles using a stippling motion (vertical tapping).
- 14.2** Break up large aggregates between fingers and continue sieving.
- 15** Once sieving is completed collect all items remaining on the sieve into a whirlpak bag labeled with the site code, removed, and depth fraction for example (1000S_PRS1_Removed_ TOP). Record the mass of the bag and the total mass (bag + removed rocks/roots) on the corresponding section of the Processing file.
- 15.1** Do not tare out the whirlpak prior to adding soil.
- 16** To measure gravimetric water content (GWC) on the top and bottom core sections record the mass of the empty tin and the mass of the tin +  10 g of soil. Place the tins in the  60 °C degree oven for  48:00:00 to  72:00:00 or a  100 °C degree oven for  24:00:00 . GWC is collected in triplicate and make sure the tins are labelled with the 1000S_SITE_GWC_DEPTH_1-3. 6d
- 16.1** Make sure to record the wet mass and dry mass (after  24:00:00 or  48:00:00) in the processing file. 3d

- 17** Collect a minimum of  10 g of soil for DNA into a whirlpak bag for the top and bottom sections. These samples will be labelled with the site code, DNA, and depth fraction (Top and Bottom). For example: 1000S_PRS1_DNA_Top
- 18** Collect  20 g (within  0.1 g) of soil for pH into a  50 mL falcon tube for the top and bottom sections. The samples should be labelled with the site code, pH, and depth fraction. For example: 1000S_PRS1_pH_Top.
- 19** Collect ~  12 g of soil for phosphorus extraction into a  50 mL falcon tube for the top and bottom sections. The samples should be labelled with the site code, P, and depth fraction. For example: PROS_P_Top.
- 19.1** Following pH samples with a pH less than 7 will be extracted using the Bray Method, while samples with pH greater than 7 will be extracted using the Olsen Method.
- 20** To measure microbial biomass and N extractions place exactly  8 g (within  0.02 g) of soil into six  50 mL falcon tubes for the top and bottom sections. Record the mass onto the MicrobialBiomass_N tab. The 6 tubes will be labelled with three for microbial biomass and three for N extraction. For example: PRS1_N_Top_1-3 & PRS1_MicB_TOP_1-3
- 21** The remaining sieved soil will be placed in a labelled clean aluminum pie pan and covered with foil and stored at  4 °C .
- 21.1** Soils are air dried in the BSC hood in 1521. Depending on the mineral composition soils can take longer to air dry.
- 21.2** Book time for the BSC hood using this calendar: [1521 BSC Hood Calendar](#)
- 22** Repeat the sieving and subsampling steps for the bottom section and sieve and store the middle



section.

22.1 Reminder Core A and Core C4 are not processed at this time.




23 Store samples in the appropriate location:

23.1 Store the DNA samples in a zip-top bag labelled 1000 soils and in the  -80 °C .
These samples are typically on the middle shelf of the furthest east  -80 °C in the corridor behind 1521.

23.2 Refrigerate the pH and P extraction samples until done processing cores.

23.3 Place GWC tins in the drying oven in 1446 at 60C for  48:00:00 to  72:00:00 .

5d

23.4 Dry the remaining soils in the BSC hood for  48:00:00 to  72:00:00 or until dry. After soils are dry collect a  50 mL Olympus centrifuge tube of OM analysis and seal the remainder in labeled Polypropylene jars for archive.

5d