



Sep 08, 2020

SARS-CoV-2 Milk ELISA

Powell RL¹¹Icahn School of Medicine at Mount Sinai**1** Works for me dx.doi.org/10.17504/protocols.io.bk2skye[XPRIZE Rapid Covid Testing](#) [Powell Lab/Lactiga Collab](#) rlp285

ABSTRACT

Currently, the COVID-19 diagnostics sector includes Point-of-Care and ELISA tests for IgG and IgM that can assess saliva or blood as described by the CDC [Centers for Disease Control and Prevention]. Similarly, the FDA webpage for EUA Authorized Serology Test Performance makes no reference to IgA COVID-19 assays [Food and Drug Administration]. Taken together, there are no direct competitors in the mucosal IgA quantitative diagnostics space.

To address this gap in diagnostic coverage we have generated a novel proof-of-concept dataset on the use of COVID-specific IgA antibodies as a marker of infection. As described below, our assay dataset demonstrates utility across several parameters:

- Highly specific
- Highly sensitive
- Works with undiluted human milk, simplifying sample collection and processing
- Defines a novel correlation between COVID-specific IgA and Secretory Component (SC), where SC is a known surrogate for dimeric IgA, the predominant IgA type in mucosal secretions
- Establishes that dimeric IgA antibodies detected by the assay are capable of neutralizing SARS-CoV-2

THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

<https://www.medrxiv.org/content/10.1101/2020.05.04.20089995v1>

DOI

dx.doi.org/10.17504/protocols.io.bk2skye

PROTOCOL CITATION

Powell RL 2020. SARS-CoV-2 Milk ELISA. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bk2skye>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

<https://www.medrxiv.org/content/10.1101/2020.05.04.20089995v1>

KEYWORDS

SARS-CoV-2, ELISA, milk, IgA, sIgA

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 08, 2020

LAST MODIFIED

Sep 08, 2020

PROTOCOL INTEGER ID

41778

GUIDELINES

keep milk at room temperature for no more than 5 hours, or else keep at 4C for up to 5 days or freeze.

MATERIALS

NAME	CATALOG #	VENDOR
Sulfuric acid	A300C-212	Fisher Scientific
Corning® 96-Well High-Binding Flat-Bottom Microplates 6 Plates/Pack	38019	Stemcell Technologies
Gibco™ (Phosphate Buffered Saline) Solution, pH 7.4 (PBS)	10010-049	Fisher Scientific
Lab Vision™ Tween™ 20 Detergent	cat# TA-125-TW	Thermo Fisher Scientific
BSA	30-AB79	Fitzgerald
Skim Milk Powder		
1-Step™ Ultra TMB-ELISA Substrate Solution	34028	Thermo Fisher
Normal Goat Serum	PCN5000	Thermo Fisher
Falcon 50mL Conical Centrifuge Tubes	14-432-22	Fisher Scientific
MilliporeSigma™ Goat anti-Human IgA alpha-Chain Specific Peroxidase Secondary Antibody	40-113-52ML	Millipore Sigma
Spike protein	Z03501	Genscript

STEPS MATERIALS

NAME	CATALOG #	VENDOR
Spike protein	Z03501	Genscript
Tween 20	P7949	Sigma Aldrich
PBS pH 7.4		
Normal Goat Serum	PCN5000	Thermo Fisher
Skim Milk Powder		
BSA		Sigma Aldrich
MilliporeSigma™ Goat anti-Human IgA alpha-Chain Specific Peroxidase Secondary Antibody	40-113-52ML	Millipore Sigma
1-Step™ Ultra TMB-ELISA Substrate Solution	34028	Thermo Fisher
Sulfuric acid	258105	Sigma – Aldrich

MATERIALS TEXT

Reagent or Consumable	Supplier	Catalog #
Text	Text	Numeric
Thermo Scientific 50mL Conical Sterile Polypropylene Centrifuge Tubes	Fisher	12-565-270
Corning™ 96-Well, High Binding, Flat-Bottom, Half-Area Microplate	Fisher	07-200-37
Gibco™ PBS, pH 7.4	Fisher	10-010-023
Tween™ 20, Fisher BioReagents	Fisher	BP337-100
MilliporeSigma™ Calbiochem™ BSA, 30% Sterile-Filtered Aqueous Solution, Preservative-Free	Fisher	12-662-650ML
Gibco™ Normal Goat Serum	Fisher	PCN5000
MilliporeSigma™ Calbiochem™ OmniPur™ Nonfat Powdered Milk	Fisher	62-501-0GM
Pierce 1-Step Ultra TMB-ELISA Substrate Solution	Fisher	PI34028

Sulfuric Acid, Fisher Chemical	Fisher	A510-P500
MilliporeSigma™ Goat anti-Human IgA, alpha-Chain Specific, Peroxidase, Secondary Antibody	Fisher	40-113-52ML
Spike protein (1mg = 1000000ng)	Genscript	Z03501

EQUIPMENT

NAME	CATALOG #	VENDOR
Centrifuge	-	
405 TS Washer	405 TS	Biotek
800 TS Absorbance Reader	800 TS	Biotek


SAFETY WARNINGS

Handle biological samples in a biological containment hood

Milk Processing 1h

- 1 ⌚ 00:15:00
⚙️ 800 x g, 22°C, 00:15:00


Pour fresh milk into 50mL conical tube and centrifuge at 800g for 15min at room temperature

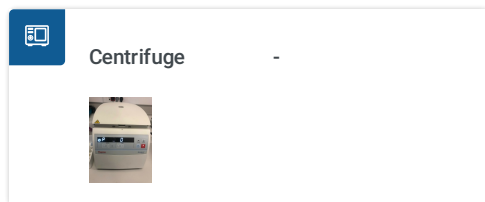

- 2 remove fat using a sterile spatula
- 3 carefully pour skim milk into a clean 50mL tube
- 4 repeat centrifugation and fat removal 2x

Day 1 5m

- 5 Coat half-area 96-well plates with 50ng SARS-CoV-2 spike protein in 50ul PBS [M]0.001 mg/ml
🧴 50 µl

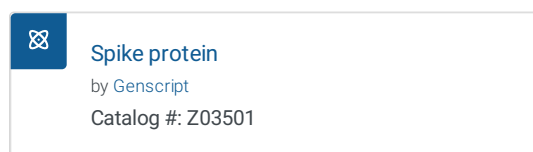


6 over and store overnight at 4C 🕒 **Overnight**

🌡 **4 °C**


Day 2 5h

7 Wash plates with 100ul 0.1% Tween 20/PBS (PBS-T) 6x 📄 **100 µl**

[M] **0.1 % volume**



405 TS Washer
plate washer
Biotek 405 TS



Tween 20
by Sigma Aldrich
Catalog #: P7949




PBS pH 7.4


8 Block plates with 100ul 3% goat serum/0.5% milk powder/3.5% PBS-T for 1h at room temperature 📄 **100 µl**

🕒 **01:00:00**

🌡 **Room temperature**



Normal Goat Serum
by Thermo Fisher
Catalog #: PCN5000



Skim Milk Powder
[View](#)




Tween 20

by Sigma Aldrich

Catalog #: P7949

9 Wash as above

10 Add 50ul milk undiluted or diluted with 1%BSA/PBS to each well in duplicate  50 µl

 1 % volume




BSA

by Sigma Aldrich

11 Incubate at room temperature for 2h  02:00:00

 Room temperature

12 Wash as above

13 Add 50ul  50 µl horseradish peroxidase (HRP)-conjugated secondary antibody: goat anti-human-IgA at 1/3000 in 1%BSA/PBS.



MilliporeSigma™ Goat anti-Human IgA alpha-Chain Specific Peroxidase Secondary Antibody



by Millipore Sigma

Catalog #: 40-113-52ML




BSA

by Sigma Aldrich

14 Incubate at room temperature for 1h  01:00:00  Room temperature


15 Wash as above

16 Add 50u Pierce 1-Step Ultra TMB-ELISA Substrate Solution



1-Step™ Ultra TMB-ELISA Substrate Solution
by Thermo Fisher
Catalog #: 34028


 50 µl



Sulfuric acid
by Sigma – Aldrich
Catalog #: 258105

followed by 50ul 2N hydrochloric acid (HCl)

17 Measure OD at 450nm on an ELISA plate reader



800 TS Absorbance Reader
Plate reader
BioTek 800 TS