



APR 05, 2024

Immunofluorescence Staining and High Content Imaging

David Sirkin¹, Gregory Tracy¹, Hanwen Zhang¹, Lilia Peyton¹, Ada McCarroll¹, Jubao Duan²

¹Endeavor Health/NorthShore University HealthSystem;

²NorthShore University HealthSystem/University of Chicago



David Sirkin
NorthShore

ABSTRACT

This protocol offers a description of immunofluorescence staining and high content imaging of iPSC induced neurons using the ImageXpress Micro Confocal High-Content Imaging System (Molecular Devices, San Jose, CA).

OPEN  ACCESS



Protocol Citation: David Sirkin, Gregory Tracy, Hanwen Zhang, Lilia Peyton, Ada McCarroll, Jubao Duan 2024.

Immunofluorescence Staining and High Content Imaging.

protocols.io

<https://protocols.io/view/immunofluorescence-staining-and-high-content-imagi-dbn22mge>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: Apr 03, 2024

Last Modified: Apr 05, 2024

MATERIALS

96 well optical bottom plate with a polymer base (Fisher Scientific, 12-566-70).

1x PBS, pH 7.4 (Fisher Scientific, 10-010-023)

Sodium Azide 5% (Avantor, BDH7465-2)

Triton X-100 (Sigma, SLCD3244)

Bovine Albumin Fraction V (7.5% solution) (Fisher Scientific, 50-121-5315)

DAPI (1mg/ml) (Fisher Scientific, EN62248)


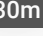
Primary Antibodies

Antibody	Catalogue Number	Concentration
Mouse anti-Synapsin 1	SySy, 106011	1:1,000
Goat anti-tdTomato	Thermofisher, TA150129	1ug/ml
Chicken anti-MAP2	Sigma, AB5543	1:5,000




Secondary Antibodies




Antibody	Catalogue Number	Concentration
Donkey anti mouse Alexa 488	Fisher Scientific, A211202	1:1,000
Donkey anti-chicken Alexa 647	Fisher Scientific, EN62248	1:1,000
Donkey anti-goat Alexa 568	Fisher Scientific, A11057	1:1,000





Immunofluorescence Staining



- 1 Wash iPSC derived neurons twice with 1X PBS and fix with 4% PFA for  00:30:00 in a 96 well optical  bottom plate with a polymer base.




2 Store fixed neurons at 4C in 1x PBS with 0.02% sodium azide.



3 Aspirate sodium azide with pipette and permeabilize neurons in  150 µL of 1x PBS with 0.5% Triton X 15m
100 for  00:15:00 at  Room temperature (RT) without shaking.

4 Aspirate PBST and block neurons with  300 µL of 3% BSA and 0.1% Triton X-100 in 1x PBS for 1h
 01:00:00 at  Room temperature .

5 Aspirate blocking buffer and stain neurons with primary antibodies mouse anti-Synapsin 1 (1:500), goat anti-tdTomato (1ug/ml), and chicken anti-MAP2 (1:5000), in blocking buffer by adding  100 µL of blocking buffer for  01:30:00 at  Room temperature or overnight at  4 °C .

6 Aspirate primary antibody solution and wash samples three times for  00:05:00 with  200 µL 1x 5m
PBS with 0.1% Triton X-100 (0.1% PBST).

7 Aspirate 0.1% PBST and incubate samples with secondary antibodies in blocking buffer by adding  100 µL for  01:00:00 at  Room temperature in the dark. 1h

8 Aspirate secondary antibody solution and wash neurons twice with  200 µL of 0.1% PBST for 5m
 00:05:00 .



9 Aspirate 0.1% PBST and incubate neurons with DAPI (0.5 mg/ml) by adding 100 μ L at

10m



Room temperature for 00:10:00 .

10 Aspirate DAPI and add 250 μ L of 1x PBS containing 0.05% sodium azide.



11 Store plate at 4 $^{\circ}$ C until use.



Image Acquisition

12 Allow plate to warm to RT before imaging.

13 Neurons are imaged using Molecular Devices (San Jose, CA) ImageXpress Micro Confocal High-Content Imaging System at both 20x and 40x.

14 The laser wavelengths used were DAPI, FITC, Texas Red, and Cy5.

15 Each well in the 96 well plate was imaged at 8 sites for 40x and 9 sites for 20x with 8-10 z stacks at 1 μ m step size. For the 40x objective the pixel size is 0.3438 μ m² with a pinhole of 60 μ m, 20x objective pixel size is 0.6842 μ m² also with a pinhole of 60 μ m.

20x Image Analysis

- 16** Acquired images are analyzed as 2D maximum projections. Analyze the first two morphometrics, mean number of branches per cell and mean length of neurite outgrowth per cell, with the built-in Neurite Outgrowth Application Module within the MetaXPress 6 software, version 6.7.2.290
- 17** Both mean number of neurite branches per cell and the mean length of neurite outgrowth per cell were calculated based on using the DAPI stain as a nuclear marker and the tdTomato stain, which labels excitatory neurons, as the neurite and cell body marker. The values used for subsequent analysis can be optimized depending on staining quality and cell density.
 - 17.1** Define cell bodies to have an approximate maximum width of 30 μm , a minimum area 300 μm^2 , and a pixel value of at least 1500 above local background level.
 - 17.2** Define nuclei to have an approximate minimum width of 8 μm , an approximate maximum width of 20 μm , and a pixel value of at least 1500 above local background level.
 - 17.3** Define neurite outgrowth to have a maximum width of 2 μm , a minimum projection length of 15 μm from the cell body, and a pixel value of at least 500 above local background level.

40x Image Analysis

- 18** Acquired images are analyzed as 2D maximum projections. Analyze the third morphometric, excitatory synapse density, using a custom synaptic assay module with MetaXpress 6 software. The values used for subsequent analysis can be optimized depending on staining quality and cell density.
 - 18.1** Use Synapsin1 staining to identify puncta with an approximate minimum width of 0.5 μm , an approximate maximum width of 2 μm , and a minimum pixel value of 2500 above local background level.
- 19** The custom module will generate the number and area of Synapsin1 positive puncta within the colocalized MAP2+ and tdTomato+ signals.

19.1

Puncta density is generated by dividing the total area of puncta within the colocalized MAP2 and tdTomato staining by the area of MAP2+ and tdTomato+ signal within the neurites.