



OCT 17, 2023


## 🌐 Masson-Fontana Staining

Nuriapenuelas<sup>1</sup>

<sup>1</sup>Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

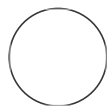

Nuriapenuelas

### ABSTRACT

Masson-Fontana Staining for melanin detection

OPEN  ACCESS



#### DOI:

[dx.doi.org/10.17504/protocols.io.ewov1q3opgr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1q3opgr2/v1)

**Protocol Citation:** Nuriapenuelas 2023. Masson-Fontana Staining. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.ewov1q3opgr2/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Oct 17, 2023

**Last Modified:** Oct 17, 2023

- 1 Heat paraffin tissue sections at 60°C for 10min.
- 2 Incubate sections in xylene (5 min, 3 times) and in ethanol serial washes (100-95-70%-H<sub>2</sub>O, 5 min each).
- 3 Incubate the sections with ammoniac solution for 40 min at 56°C. Rinse sections in distilled water.
- 4 Incubate the sections with sodium thiosulphate for 2 min. Rinse sections in distilled water.
- 5 Finally, counterstain with Kernechtrot for 7 min. Rinse sections in distilled water.