



AUG 10, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.n2bvj384plk5/v1

Protocol Citation: rund.tawfiq 2023. XRF Soil Sample Preparation. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n2bvj384plk5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working


Created: Aug 10, 2023

Last Modified: Aug 10, 2023

XRF Soil Sample Preparation

rund.tawfiq¹

¹CBRC, KAUST

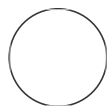

rund.tawfiq

ABSTRACT

This protocol outlines a procedure for the preparation of both wet and dry soil samples for X-ray Fluorescence (XRF) inorganic chemical analysis. The accurate analysis of soil samples is essential for understanding elemental composition and its implications for various environmental studies. Proper preparation ensures reliable and reproducible results.

For wet soil samples, the protocol encompasses drying to remove excess moisture, preventing potential interference with subsequent XRF analysis.

For wet and dry soil samples, the protocol outlines grinding to ensure homogeneity and fine particle size. Proper grinding enhances the accuracy of XRF analysis by promoting even distribution and facilitating the release of elemental components.

1 Weigh ~500 mg of soil

30s

Step 1 includes a Step case.

Dry Soil

Wet Soil

step case

Dry Soil

Dry soil is loose and crumbly to the touch and doesn't hold its shape when squeezed.

2 Grind with a mortar and pestle until the soil turns to powder (1-2 mins)