

Jul 25, 2024 Version 1

T cell purification and activation V.1

DOI

dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v1

Moustafa Nouh Elemeery¹, Salix Boulet²

¹Department of neurosciences, Faculty of Medicine, Université de Montréal, Canada. Medical Biotechnology Department, National Research Centre, Dokki, Cairo, Egypt;

²Centre de recherche de l'hôpital Maisonneuve-Rosemont (CRHMR), Department of medicine, Faculty of Medicine, Université de Montréal

ASAP Collaborative Rese...



Lilia Rodriguez

Université de Montréal

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v1

Protocol Citation: Moustafa Nouh Elemeery, Salix Boulet 2024. T cell purification and activation. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.81wgbz431gpk/v1>

Manuscript citation:

Elemeery, M. N., Tchung, A., Boulet, S., Mukherjee, S., Giguere, N., Daudelin, J. F., ... & Trudeau, L. E. (2024). Adoptive transfer of mitochondrial antigen-specific CD8+ T-cells in mice causes parkinsonism and compromises the dopamine system. *bioRxiv*, 2024-02.

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: July 02, 2024

Last Modified: July 25, 2024

Protocol Integer ID: 103033



Keywords: ASAPCRN

Funders Acknowledgement:
Aligning Science Across
Parkinson's (ASAP)
Grant ID: ASAP 000525

Abstract

This protocol details the purification and activation of Mouse Naïve CD8+ T Cell using an isolation kit from STEMCELL. This kit is designed to isolate naïve CD62L+CD44-CD8+ T cells from single-cell suspensions of splenocytes by negative selection. Unwanted cells are targeted for removal with biotinylated antibodies that are directed against non-naïve CD8+ T cells (CD4, CD11b, CD11c, CD19, CD44, CD45R/B220, CD49b, TCR γ/δ , TER119) and streptavidin-coated magnetic particles. Labeled cells are separated using an EasySepTM magnet without the use of columns.

Materials

Reagents and solution:

A	B	C
Reagent or solution	Supplier	Catalogue #
Anti-mouse CD3 (Clone: 145-2C11)	BioXcell	BE0001
Anti-mouse CD28 (clone: 37.51)	Leinco	C379-5.0 mg
RPMI media	VWR (Corning)	CA45000-396
Phosphate Buffer saline (PBS)	Gibco	14190144
EasySep buffer	STEMCELL	20144
EasySep™ Mouse CD8+ T Cell Isolation Kit	STEMCELL	19858
ACK (Ammonium-Chloride-Potassium) Lysing Buffer	ThermoFisher	A1049201
L-glutamine 200 mM	VWR	CA45000-676
HEPES 1M	Fisher	MT25060CI
Sodium Pyruvate	VWR	CA45000-710
2-Mercaptoethanol	ThermoFisher	21985023
Non-Essential Amino acids	VWR (Corning)	CA45000-700
Fetal Bovin Serum	Gibco	12483020
24 wells plate flat bottom suspension plates	Sarstedt	83.1836.500
96 wells plate flat bottom suspension plates	Sarstedt	82.1581.001
Cell strainers (70um)	Fisher	08-771-2

⊗ InVivoMAb anti-mouse CD3 **BioXcell Catalog #BE0001**

⊗ Anti-Mouse CD28 (Clone 37.51) **Leinco Catalog #C379**

⊗ RPMI **VWR International Catalog #45000-396**

⊗ Phosphate buffered saline (PBS) without Ca/Mg **Thermo Fisher Scientific Catalog #14190144**

⊗ EasySep™ Buffer **STEMCELL Technologies Inc. Catalog #20144**

⊗ EasySep™ Mouse Naïve CD8+ T Cell Isolation Kit **STEMCELL Technologies Inc. Catalog #19858**

⊗ ACK Lysing Buffer **Thermo Fisher Scientific Catalog #A1049201**

⊗ L(+)-Glutamine solution 200 mM **VWR International Catalog #45000-676**

⊗ Sodium pyruvate solution 100 mM **VWR International Catalog #45000-710**

⊗ 2-mercaptoethanol **Gibco - Thermo Fisher Catalog #21985023**

⊗ Non essential amino acids **VWR International Catalog #45000-700**

⊗ Fetal Bovine Serum, qualified, Canada **Thermo Fisher Catalog #12483020**

⊗ Falcon™ Cell Strainers **Fisher Scientific Catalog #08-771-2**

**RPMI complete (RPMIc):**

A	B
RPMI	500 mL
FBS 10%	50 mL (decomplemented)
L-Glutamine	5 mL
Sodium pyruvate	5 mL
Antibiotic (Pen-Strep)	5 mL
Non-essential amino acids	5 mL
2-Mercaptoethanol	50 mmol/L (final)


















Note

Note! Very important, 2-Mercaptoethanol is an essential growth factor for mouse T-lymphocytes.







Purification and activation

4d 0h 30m

- 1 Dilute CD3 antibody (145-2C11) (clone KT3 can also be used) to  1 undetermined in PBS.
- 2 Coat plates with anti-CD3 antibody.
- 3 Use 24 or 96 wells flat bottom suspension plates. If these plates are not used, there is a risk of partial stimulation due to low absorbance of the antibody on the plate.
- 3.1 Add  100 μL of antibody/well if 96 wells or  1 mL per well if 24 wells. 
- 3.2 Incubate the plate for  24:00:00 at  4 °C or  01:00:00 at  37 °C . 
1d
- 4 Remove the antibody (aspirate) and wash 2 times with PBS. 
- 5 Add  100 μL of RPMIc if 96 wells or  1 mL if 24 wells and incubate the plate at  37 °C for  00:15:00 . 
15m
- 6 Collect spleen from mice 6-8 weeks in complete RPMI media (RPMIc), sex matched with recipient mice.
- 7 Purify CD8+ T cells using EasySep mouse naïve purification kit (STEMCELL, catalo # 19858) as follow:
- 7.1 Use a frosted microscope slide to homogenize spleens in PBS or Hanks' Balanced Salt Solution (HBSS) containing 2% fetal bovine serum (FBS).
- 7.2 Remove aggregates and debris by passing cell suspension through a 70 μm mesh nylon strainer. Collect cells in a 15 mL tube.
- 7.3 Centrifuge at  1300 rpm, 00:05:00 and discard the supernatant. 
5m




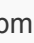
7.4 Red blood cells lysis is done by adding  5 mL /spleen of 0.83% ammonium chloride and incubate for  00:05:00 (or  00:02:00 for ACK lysing buffer) at  Room temperature while continuously shaking tubes.



5m



7.5 Quench by filling up tube with RPMIc.

7.6 Centrifuge at  1300 rpm,  00:05:00 and discard supernatant.


5m





7.7 Resuspend in EasySep buffer at 1×10^8 cells/ml (Easysep mouse naïve CD8+ T cell isolation kit (#19858A) and follow the protocol provided by STEMCELL.


8 Wash cells with RPMIc medium.





9 Resuspend cells in  5 mL of RPMIc, count cells and adjust the concentration to 2×10^6 cells/mL

10 Add  100 μ L of cells (2×10^5 cells) per well if 96 wells or  1 mL of cells (2×10^6 cells) per well if 24 wells).



11 Add purified anti-CD28 antibodies to reach a concentration of  5 undetermined .





12 Incubate the plate at  37 °C and 5% CO₂ for  24:00:00 .

1d



13 Add 20 Unit/ml of IL-2.



14 Incubate the plate at  37 °C and 5% CO₂ for another  48:00:00 .

2d



**Note******* Don't forget the non-stimulated controls.**

- Resuspend the cells well throughout the experiment (the cells quickly settle to the bottom of the tube)
- Check purity of CD8+ T cells after STEMCELL isolation by staining with anti-CD8 antibody followed by flow cytometry.
- Check activation after stimulation using anti-CD8 and anti-CD44 staining followed by flow cytometry.