




APR 02, 2023

Alpha-synuclein immunohistochemistry on STC-1 cells using DAB

Michael J Hurley^{1,2}

¹Department of Clinical and Movement Neurosciences, UCL Queen Square Institute of Neurology, London, UK;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, Maryland, USA.

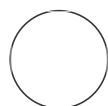

Michael J Hurley

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.kqdg39yjpg25/v1

Protocol Citation: Michael J Hurley 2023. Alpha-synuclein immunohistochemistry on STC-1 cells using DAB. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kqdg39yjpg25/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Feb 27, 2023

Last Modified: Apr 02, 2023

PROTOCOL integer ID:
 77692

Keywords: ASAPCRN, STC-1, immunohistochemistry, alpha-synuclein

ABSTRACT

This protocol describes how to visualise alpha-synuclein in STC-1 cells by DAB immunohistochemistry. It also works for other antibodies (e.g. 5-HT, CCK, GLP).

MATERIALS

Rabbit anti-alpha-synuclein	ab212184	Abcam
Rabbit anti-beta-actin	ab241153	Abcam
Goat anti-rabbit IgG biotin	111-065-144-JIR	Strattech
Streptavidin	POD11089153001	Roche

Alpha-synuclein immunohistochemistry on STC-1 cells using DAB 1d 0h 58m

1 Grow STC-1 cells on 12 mm glass coverslips treated with poly-D-lysine or Geltrex™.

2 Fix cells with 4 % paraformaldehyde ⌚ 00:15:00

15m

Safety information

Paraformaldehyde is toxic by inhalation

[158127 \(sigmaaldrich.com\)](https://www.sigmaaldrich.com)

3 Wash cells with PBS 2 x ⌚ 00:05:00

5m

4 Aspirate and discard PBS and freeze cells < 🌡 -20 °C

5 Thaw frozen 24-well plates containing fixed cells grown on glass coverslips

6 Wash cells with PBS ⌚ 00:05:00









5m


7 Quench with 0.3% H₂O₂ ⌚ 00:05:00

5m

8 Wash cells with PBS 2 x ⌚ 00:05:00

5m

- 9 Block with 10 % goat serum in PBS containing 0.005 % Triton™ X-100 and 0.05 % thimerosal  01:00:00 1h
- 10 Incubate with rabbit anti-alpha-synuclein (1:500) (Abcam, ab212184) in 10 % goat serum containing 0.005 % Triton™ X-100 and 0.05 % thimerosal (block buffer)  20:00:00 20h
- 11 Wash with PBS 4 x  00:05:00 5m
- 12 Incubate with biotinylated goat anti-rabbit IgG (1:500) (Stratech, 111-065-144-JIR) in 2 % goat serum containing 0.005 % Triton™ X-100 and 0.05 % thimerosal (antibody buffer)  02:00:00 2h
- 13 Wash with PBS 4 x  00:05:00 5m
- 14 Incubate with 1:250 streptavidin-HRP (Roche, 11089153001) for 1 hour in antibody buffer  01:00:00 1h
- 15 Wash with PBS 6 x  00:05:00 5m
- 16 Visualise staining by incubation with PBS containing 0.05 % (w/v) 3,3'-Diaminobenzidine/0.015 % (v/v) H₂O₂/0.05 % (w/v) nickel ammonium sulphate for ≤  00:03:00 3m

17 Rinse for >  00:05:00 in running tap water

5m

18 Remove coverslips from wells and mount on microscope slides with Hydromount™ (National Diagnostics)