



NOV 02, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.4r3l22er4l1y/v1](https://dx.doi.org/10.17504/protocols.io.4r3l22er4l1y/v1)

**Protocol Citation:** Alexander Röntgen 2023. SDS-PAGE.

**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.4r3l22er4l1y/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working


**Created:** Nov 02, 2023

**Last Modified:** Nov 02, 2023

## SDS-PAGE

Alexander  
 Röntgen<sup>1</sup>

<sup>1</sup>University of Cambridge

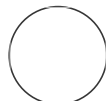

Alexander Röntgen

### ABSTRACT

Protocol for running an SDS-PAGE.

### MATERIALS

- NuPAGE Sample Buffer 4X Concentrate (ThermoFisher Scientific)
- Heating block
- XCell *SureLock*<sup>TM</sup> Mini-Cell (ThermoFisher Scientific)
- 4-12% Bis-Tris gels (ThermoFisher Scientific)
- NuPAGE MES SDS Running Buffer (ThermoFisher Scientific)
- SeeBlue Plus2 Pre-stained Protein Standard (ThermoFisher Scientific)
- Instant Blue Protein stain (Abcam)
- ChemiDoc MP Imaging system (Bio-Rad)

## Sample preparation

- 1 Dilute NuPAGE Sample Buffer 4X Concentrate into protein samples
- 2 Heat for 10 min at 95 °C

## Gel

- 3 Run gel at ~150 V for 35 min
- 4 Incubate in staining solution until bands appear
- 5 Image using ChemiDoc MP Imaging system