

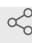


Oct 31, 2022

Immunocytochemical analysis

 In 1 collection**Felix Kraus**^{1,2}¹Department of Cell Biology, Blavatnik Institute, Harvard Medical School, 240 Longwood Ave, Boston MA 02115, USA;²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD 20815, USA

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.bp2l69bmlqe/v1 Felix Kraus

ABSTRACT

Protocol for the immunocytochemical analysis of cell culture cells.

DOI

dx.doi.org/10.17504/protocols.io.bp2l69bmlqe/v1

PROTOCOL CITATION

Felix Kraus 2022. Immunocytochemical analysis. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bp2l69bmlqe/v1>



COLLECTIONS

 **Kraus et al., 2022 FBX07 /Park15**

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 11, 2022

LAST MODIFIED

Oct 31, 2022

PROTOCOL INTEGER ID

71172


Seeding of HeLa

- 1 Wash HeLa cells with 1x PBS
- 2 Add Trypsin to cells for 5 min and incubate at 37°C to dissociate cells from plastic well
- 3 Resuspend cells in 1 mL DMEM media
- 4 Count cells
- 5 Seed appropriate number of cells into 6-/12-/24-well glass bottom dish
- 6 Top up glass bottom dish with either 3-/1.5/1 mL DMEM and place cells back into incubator
- 7 The next day exchange DMEM with DMEM + 2µg/ml doxycycline for 18h to induce Parkin expression.
- 8 Induce mitophagy using Antimycin A / Oligomycin A for the desired time.

Staining

- 9 Aspirate DMEM and fix cells in 1 ml pre-warmed 4% PFA for 30 min.

- 10 Aspirate PFA solution and wash wells 3x with PBST (1x PBS, 0.02% Tween 20)
- 11 Permeabilize the cells by adding 0.2% Triton X-100 in PBS.
- 12 Remove the detergent solution by aspiration. Wash wells 3x with PBST (1x PBS, 0.02% Tween 20). Drain well.
- 13 Block cells for 10 min with 3% BSA – 1x PBS.
- 14 Remove BSA solution by aspiration. Wash wells 3x with PBST (1x PBS, 0.02% Tween 20). Drain well.
- 15 Incubate with primary antibodies in 3% BSA - 1x PBS for 3h at RT or overnight at 4°C with gentle shaking.
- 16 Wash wells 3x with PBST (1x PBS, 0.02% Tween 20). Drain well.
- 17 Incubate with secondary antibodies in 3% BSA - 1x PBS for 45 min – 1h.
- 18 Wash wells 3x with PBST (1x PBS, 0.02% Tween 20). Drain well.
- 19 If required, add Hoechst33342 or DAPI 1:2000 to wells for 5 min with gentle shaking.
- 20 Wash wells 3x with PBST (1x PBS, 0.02% Tween 20). Drain well.

- 
- 21 Exchange PBST with 1x PBS and keep cells at 4°C until imaging. Image within the next few days.