



VERSION 2
DEC 15, 2022

WORKS FOR ME

1

Determining biofilm growth amount (absorbance) V.2

COMMENTS 0

DOI

dx.doi.org/10.17504/protocols.io.n2bvj8x95gk5/v2

An.Huang¹

¹XJTLU



An.Huang

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

ABSTRACT

This protocol describes a method to determine the growth amount of biofilm at the early stage of biofilm formation by measuring absorbance. Here we use our own engineered bacteria, and it requires induction of IPTG and cultured with silver ion. You may use a different kind of bacteria and different culturing medium.

DOI

dx.doi.org/10.17504/protocols.io.n2bvj8x95gk5/v2

PROTOCOL CITATION

An.Huang 2022. Determining biofilm growth amount (absorbance) . **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n2bvj8x95gk5/v2>
Version created by An.Huang

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Dec 15, 2022

LAST MODIFIED

Dec 15, 2022

PROTOCOL INTEGER ID

74006

MATERIALS TEXT

(Optional) IPTG, silver nitrate, *Escherichia coli*, LB broth medium.

Droppers, forceps, flasks with caps, 50 mL centrifuge tubes, MBBR carrier K1, oven, 1X phosphate-buffered saline (PBS), centrifuge, spectrometer.

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

IPTG induction

- 1 *Escherichia coli* grown overnight was diluted by LB to OD₆₀₀=0.4-0.6.
- 2 IPTG was added to cell culture to 1mM IPTG finally, and incubated 3h at 171 rpm, 37°C in orbital shaking incubator.

Sample preparation

- 3 Preparing several 100mL flasks by filling the flasks with MBBR carrier K1. Autoclave all the flasks and dry them in an oven at 60s°C.

Note

- 4 Cells cultures and silver nitrate (to 6 μ M) solution were mixed in advance. Add the mixture into each bottle.

Note

- 5 Incubate all bottles in a biochemical incubator at 37 °C .

Absorbance measurement

10m

- 6 At the day of measuring the absorbance, take the flasks to be tested out from incubator.

- 7 Using a pair of forceps, carefully pick out the carriers that has not been submerged into the culture and discard.

Note

- 8 Randomly pick five carriers out from culturing flask and put into a 50 mL centrifuge tube. Label this tube as "Washing tube".

Note

- 9 Wash these five carriers using 1X phosphate-buffered saline (PBS) for three times. Collect all eluents.

Note

- 10 Use a dropper, add 2-3 mL 1X PBS into the Washing tube. Cap the tube and rinse every corner of the inside wall of tube using the PBS added by rotating the tube. Collecting the eluents by pouring these PBS into Eluent tube. Repeat these steps for three to five times.

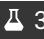
Note


- 11 Centrifuge the Eluent tube at  10000 rpm, Room temperature, 00:10:00 .

10m

- 12 Discard all supernatants.

Note

- 13 Add  3 μ L 1X PBS into Eluent tube. Resuspend the pellet.

- 14 Measure the absorbance of the liquid  13 at OD600. Record the absorbance.
- 15 Repeat steps from Step 7 to Step 14 for three times. Calculate the mean value for the three measurements. A growing curve of biofilm can be generated from the data collected.