




APR 03, 2023

Endogenous coimmunoprecipitation

In 2 collections

michela.deleidi¹, Pascale Baden¹

¹German Center for Neurodegenerative Diseases (DZNE), Tübingen, 72076 Germany

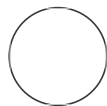

Federico Bertoli

ABSTRACT

Protocol used for immunoprecipitation of HSP60 and LONP1 in HEK cells to show the interaction with V5-Flag-tagged WT-GCase

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.eq2ly7bqmlx9/v1

Protocol Citation: michela.deleidi, Pascale Baden 2023. Endogenous coimmunoprecipitation. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.eq2ly7bqmlx9/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Apr 01, 2023

Last Modified: Apr 03, 2023

PROTOCOL integer ID:
 79871

Endogenous coimmunoprecipitation

1 Wash Protein G agarose fast-flow beads in TBS 1X + 0.05% NP40. To this

end, the beads were incubated with lysis buffer on the spinning wheel (25 RPM) for 2 –5 min at 4 °C, followed by a 1 min centrifugation at 500 RPM 4 °C.

- 2 Repeat step 1 for a total of 3 washes.
- 3 Incubate 20ul of pure beads with 5ug of anti-LONP1 antibody or 5ug of normal rabbit IgG as control (or 3ug of anti-Hsp60 and 3ug of mouse IgG as control) and add 300ul of washing buffer to each tube.
- 4 Incubate 2h at 20RPM on a rotating wheel at 4°C
- 5 Meanwhile, detach HEK cells using Accutase for 5 minutes at 37°C and collect them.
- 6 Detach cells using Accutase for 5 minutes at 37°C and collect them.
- 7 Spin cells in a centrifuge at 250g for 5 minutes at room temperature.
- 8 Remove the supernatant and wash cells in PBS.
- 9 Repeat steps 7 and 8 for a total of 2 washes.

- 10** Lyse cells in 1% TBS + 0.5% NP40 + PI/PHI (Pierce #A32959)
- 11** After 2h incubation, wash beads again 3 times in washing buffer.
- 12** Incubate antibody-coated beads with 3.7 ug of lysate on the spinning wheel for 2h.
- 13** After 2h incubation, wash beads again 3 times in washing buffer.
- 14** Elute by boiling the beads twice with 2x Laemmli buffer at 95°C for 8min at 300rpm in a thermoblock.
- 15** Spin beads for 1 minute at 10000 RCF.
- 16** Collect supernatant and proceed with western blot analysis.