



SEP 20, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.4r3l229ojl1y/v1](https://dx.doi.org/10.17504/protocols.io.4r3l229ojl1y/v1)

**Protocol Citation:** ghanem.elkassem, micboe 2023. Rapid and robust cloning of sgRNA expression plasmids. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.4r3l229ojl1y/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Sep 12, 2023

## Rapid and robust cloning of sgRNA expression plasmids

ghanem.elkassem<sup>1</sup>, micboe<sup>1</sup>

<sup>1</sup>Medical Faculty of the Martin Luther University Halle-Wittenberg

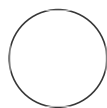

ghanem.elkassem

### ABSTRACT

In this lab protocol, we will outline a step-by-step procedure for cloning of spCas9 sgRNAs using oligo annealing and T4 ligation. This protocol will guide you through the crucial steps of designing, annealing, ligating oligonucleotides and digested plasmids, as well as heat shock transformation procedure to clone sgRNAs expression vectors.

### ATTACHMENTS

[sgLenti vector for CRISPR](#)  
[sgRNA expression.gbk](#)

### PROTOCOL MATERIALS

✕ Macherey-Nagel™ NucleoSpin™ Gel and PCR Clean-up Kit Fischer  
 Scientific Catalog #11992242

Step 10

✕ NucleoSpin Plasmid Transfection-grade, Mini kit for ultrapure plasmid  
 DNA Macherey-Nagal Catalog #740490.250

Step 26

✕ T4 DNA Ligase Buffer (10X) Thermo  
 Fisher Catalog #B69

Step 2

✕ AarI (2 U/μL) Thermo  
 Fisher Catalog #ER1581

Step 4

✕ T4 DNA Ligase (5 U/μL) Thermo  
 Fisher Catalog #EL0011

Step 13

✕ DNA Gel Loading Dye (6X) Thermo Fisher  
 Scientific Catalog #R0611

Step 7

### SAFETY WARNINGS



Wear gloves and a UV filter mask and ensure sleeves are completely covering wrists to avoid direct UV light exposure

## Oligo Design

- 1 Order **TOP** and **BOTTOM** strand oligos, where the TOP sequence is the desired 20 nt sgRNA 'spacer' sequence and the BOTTOM sequence is the reverse complement of the TOP strand. After annealing of both oligos, the 5' overhangs should result in compatible sticky ends for T4 ligation into the target vector. Shown are example sticky ends to clone into the vector sgLenti (Addgene #105996). Alternative sgRNA expression vectors may require different sticky ends.

5'-**TTGG**NNNNNNNNNNNNNNNNNNNNNN -3'  
3'-NNNNNNNNNNNNNNNNNNNNNN**CAAA**-5'

**TOP** strand oligo: 5'-**TTGG**NNNNNNNNNNNNNNNNNNNNNN-3'

**BOTTOM** strand oligo: 5'-**AAAC**NNNNNNNNNNNNNNNNNNNNNN-3'

## Oligo Annealing



1h 40m

- 2 Set up the annealing reaction:

10m

1  1  $\mu$ L Top strand oligo 1M 100 micromolar ( $\mu$ M)

1  1  $\mu$ L Bottom strand oligo 1M 100 micromolar ( $\mu$ M)

1  1  $\mu$ L  T4 DNA Ligase Buffer (10X) Thermo  
Fisher Catalog #B69

1  7  $\mu$ L ddH<sub>2</sub>O

- 3 Anneal oligos by running the following thermocycler program:

1h 30m

1.  95 °C for 5 min

2. Ramp down at  0.1 °C /sec from  95 °C to  25 °C

## Plasmid Restriction Digest


1d

- 4 Set up AarI plasmid digest:

10m

1  1  $\mu$ g vector for spCas9 sgRNA expression

🧪 2 µL 10x AarI Buffer

🧪 1 µL  AarI (2 U/µL) Thermo  
Fisher Catalog #ER1581

🧪 0.4 µL 50x Oligo  
ddH<sub>2</sub>O to 🧪 20 µL

5 Incubate in a Thermocycler at 🔥 50 °C ⌚ Overnight


16h



## Digested Plasmid Purification by Agarose Gel Electrophoresis 3h

6 Prepare 🧪 1 % Agarose gel

30m

7 Add 🧪 4 µL of 6x  DNA Gel Loading Dye (6X) Thermo Fisher  
Scientific Catalog #R0611  
to 🧪 20 µL of the restriction digestion mix.

10m

8 Load samples on the 🧪 1 % Agarose gel and run the electrophoresis at 120 V for 50 min

1h

9 Place the gel in a UV light box


10m

### Safety information

Wear gloves and a UV filter mask and ensure sleeves are completely covering wrists to avoid direct UV light exposure

10 Cut the digested plasmid band from the gel and purify the DNA using the

1h

 Macherey-Nagel™ NucleoSpin™ Gel and PCR Clean-up Kit Fischer  
Scientific Catalog #11992242

11 Determine the concentration of the purified plasmid DNA (e.g. NanoDrop One)

10m

## Ligation of Plasmid and Annealed Oligos

30m

12 Dilute the annealed oligos 1:200 in ddH<sub>2</sub>O

10m

13 Set up the ligation reaction:

10m

50 ng purified vector

1 µL diluted oligos

1 µL 10x T4 Ligation buffer

1 µL T4 DNA Ligase (5 U/µL) Thermo  
Fisher Catalog #EL0011

ddH<sub>2</sub>O to 10 µL

### Note

As a control set up a ligation reaction without oligos

14 Incubate for 10 min at Room temperature

10m













## Heat Shock Transformation of the Ligation Product into Chem...

1d

15 Take DH5a chemocompetent cells from -80 °C storage and thaw on ice for 30 min. Use 1



30m

50 µL aliquot for the ligation reaction and 1 50 µL aliquot for the control reaction

- 16 Add  2  $\mu\text{L}$  of the T4 Ligation reaction and the control reaction to the cells and mix gently by flicking the tube
- 17 Incubate on ice for 30 min 30m
- 18 Incubate the tube containing the cells in a water bath at  42 °C for 45 sec and place back on ice for 2 min 2m 45s
- 19 Add  0.5 mL of LB medium to the cells and incubate in a thermo-block at  37 °C with shaking at 250 rpm for 1 hr 1h
- 20 Plate  100  $\mu\text{L}$  of the cells on  Room temperature LB agar plates with  100  $\mu\text{g/mL}$  Carbenicillin 10m
- 21 Incubate at  37 °C  Overnight 16h
-   

- 22 After the overnight incubation, count the colonies on the ligation reaction and the control plates. The ligation reaction plate should have a minimum of 10x more colonies. 10m
- 23 Seal the plates with parafilm and store at  4 °C for up to 2 weeks

## Plasmid DNA Purification


1d

- 24 Pick up to 3 colonies from the plate using a pipette tip and inoculate  5 mL LB medium with  100 µg/mL Carbenicillin in 15 mL cell culture tubes

- 25 Incubate the cell culture tubes  Overnight at  37 °C and 225 rpm in a shaking incubator

16h



- 26 Extract the plasmid DNA using the  NucleoSpin Plasmid Transfection-grade, Mini kit for ultrapure plasmid DNA Macherey-Nagel Catalog #740490.250

1h 30m

- 27 Measure the concentration of the purified plasmid DNA (e.g. NanoDrop One)

- 28 Send the plasmids for Sanger sequencing with the following primer to confirm the correct sequence:

GGCTTGGATTTCTATAACTTCGTATAGCA

### Note

Provided primer is complementary to the 'sgLenti vector for CRISPR sgRNA expression' provided in the attachments and needs to be adjusted for alternative vectors.