




MAY 17, 2023

Adding 1-hydroxyphenazine to L4 *C. elegans* in liquid culture

Man
Muhammad Zaka Asif¹, Shah¹, Yosef Smadi¹

¹Edison Lab, UGA

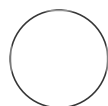

Muhammad Zaka Asif

ABSTRACT

This protocol describes day 4 of our workflow to grow worms in liquid culture to induce the production of natural products (in this case, 1-HP derivatives). In this protocol, we add the toxin to our worms grown in liquid culture to induce the production of glycosylated products.

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.yxmvm28png3p/v1

Protocol Citation: Muhammad Zaka Asif, Man Shah, Yosef Smadi 2023. Adding 1-hydroxyphenazine to L4 *C. elegans* in liquid culture.

protocols.io
<https://dx.doi.org/10.17504/protocols.io.yxmvm28png3p/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: May 17, 2023

Last Modified: May 17, 2023

PROTOCOL integer ID:
82036

- 1 Follow steps 1-6 from the protocol for transferring *C. elegans* to S-buffer, but centrifuge at 525 RCF for 1 minute

- 2 Dilute to a final concentration of 30,000 worms/mL in S-basal and final toxin concentration of 22.3 μ M 1-HP.
- 3 Repeat steps 1-2 for 1.1% DMSO and bacteria-only controls.
- 4 If using negative controls, repeat steps 1-3 for them.
- 5 Put the flasks back in the shaker and allow to incubate for 24 hours.