

MAR 29, 2024

# Metabarcoding Fecal Swabs or Stomach Contents for Fish and Crustaceans using 2-PCR protocol and Illumina MiSeq

Eldridge Wisely<sup>1</sup>

<sup>1</sup>University of Arizona



Eldridge Wisely  
University of Arizona

## ABSTRACT

This protocol describes a method to metabarcode a 170bp region of the mitochondrial 16S rRNA gene of crustaceans and a 163-185bp region of the mitochondrial 12S rRNA gene of fishes. These regions are subjected to PCR separately in multiple replicates and the resulting PCR products are pooled by sample and then indexed for sequencing on an Illumina MiSeq platform.

OPEN ACCESS



DOI:

[dx.doi.org/10.17504/protocols.io.ewov1qxokgr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1qxokgr2/v1)

**Protocol Citation:** Eldridge Wisely 2024. Metabarcoding Fecal Swabs or Stomach Contents for Fish and Crustaceans using 2-PCR protocol and Illumina MiSeq. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.ewov1qxokgr2/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Nov 19, 2023

**Last Modified:** Mar 29, 2024

**PROTOCOL integer ID:** 91145

## PROTOCOL REFERENCES

Berry, Tina E., Sylvia K. Osterrieder, Dáithí C. Murray, Megan L. Coghlan, Anthony J. Richardson, Alicia K. Grealy, Michael Stat, Lars Bejder, and Michael Bunce. 2017. "DNA Metabarcoding for Diet Analysis and Biodiversity: A Case Study Using the Endangered Australian Sea Lion (*Neophoca cinerea*)." *Ecology and Evolution* 7 (14): 5435–53. <https://doi.org/10.1002/ece3.3123>.

Miya, M., Y. Sato, T. Fukunaga, T. Sado, J. Y. Poulsen, K. Sato, T. Minamoto, et al. 2015. "MiFish, a Set of Universal PCR Primers for Metabarcoding Environmental DNA from Fishes: Detection of More than 230 Subtropical Marine Species." *Royal Society Open Science* 2 (7): 150088. <https://doi.org/10.1098/rsos.150088>.

Glenn, Travis C., Roger A. Nilsen, Troy J. Kieran, Jon G. Sanders, Natalia J. Bayona-Vásquez, John W. Finger, Todd W. Pierson, et al. 2019. "Adapterama I: Universal Stubs and Primers for 384 Unique Dual-Indexed or 147,456 Combinatorially-Indexed Illumina Libraries (iTru & iNext)." *PeerJ* 2019 (10). <https://doi.org/10.7717/peerj.7755>.

16S Metagenomic Sequencing Library Preparation." 2013. Illumina. [https://support.illumina.com/content/dam/illumina-support/documents/documentation/chemistry\\_documentation/16s/16s-metagenomic-library-prep-guide-15044223-b.pdf](https://support.illumina.com/content/dam/illumina-support/documents/documentation/chemistry_documentation/16s/16s-metagenomic-library-prep-guide-15044223-b.pdf)

## IMAGE ATTRIBUTION

Haley Capone

## GUIDELINES

The PCR conditions described here are different from the PCR conditions described by Miya et al., and Berry et al. in their respective publications introducing the primers used here. This difference is due to the use of the Takara High Fidelity PCR EcoDry Premix in this protocol.

MATERIALS

96-well PCR plates  
 Adhesive foil PCR plate covers

1.5mL tubes

Glenn et al. Adapterama I iNext indexing primers A-H and 1-12.

PCR machine

Equipment to run gels  
 optionally: equipment for fluorometric quantification

| Equipment   |       |
|---|-------|
| 96-well Magnetic Rack Separator   | NAME  |
| Magnetic Rack Separator   | TYPE  |
| Sergi Lab Supplies  | BRAND |
| B08134P9RT  | SKU   |
| <a href="https://www.amazon.com/Magnetic-Separator-Protein-Purification-Format/dp/B08134P9RT/ref=asc_df_B08134P9RT/?tag=&amp;linkCode=df0&amp;hvadid=416872221972&amp;hvpos=&amp;hvnetw=g&amp;hvrnd=12953200023550024012&amp;hvpone=&amp;hvptwo=&amp;hvqmt=&amp;hvdev=c&amp;hvdvcmdl=&amp;hvlocint=&amp;hvl ocphy=903024">https://www.amazon.com/Magnetic-Separator-Protein-Purification-Format/dp/B08134P9RT/ref=asc_df_B08134P9RT/?tag=&amp;linkCode=df0&amp;hvadid=416872221972&amp;hvpos=&amp;hvnetw=g&amp;hvrnd=12953200023550024012&amp;hvpone=&amp;hvptwo=&amp;hvqmt=&amp;hvdev=c&amp;hvdvcmdl=&amp;hvlocint=&amp;hvl ocphy=903024</a> | LINK  |

## Equipment

### Magnetic Rack for for 1.5 mL Tubes

NAME

Magnetic Rack for DNA, RNA Purification; for 1.5 mL centrifuge Tubes

TYPE

Sergi Lab Supplies

BRAND

B0BZWXZM22

SKU

[https://www.amazon.com/Magnetic-Rack-Purification-centrifuge-Tubes/dp/B0BZWXZM22/ref=asc\\_df\\_B0BZWXZM22/?tag=hyprod-20&linkCode=df0&hvadid=652498086131&hvpos=&hvnetw=g&hvrnd=6716034042841103246&hvpone=&hvptwo=&hvqmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=9](https://www.amazon.com/Magnetic-Rack-Purification-centrifuge-Tubes/dp/B0BZWXZM22/ref=asc_df_B0BZWXZM22/?tag=hyprod-20&linkCode=df0&hvadid=652498086131&hvpos=&hvnetw=g&hvrnd=6716034042841103246&hvpone=&hvptwo=&hvqmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=9)

L  
I  
N  
K

## PROTOCOL MATERIALS

⊗ Crustacean16S-F **Integrated DNA Technologies, Inc. (IDT) Catalog #custom**

In 2 steps

⊗ Takara High Fidelity PCR EcoDry Premix **Takara Bio Inc. Catalog #639280**

In 2 steps

⊗ Agencourt AMPure XP **Beckman Coulter Catalog #A63880** Step 14.1

⊗ 2x Kapa HiFi Hotstart Readymix **Kapa Biosystems Catalog #KK2602** Step 18

⊗ Nuclease-free water **Integrated DNA Technologies, Inc. (IDT) Catalog #11-05-01-14**

In 4 steps

⊗ MiFish-R **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** Step 1

⊗ Crustacean16S-R **Integrated DNA Technologies, Inc. (IDT) Catalog #custom**

In 2 steps

⊗ Buffer EB **Qiagen Catalog #19086** In 2 steps

⊗ MiFish-F **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** Step 1


## BEFORE START INSTRUCTIONS


Work in a pre-PCR lab, as separated as possible from post-PCR products.


Clean work area with 10% bleach solution before beginning work for the day, then change gloves so that no bleach carryover to your samples or reactions occurs.


### Prepare Primers

- Order metabarcoding primers with diversity spacers and Illumina overhang sequences (Illumina, 2013):  


 MiFish-F **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** (Miya et al., 2015):  
TCGTCGGCAGCGTCAGATGTGTATAAGAGACAGNNNGTCGGTAAACTCGTGCCAGC



 MiFish-R **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** (Miya et al., 2015):  
GTCTCGTGGGCTCGGAGATGTGTATAAGAGACAGNNNCATAGTGGGGTATCTAATCCCAGTTTG



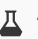
 Crustacean16S-F **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** Berry et al., 2017):  
TCGTCGGCAGCGTCAGATGTGTATAAGAGACAGNNNGGACGATAAGACCCTATA

 Crustacean16S-R **Integrated DNA Technologies, Inc. (IDT) Catalog #custom** (Berry et al., 2017):  
GTCTCGTGGGCTCGGAGATGTGTATAAGAGACAGNNNATTACGCTGTTATCCCTAAAG

We got ours from <https://www.idtdna.com/> as custom oligos at 25nm scale, with standard desalting.

- Reconstitute primers to [M] 100 micromolar ( $\mu\text{M}$ ) stock solutions by adding  40  $\mu\text{L}$  of  

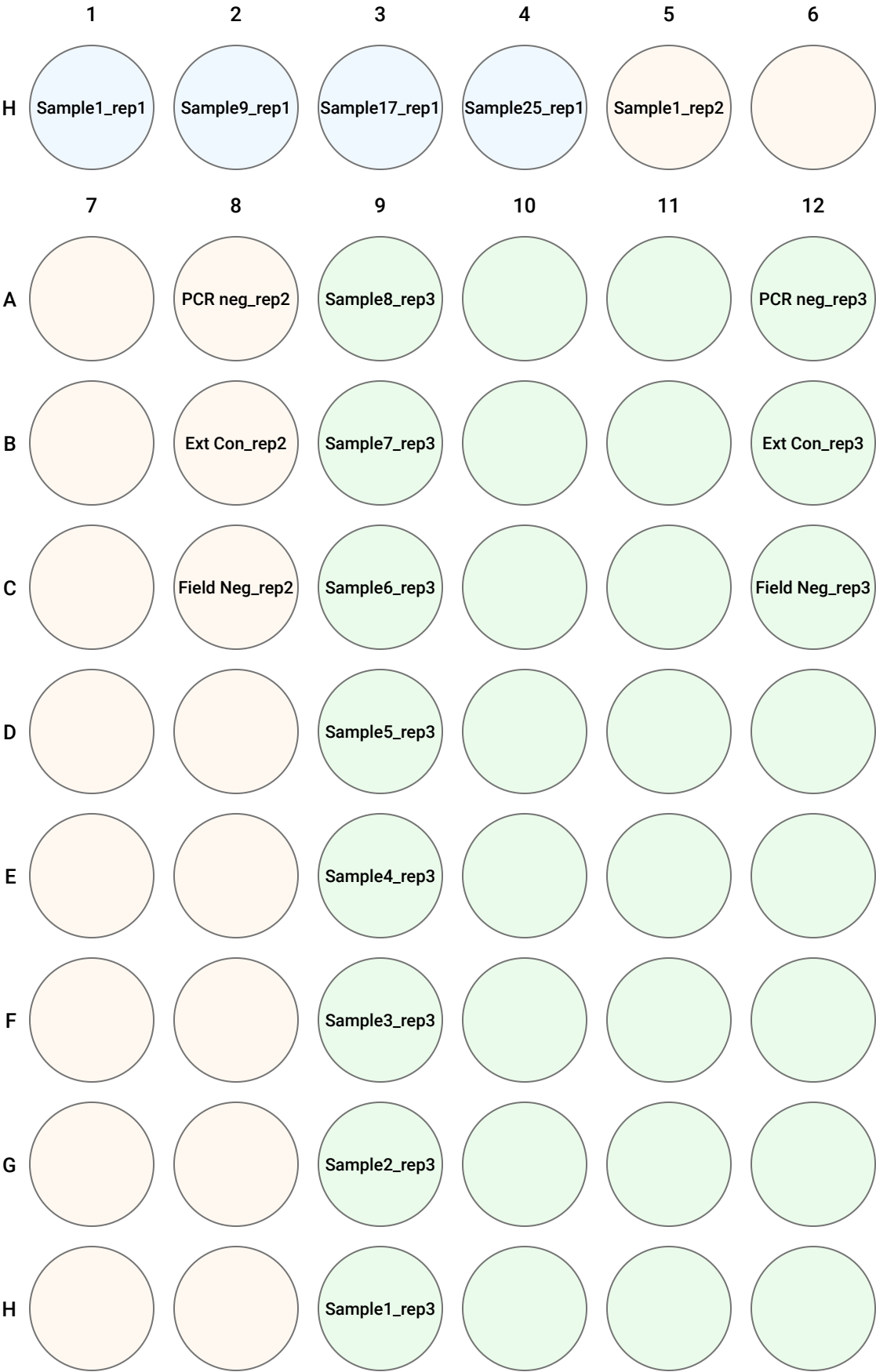
 Nuclease-free water **Integrated DNA Technologies, Inc. (IDT) Catalog #11-05-01-14**
- Make [M] 5 micromolar ( $\mu\text{M}$ ) working solutions of each primer by adding  95  $\mu\text{L}$  of  

 Nuclease-free water **Integrated DNA Technologies, Inc. (IDT) Catalog #11-05-01-14** and  
 5  $\mu\text{L}$  of primer stock solution for each  100  $\mu\text{L}$  of primer that you intend to use within the next week or so.

## Create Plate Map



- 4 Determine which sample will go into each well. This should be the same for each primer set and each replicate. Include at least one extraction control (you can combine aliquots of the extraction controls from each round of DNA extraction into one tube, and use that as your single extraction control), and include a PCR negative control for each plate of PCR. See example below of 21 samples, a field negative sample, a combined extraction control, and a PCR negative.


|   | 1            | 2             | 3             | 4              | 5            | 6 |
|---|--------------|---------------|---------------|----------------|--------------|---|
| A | Sample8_rep1 | Sample16_rep1 | Sample24_rep1 | PCR neg_rep1   | Sample8_rep2 |   |
| B | Sample7_rep1 | Sample15_rep1 | Sample23_rep1 | Ext Con_rep1   | Sample7_rep2 |   |
| C | Sample6_rep1 | Sample14_rep1 | Sample22_rep1 | Field Neg_rep1 | Sample6_rep2 |   |
| D | Sample5_rep1 | Sample13_rep1 | Sample21_rep1 | Sample29_rep1  | Sample5_rep2 |   |
| E | Sample4_rep1 | Sample12_rep1 | Sample20_rep1 | Sample28_rep1  | Sample4_rep2 |   |
| F | Sample3_rep1 | Sample11_rep1 | Sample19_rep1 | Sample27_rep1  | Sample3_rep2 |   |
| G | Sample2_rep1 | Sample10_rep1 | Sample18_rep1 | Sample26_rep1  | Sample2_rep2 |   |



- 4.1** Do not mix sample types between invasively sampled methods (fecal swabs, or stomach contents) and non-invasively sampled methods (eDNA from water or sediment) in the same PCR procedure. And don't plan to sequence both types in the same sequencing run with the combinatorial indexing scheme used here. The potential for contamination of the lower quantity eDNA samples by the higher quantity fDNA samples is too high.

## MiFish Takara PCR Recipe



- 5** Add  24 µL of your MiFish metabarcoding mastermix to each well of  Takara High Fidelity PCR EcoDry Premix **Takara Bio Inc. Catalog #639280**

- 5.1** Add  1 µL DNA extracted from stomach contents or fecal swabs.

- 5.2** Mix and stir together with pipette tip, swirling to make sure the liquid is in the bottom, and bringing any bubbles to the surface of each reaction.


- 5.3** Cap each row of reaction tightly before beginning any other PCR reaction in the same room.

## MiFish Takara PCR Conditions

- 6**  95 °C for  00:01:00


3m 30s

35 cycles of:

 95 °C for  00:00:30

 66 °C for  00:01:00

followed by:

 68 °C for  00:01:00


Hold at  4 °C





## Crustacean\_16S Takara PCR Recipe

### 7 Make your Crustacean\_16S Mastermix:

For each **PCR replicate of each sample** you intend to process (+10% overage), mix:

 2  $\mu$ L  5 micromolar ( $\mu$ M)

 Crustacean16S-F **Integrated DNA Technologies, Inc. (IDT) Catalog #custom**

 2  $\mu$ L  5 micromolar ( $\mu$ M)

 Crustacean16S-R **Integrated DNA Technologies, Inc. (IDT) Catalog #custom**

 20  $\mu$ L  Nuclease-free water **Integrated DNA Technologies, Inc. (IDT) Catalog #11-05-01-14**

For a full plate of 96 reactions, multiply 105.6\*the per-sample volumes in the recipe to make the mastermix.

### 8 Add 24 $\mu$ L of your Crustacean\_16S metabarcoding mastermix to each well of

 Takara High Fidelity PCR EcoDry Premix **Takara Bio Inc. Catalog #639280**



#### 8.1 Add 1 $\mu$ L DNA extract

**8.2** Mix and stir together with pipette tip, swirling to make sure the liquid is in the bottom, and bringing any bubbles to the surface of each reaction.

**8.3** Cap each row of reaction tightly before beginning any other PCR reaction in the same room.


## Crustacean\_16S Takara PCR Conditions



4m



9  95 °C for  00:01:00

4m


35 cycles of:


 95 °C for  00:00:30

 50 °C for  00:01:00

 68 °C for  00:00:30

followed by:

 68 °C for  00:01:00

then hold at  4 °C


## Visualize PCR Products

- 10 Make a 1.7% to 2% agarose gel and run a representative sample of reactions on it to make sure the PCRs worked, producing bands in the 250-300bp range. Check some PCR negatives to see that they don't have bands. Be very careful opening the PCR plate wells at this point to avoid cross-contamination.

## Prepare EtOH for bead cleanup, and bring beads to room temperature

12m 30s

- 11 Get AmpureXP beads out of the refrigerator, and bring to room temp, swirl to mix occasionally, or use a rocking platform.

- 12 Make fresh 80% EtOH so that you will have at least  200 µL of EtOH per well of the combined plate.

- 13 Get 2 sterile DNAase/RNAse free 96-well PCR plates out of their packaging and immediately cover with adhesive foil.



15m

UV clean the plates for  00:15:00

One plate will be for the bead-cleanup steps, and the other will be for the final, cleaned reactions.

## Perform a 1.5x bead cleanup with Ampure XP beads.


**14** in the bead-cleanup plate, do the following steps for one 8-sample row of the plate at a time, pulling back the foil cover for each row after the previous one has been completed.

**14.1** pipette mix  10  $\mu$ L combined PCR product with  15  $\mu$ L

5m

 Agencourt AMPure XP **Beckman Coulter Catalog #A63880**

Incubate  00:05:00 at room temperature.

**15** After the  00:05:00 incubation, place 96-well plate on a

7m

### Equipment

**96-well Magnetic Rack Separator**

NAME

Magnetic Rack Separator

TYPE

Sergi Lab Supplies


BRAND

B08134P9RT
















SKU

[https://www.amazon.com/Magnetic-Separator-Protein-Purification-Format/dp/B08134P9RT/ref=asc\\_df\\_B08134P9RT/?tag=&linkCode=df0&hvadid=416872221972&hvpos=&hvnetw=g&hvrnd=12953200023550024012&hvpone=&hvptwo=&hvqmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=903024](https://www.amazon.com/Magnetic-Separator-Protein-Purification-Format/dp/B08134P9RT/ref=asc_df_B08134P9RT/?tag=&linkCode=df0&hvadid=416872221972&hvpos=&hvnetw=g&hvrnd=12953200023550024012&hvpone=&hvptwo=&hvqmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=903024)

L  
I  
N  
K

for  00:02:00 or until liquid is clear.

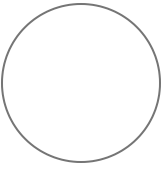
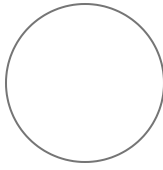
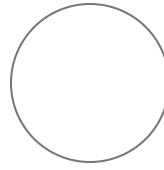
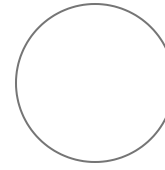
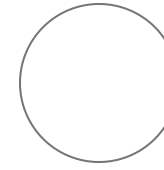
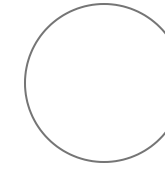
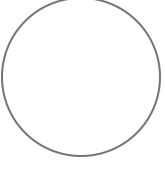
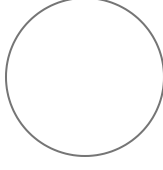
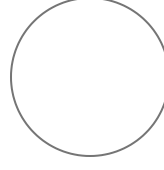
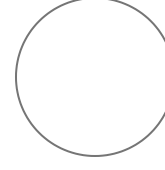
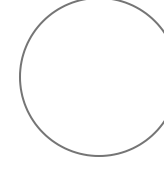
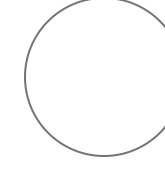
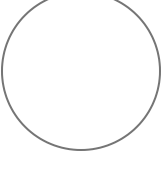
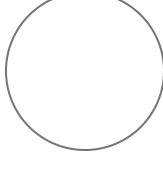
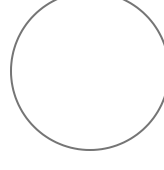
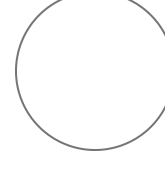
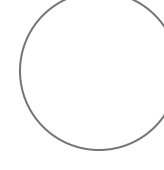
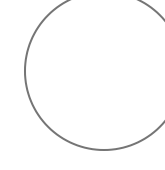
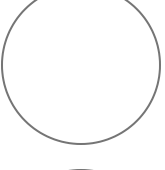
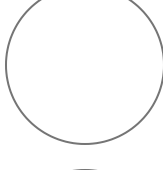
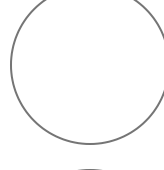
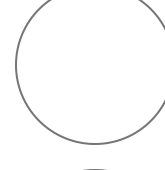
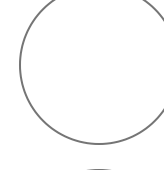
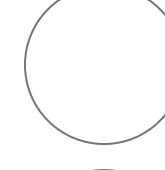
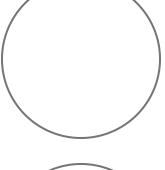
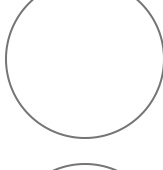
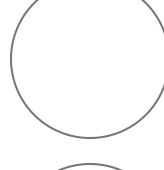
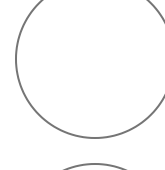
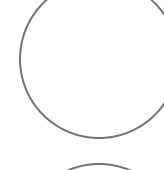
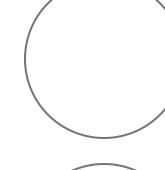
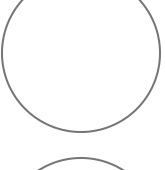
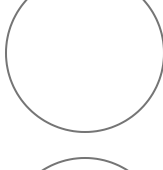
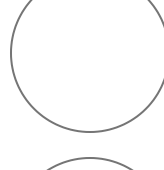
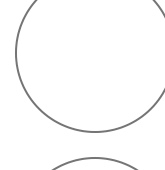
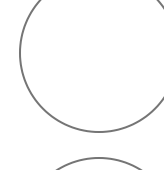
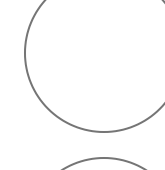
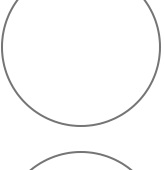
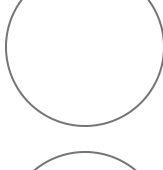
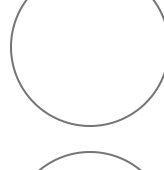
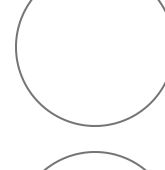
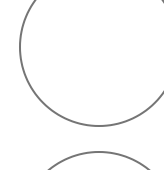

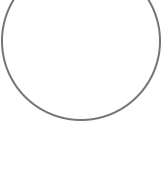
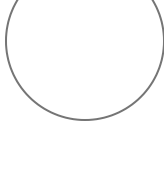
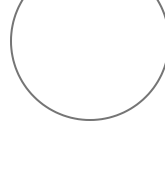
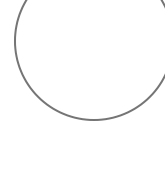
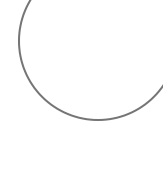
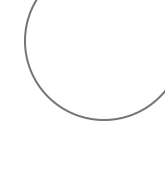
**16** remove and discard liquid from the row, being careful not to touch the beads with the pipette or to let the beads dry for more than 30 seconds.

- 16.1 Add  100  $\mu$ L of 80%EtOH to each well of beads. Incubate at  Room temperature for  00:00:30 30s
- 16.2 Remove the EtOH, then immediately add another  100  $\mu$ L of 80% EtOH to the wells, incubate for  00:00:30  Room temperature . 30s
- 16.3 Remove ALL EtOH, and let the row of beads dry just enough to lose some shine but not enough to start cracking. This should be approximately  00:00:30 to  00:01:00 . 1m 30s
- 16.4 Remove the plate with cleaned beads from the magnetic plate, and add  30  $\mu$ L of  Buffer EB **Qiagen Catalog #19086** to each well of beads, pipette mixing each well thoroughly. Incubate  00:05:00 at  Room temperature 5m
- 16.5 Place back on the magnetic rack for  00:01:00 until liquid is clear again. 1m
- 16.6 Roll back the foil on the final cleaned reactions plate for the appropriate row. Remove the  30  $\mu$ L clear eluate from the bead-cleanup plate, and place in the appropriate wells of the final cleaned reactions plate. Immediately cover this cleaned PCR product with either 8-strip caps.
- 16.7 uncover the next row of samples for cleaning and  [go to step #14](#) until all rows are cleaned.

## Prepare Indexing PCR

- 17 Create an indexing plate map and make sure your chosen indexes (iNext indexes) are color balanced if you aren't doing full 96-well plates at one time.

|   | 1                         | 2                          | 3    | 4                           | 5                            | 6 |
|---|---------------------------|----------------------------|------|-----------------------------|------------------------------|---|
| A | S8 iNextA-F +<br>iNext1-R | S16 iNextA-F +<br>iNext2-R | etc. | MiFishPCRneg i<br>NextAF+4R | Crust..PCRneg i<br>NextAF+5R |   |
| B | S7 iNextB-F +<br>iNext1-R | S15 iNextB-F +<br>iNext2-R |      |                             |                              |   |
| C | S6 iNextC-F +<br>iNext1-R | S14 iNextC-F +<br>iNext2-R |      |                             |                              |   |
| D |                           | etc.                       |      |                             |                              |   |
| E |                           |                            |      |                             |                              |   |
| F |                           |                            |      |                             |                              |   |
| G |                           |                            |      |                             |                              |   |
| H |                           |                            |      |                             |                              |   |

|   | 7   | 8   | 9   | 10   | 11  | 12  |
|---|---|---|---|--|---|---|
| A |    |    |    |    |    |    |
| B |    |    |    |    |    |    |
| C |    |    |    |    |    |    |
| D |    |    |    |    |    |    |
| E |   |   |   |   |   |   |
| F |  |  |  |  |  |  |
| G |  |  |  |  |  |  |
| H |  |  |  |  |  |  |

See: Glenn, Travis C., Roger A. Nilsen, Troy J. Kieran, Jon G. Sanders, Natalia J. Bayona-Vásquez, John W. Finger, Todd W. Pierson, et al. 2019. "Adapterama I: Universal Stubs and Primers for 384 Unique Dual-Indexed or 147,456 Combinatorially-Indexed Illumina Libraries (iTru & iNext)." *PeerJ* 2019 (10). <https://doi.org/10.7717/peerj.7755>. Supplemental file S10

Prepare working solutions of [M] 5 micromolar ( $\mu\text{M}$ ) of each indexing primer you intend to use.

## 18 Indexing PCR Mastermix Recipe:

🧪 6  $\mu\text{L}$  🧬 2x Kapa HiFi Hotstart Readymix **Kapa Biosystems Catalog #KK2602**

🧪 2.1  $\mu\text{L}$  🧬 Nuclease-free water **Integrated DNA Technologies, Inc. (IDT) Catalog #11-05-01-14**  
per sample.

Multiply by number of wells \*10% as explained above, to create master mix.

## 19 In a new, clean 96-well plate (UV before use if possible and prepare in a pre-PCR space):

Add 🧪 8.1  $\mu\text{L}$  Indexing Mastermix to each well that will be used and add 🧪 0.7  $\mu\text{L}$  of the [M] 5 micromolar ( $\mu\text{M}$ ) iNext forward indexed primer for each horizontal row of the plate (8 letters), and 🧪 0.7  $\mu\text{L}$  [M] 5 micromolar ( $\mu\text{M}$ ) of the iNext reverse indexed primer for each vertical column of the plate (12 numbers) according to the indexing plate map.

Take the prepared indexing reactions to the post-PCR space to add the cleaned PCR product.

## 20 In the post-PCR area, add 2.5uL of cleaned PCR 1 product to their associated wells from the indexing plate map.

### Indexing PCR Conditions

21 🌡️ 95 °C ⌚ 00:03:00

4m 35s

8 cycles of:

🌡️ 98 °C ⌚ 00:00:20

🌡️ 65 °C ⌚ 00:00:15

🌡️ 72 °C ⌚ 00:01:00

then hold 🌡️ 4 °C

### Optional gel to check Indexing PCR

22 Optional: visualize PCR products in a 1.7-2% gel. Bands should be around 350-400bp.

### Combine and Clean all indexed samples from each plate

23 Combine 10uL of up to 70 indexed samples (library) into a single 1.5mL tube. If there are more than 70 samples, you will need another tube.

24 Multiply the volume of the pooled libraries in each tube by 0.9 to get the volume of Ampure XP beads needed to clean up the reactions.

### Perform a 0.9x bead cleanup with Ampure XP beads

28m

25 In the 1.5mL tube of pooled libraries, add 0.9x volume of Ampure XP beads and pipette mix well. incubate 10m

🌡 Room temperature for ⌚ 00:10:00

26 Make enough fresh 80% EtOH to have 2x the total volume of the beads+library pool plus a bit extra.

27 Place 1.5mL tube into a magnetic rack

5m



## Equipment

### Magnetic Rack for for 1.5 mL Tubes

NAME

Magnetic Rack for DNA, RNA Purification; for 1.5 mL centrifuge Tubes

TYPE

Sergi Lab Supplies



BRAND









B0BZWXZM22

SKU

[https://www.amazon.com/Magnetic-Rack-Purification-centrifuge-Tubes/dp/B0BZWXZM22/ref=asc\\_df\\_B0BZWXZM22/?tag=hyprod-20&linkCode=df0&hvadid=652498086131&hvpos=&hvnetw=g&hvrnd=6716034042841103246&hvpone=&hvpstwo=&hvgmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=9](https://www.amazon.com/Magnetic-Rack-Purification-centrifuge-Tubes/dp/B0BZWXZM22/ref=asc_df_B0BZWXZM22/?tag=hyprod-20&linkCode=df0&hvadid=652498086131&hvpos=&hvnetw=g&hvrnd=6716034042841103246&hvpone=&hvpstwo=&hvgmt=&hvdev=c&hvdvcmdl=&hvlocint=&hvlocphy=9)

LINK

and incubate  Room temperature for  00:05:00

- 28 Discard liquid and add an equal or greater volume of 80% EtOH. Incubate  Room temperature for  00:01:00 1m
- 29 Repeat the ethanol wash a second time  [go to step #28](#) , then after the second 80% EtOH wash, remove all EtOH and dry the beads slightly (just until no longer wet-looking but not cracking either).
- 30 Resuspend beads with  100  $\mu$ L  Buffer EB **Qiagen Catalog #19086** by pipette mixing thoroughly, 10m  
Incubate  Room temperature  00:10:00
- 31 Place 1.5 mL tube back on magnet rack and wait until liquid is clear, approximately  00:02:00 2m

- 32**      remove 100uL of the clear eluate from the tube with beads while on the magnet and place in a new 1.5mL tube.
- 33**      Quantify with Qubit Broad range and visualize in a gel, then send for sequencing on a lane of MiSeq.