



OCT 18, 2023


Optical densitometry of neuronal fibers

Núria
Peñuelas¹

¹Vall d'Hebron Research Institute

Vilalab Public

Vilalab Public

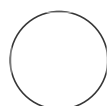

Núria Peñuelas
Vall d'Hebron Research Institute

OPEN  ACCESS



ABSTRACT

Optical densitometry of striatal neuronal fibers in the rodent brain

DOI:
dx.doi.org/10.17504/protocols.io.8epv5x3nng1b/v1

Protocol Citation: Núria Peñuelas 2023. Optical densitometry of neuronal fibers. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.8epv5x3nng1b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Oct 18, 2023

Last Modified: Oct 18, 2023

- 1 Scan the immunostained serial coronal sections covering different caudorostral levels of the brain regions (for striatum, 4 sections/animal) with an Epson Perfection v750 Pro scanner.
- 2 Open the resulting images in ImageJ (RRID:SCR_003070,<https://imagej.net/>).
- 3 With the free hand selections tool, draw the region of interest (e.g. striatum, nucleus accumbens and olfactory tubercle) and measure the optical density (pixel brightness). Then draw a region of interest in a blank area (e.g. cortex) and also measure the optical density.
- 4 Calculate the optical density or absorbance using the formula: $-\log_{10} (\text{Striatum Intensity} / \text{Cortex Intensity})$.