

JAN 12, 2024

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.ewov1q15pgr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1q15pgr2/v1)

**Protocol Citation:** Jacob Robertson, Erin Garza 2024. Preservation method for long-term storage of fluorescently labeled cells for microscopy. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.ewov1q15pgr2/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working


**Created:** Dec 18, 2023

## 🌐 Preservation method for long-term storage of fluorescently labeled cells for microscopy

Jacob Robertson<sup>1</sup>, Erin Garza<sup>2</sup>

<sup>1</sup>University of California San Diego; <sup>2</sup>J. Craig Venter Institute

JCVI West Protocols

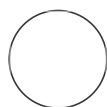

Erin Garza  
 J. Craig Venter Institute

### ABSTRACT

This method can be used to preserve fluorescently labeled bacterial cells for long-term storage before imaging on a fluorescent microscope. Using this protocol, samples can be saved at 4°C for weeks-months while maintaining strong, easily detectable fluorescence and cell integrity.

### IMAGE ATTRIBUTION

The image was produced by the authors using a Leica SP8 confocal microscope. This is a comparative photo of biofilms expressing a fluorescent protein that were preserved and imaged after 1 week vs 5 months.

### GUIDELINES

This protocol has been successfully used on liquid bacterial cultures and on biofilms on solid plastic material. Cells were imaged using a confocal microscope up to 5 months after preservation. Fluorescence may be maintained longer than 5 months, but has not been tested by the authors.

### MATERIALS

Microfuge tubes  
 Glycerol  
 Paraformaldehyde  
 PBS  
 4°C storage

### SAFETY WARNINGS



Paraformaldehyde is toxic and a skin irritant. Wear appropriate PPE when preparing and working with this solution.

**Last Modified:** Jan 12, 2024

## BEFORE START INSTRUCTIONS

**PROTOCOL integer ID:**  
92469

Prepare the preservation components: 10% glycerol (sterile), 4% paraformaldehyde in PBS pH 7.4 (sterile)

**Keywords:** Preserve,  
Microscopy, Fluorescence,  
Long-term storage, Confocal

### Sample



2m

- 1 Collect your cells (up to 500  $\mu$ l) in a sterile microfuge tube.  
- This protocol can also be used on solid material containing biofilms.

2m

### Preserve

3m

- 2 To your sample tube, add  250  $\mu$ L sterile 10% glycerol and  250  $\mu$ L sterile 4% paraformaldehyde in PBS (pH 7.4).  
- If preserving solid material, make sure the entire sample is submerged in the preservation solution.

2m

- 3 Gently Mix

30s

### Storage

1m

- 4 Store samples at  4  $^{\circ}$ C until ready to image.

1m