



VERSION 2

JUN 28, 2023

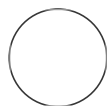
OPEN  ACCESS**DOI:**dx.doi.org/10.17504/protocols.io.261ged6xyv47/v2**Protocol Citation:** Alex N Nguyen Ba 2023. Converting ssDNA oligos to dsDNA with T4 DNA polymerase.**protocols.io**<https://dx.doi.org/10.17504/protocols.io.261ged6xyv47/v2>Version created by [Alex N Nguyen Ba](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jun 28, 2023**Last Modified:** Jun 28, 2023**PROTOCOL integer ID:**
84169

Converting ssDNA oligos to dsDNA with T4 DNA polymerase V.2

Alex N Nguyen Ba¹¹University of Toronto

Alex N Nguyen Ba

ABSTRACT

This protocol allows one to convert ssDNA to dsDNA oligos. In principle, one can buy two complementary oligos and anneal them. However, there are a few cases where randomized bases are desired, and therefore complementary oligos can not be ordered. The protocol works similarly to PCR, except that the reaction is not performed at the high temperature of typical thermocycling, because the oligos are short and their melting temperature can be in the order of the extension temperature of high-fidelity polymerases. This protocol is therefore particularly useful when creating barcodes or when short dsDNA fragments containing degenerate bases are desired.

Primer design

- 1 Design your target single-stranded DNA oligo such that the 3' end contains the reverse complement of an extending oligonucleotide. For example, appending "GTCATAGCTGTTTCCTG" to the end of your oligo will allow an oligonucleotide matching the M13 Reverse (-27) primer to extend it (5'-CAGGAAACAGCTATGAC-3').
- 2 Resuspend your oligonucleotides in 1/10th TE (10mM Tris, 0.1mM EDTA pH 8) to a final concentration of **100 micromolar (μM)** to form your stock oligonucleotide solutions.
- 3 Make working oligonucleotide solutions by diluting the stock oligonucleotides to a final concentration of **10 micromolar (μM)** in water.

Oligo annealing

- 4 Mix in a PCR tube **5 μL** of a **10 micromolar (μM)** stock of the target oligo with **10 μL** of a **10 micromolar (μM)** stock of the extending oligo. Add 27ul of molecular biology grade water.
- 5 Add **5 μL** of **NEBuffer 2.1 - 5.0 ml New England Biolabs Catalog #B7202S**. Mix well by vortexing.
- 6 Place the tube in a thermocycler and run the following protocol: 15m
 1. 95 degrees for 10 seconds
 2. Decrease by 1 degrees
 3. Repeat step 2 every 10 seconds, for 90 cycles.
 4. Hold at 4 degrees.



Oligo extension

- 7 Add **2.5 μL** **dNTPs Contributed by users** (**10 millimolar (mM)** stock) and mix well by vortexing.

8 Go to the thermocycler and prepare the reaction cycle. When ready, start the protocol and pause when the block reaches 0 degrees (when Step 1 begins).

40m

1. 5 minutes at 0 degrees C.
2. 5 minutes at 22 degrees C.
3. 30 minutes at 37 degrees C.
4. Hold at 0 degrees.


9 Add  0.5 µL of  T4 DNA Polymerase - 750 units New England Biolabs Catalog #M0203L and mix well.

10 Place the tube in the thermocycler and unpause the run.

Oligonucleotide cleanup

10m


11 As soon as the reaction is completed, add

 EDTA (0.5 M), pH 8.0 Life Technologies Catalog #AM9260G

to [M] 10 millimolar (mM) .

12 Use a PCR purification kit to purify the oligo. Ensure the kit is able to purify short oligonucleotides if needed (a kit such as

10m

 Monarch PCR and DNA Cleanup Kit - 250 preps New England Biolabs Catalog #T1030L

can

be suitable).