



VERSION 2
NOV 29, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.5qpvp6bv4o/v2

Protocol Citation: Rebecca Andrews, Joanne Lachica, Steven F. Lee, Sonia Ghandi 2023. Single-molecule Immunofluorescence Tissue Staining Protocol for Oligomer Imaging.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.5qpvp6bv4o/v2>
Version created by [Rebecca Andrews](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Single-molecule Immunofluorescence Tissue Staining Protocol for Oligomer Imaging V.2

Rebecca Andrews¹, Joanne Lachica², Steven F. Lee¹, Sonia Ghandi²

¹University of Cambridge; ²UCL

ASAP Collaborative Research Network



Rebecca Andrews
University of Cambridge

ABSTRACT

This protocol details background fluorescence quenching and immunofluorescence staining of human brain tissue for oligomer imaging.

ATTACHMENTS

[kb3ib25np.pdf](#)

GUIDELINES

- Use only clean bottles, flasks, magnetic stirrers, tweezers, weighing spatulas, measuring cylinders – everything should be cleaned, dried and covered if left on the side before next use.
- Everything should be handled with clean tweezers – gloves should not touch the samples, solutions and ideally anything placed into the solutions where the slides are.

Protocol status: Working
We use this protocol and it's working

Created: Jan 19, 2023

Last Modified: Nov 29, 2023

PROTOCOL integer ID:
75550


Keywords:
immunofluorescences,
oligomer imaging

MATERIALS

Materials and Reagents



- Microtome
- Glass slides
- Xylene solution
- 100% alcohol
- Methanol
- Hydrogen peroxide (H₂O₂) solution
- Citrate buffer pH6
- Milli Q water
- Pressure cooker
- PBS
- Goat Serum 10%
- AlexaFluor antibody
- 0.1% Sudan black solution
- Vectashield
- Overslip

Immunofluorescences staining protocol for oligomer imaging 10m

1 Cut  8 µm tissue sections on a microtome and load onto glass slides.


2 Dry slides  Overnight at  37 °C – cover over the top. 10m




3 Before staining commences keep slides for a few hours but ideally  Overnight at  60 °C . 10m




4 De-wax sections through three pots of xylene solution. Use each fresh pots of xylene each time.

4.1 De-wax section through pot of xylene solution for  00:02:00 . (1/3)

2m

4.2 De-wax section through pot of xylene solution for  00:02:00 . (/23)

2m


4.3 De-wax section through pot of xylene solution for  00:02:00 . (3/3)

2m


5 Take sections through two pots of 100% alcohol.

Note



Use fresh pots each time – methylated spirits.

5.1 Take sections through pot of 100% alcohol for  00:02:00 . (1/2)

2m

5.2 Take sections through pot of 100% alcohol for  00:02:00 . (2/2)

2m

6 Put slides into methanol + hydrogen peroxide (H₂O₂) solution (100 ml: 1 ml) for  00:10:00 at  Room temperature **fresh pot each time.**

10m

Note



This process will block any staining of endogenous peroxidase in the tissue sections.

7 Perform necessary antigen retrieval pretreatments by pressure cooking in citrate buffer.

7.1 Pressure cook sections in citrate buffer  6 for  00:10:00 at pressure (wait for it to release high pressure air) in **cleaned pressure cooker**. 

7.2 Cool under running **Milli Q** water – *never* under tap water.



8 Block non-specific antigen/antibody binding by placing sections in **PBS and Goat Serum 10%** (taken  the **DNA PAINT protocol**) for  00:30:00 .

9 Apply primary antibody for  01:00:00 at  Room temperature . 





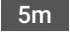

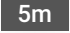

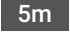
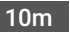



10 Wash in PBS with fresh buffer (**at least filtered if not cell culture grade**).



10.1 Wash  00:05:00 in PBS clean squirty bottle with fresh buffer. (1/3) 

10.2 Wash  00:05:00 in PBS clean squirty bottle with fresh buffer. (2/3) 

10.3 Wash  00:05:00 in PBS clean squirty bottle with fresh buffer. (3/3) 

- 11 Apply secondary AlexaFluor antibody for  01:00:00 at  Room temperature in the dark. 
- 12 Wash in PBS.

- 12.1 Wash  00:05:00 in PBS in the dark. (1/3) 
- 12.2 Wash  00:05:00 in PBS in the dark. (2/3) 
- 12.3 Wash  00:05:00 in PBS in the dark. (3/3) 
- 13 Add **filtered (0.22 μ m)** 0.1% Sudan black solution (0.1% sudan black/70% ethanol) for  00:10:00 
  Room temperature in the dark.
- 14 Wash 2-3 times in 30% ethanol.

- 15 Mount section with Vectashield and coverslip (**Plasma cleaned slides**).

Take for imaging

