



Jan 18, 2022

# Protocols from: Evolutionary analyses of visual opsin genes in frogs and toads: diversity, duplication, and positive selection

Ryan K Schott<sup>1</sup>, Leah Perez<sup>2</sup>, Matthew A Kwiatkowski<sup>2</sup>, Vance Imhoff<sup>3</sup>, Jennifer M Gumm<sup>3</sup>

<sup>1</sup>York University; <sup>2</sup>Stephen F. Austin State University; <sup>3</sup>US Fish and Wildlife Service

1



[dx.doi.org/10.17504/protocols.io.b3yqqpvw](https://dx.doi.org/10.17504/protocols.io.b3yqqpvw)

Schott Lab

Ryan K Schott

Protocols used to extract mRNA from frog retinas, create cDNA libraries, and amplify opsins by PCR for sequencing at the UT core facility under their standard protocols. These protocols were used to obtain the opsin sequences in the paper: Evolutionary analyses of visual opsin genes in frogs and toads: diversity, duplication, and positive selection.

DOI

[dx.doi.org/10.17504/protocols.io.b3yqqpvw](https://dx.doi.org/10.17504/protocols.io.b3yqqpvw)

<https://doi.org/10.1002/ece3.8595>

Ryan K Schott, Leah Perez, Matthew A Kwiatkowski, Vance Imhoff, Jennifer M Gumm 2022. Protocols from: Evolutionary analyses of visual opsin genes in frogs and toads: diversity, duplication, and positive selection. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.b3yqqpvw>



protocol

Schott RK, L Perez, MA Kwiatkowski, V Imhoff, JM Gumm. 2022. Evolutionary analyses of visual opsin genes in frogs and toads: diversity, duplication, and positive selection. Ecology and Evolution (in press).

protocol ,







Jan 18, 2022












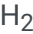

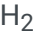
Jan 27, 2022

57072



Frog retinas, RNeasy Kit, Qias shredder.

## RNA Extraction

- 1 Transfer sample into a  **1.5 mL** collection tube.
- 2 Pipette off RNALATER.
- 3 Add  **600 µL** Buffer RLT.
- 4 Add  **6 µL** Beta-mercaptoethanol.
- 5 Disrupt tissue with sterile pestle.
- 6 Pipette into Qias shredder column. Spin  **00:02:00** @ 8,000rpm. 2m
- 7 Remove Qias shredder column; Add cap; Spin  **00:03:00** @ max speed. 3m
- 8 Add  **600 µL** 70% Ethanol to new collection tube.
- 9 Transfer lysate to the collection tube; mix lysate and 70% Ethanol by pipetting.

- 10 Transfer lysate to RNeasy column (  **700 µL** at a time). Spin  **00:00:15** @ 9,800rpm; Discard flow through. Add rest of lysate; Spin  **00:00:15** @ 9,800rpm; Discard flow through. 30s
- 11 Add  **700 µL** Buffer RWI. Spin  **00:00:15** @ 9,800rpm. 15s
- 12 Transfer RNeasy column to new collection tube.
- 13 Add  **500 µL** Buffer RPE. Spin Spin  **00:00:15** @ 9,800rpm; Discard flow through. 15s
- 14 Add  **500 µL** Buffer RPE – Spin  **00:01:00** @ 9,800rpm; Discard Flow through/ Spin  **00:02:00** @ 13,000rpm. 3m
- 15 Transfer RNeasy column to new collection tube.
- 16 Elute with  **30 µL** RNase-Free H<sub>2</sub>O. Spin  **00:01:00** @ 13,000rpm. 1m
- 17 Elute with  **30 µL** RNase-Free H<sub>2</sub>O- Spin  **00:01:00** @ 13,000rpm. 1m

#### CDNA Synthesis

- 18 Combine mRNA and RNase-free H<sub>2</sub>O to standardize all samples to aliquots containing  **0.4 µg** mRNA total in  **10 µL** .
- 19 Make 2 Master mixes:

19.1 Master Mix 1: add **1 µL** dNTP mix and **2 µL** dT primer per sample.

19.2 Master Mix 2: add **4 µL** Buffer, **2 µL** DTT and **0.5 µL** RNAase inhibitor per sample.

20 Pipette **3 µL** of Master Mix 1 into each sample.

21 Place sample on dry bath at **65 °C** for **00:05:00** . 5m

22 **6.5 µL** Put samples on ice for **00:01:00** . 1m

23 Pipette **6.5 µL** of Master Mix 2 into each sample.

24 Pipette **1 µL** Superscript into each sample.

25 **65 °C** Incubate samples at room temp for **00:10:00** . 10m

26 Incubate samples at **42 °C** for **00:50:00** . 50m

PCR 17m 30s

27 Keep all reagents on ice at all times.

28 Make a master mix. Per sample add the following:

- **2.0 µL** 10X Buffer
- **1.0 µL** [**50 millimolar (mM)**] MgSO<sub>4</sub>
- **0.5 µL** dNTP mix ([**10 micromolar (µM)**] each)
- **18.4 µL** ddH<sub>2</sub>O
- **1 µL** forward primer ([**10 micromolar (µM)**])
- **1 µL** reverse primer ([**10 micromolar (µM)**])
- **0.5 µL** Taq polymerase

29 Mix well by spinning.

30 Add **24 µL** of Master Mix to each PCR tube.

31 Add **1 µL** of sample for a total of **25 µL** per tube.

32 Program the thermocycler for the following program:

17m 30s

- **95 °C** for **00:10:00**
- **94 °C** for **00:02:00**
- REPEAT FOLLOWING 3 steps 35-50 times:
  - **94 °C** for **00:00:30**
  - **45-50 °C** for **00:01:00** \*temperature depends on primer
  - **72 °C** for **00:02:00**
- **72 °C** for **00:02:00**
- **4 °C** hold