



Version 2

Sep 07, 2020

covid 19 indirect detection thru rise of 100nm filter air resistance V.2

Ladislav Nevery

In Development

dx.doi.org/10.17504/protocols.io.bkznkx5e

Coronavirus Method Development Community

XPRIZE Rapid Covid Testing

1 more workspace



Ladislav Nevery

ABSTRACT

By passing saliva sample thru 200nm filter to remove particles larger than 120nm covid 19 virus particles and concentrating resulting fluid to very small 1 mm² area of 100nm filter. Pores in this small area of 100nm filter should be clogged up fast raising filter air resistance. Due to high virus density 5.2 log₁₀ ml saliva of infected person should clog filter much faster than healthy person. Test is setup in such way that act of compressor automatically turning off after reaching threshold PSI value will indicate high count of particles larger than 100nm of being present in measured filter area.

ATTACHMENTS

[PrototypeV8Printable.zip](#)

DOI

dx.doi.org/10.17504/protocols.io.bkznkx5e

PROTOCOL CITATION

Ladislav Nevery 2020. covid 19 indirect detection thru rise of 100nm filter air resistance . **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bkznkx5e>
 Version created by Ladislav Nevery

KEYWORDS

covid, pressure, test

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 07, 2020

LAST MODIFIED

Sep 07, 2020

PROTOCOL INTEGER ID

41742

MATERIALS

NAME	CATALOG #	VENDOR
pragopor 10		pragochema
pragopor 8		pragochema
3d filament		
rubber seal 24x20x2		

STEPS MATERIALS

NAME	CATALOG #	VENDOR
pragopor 10		pragochema
rubber seal 24x20x2		
Target2™ Nylon Syringe Filter, 0.2µm, 30mm	F2500-2	Thermo Fisher

EQUIPMENT


NAME	CATALOG #	VENDOR
Xiaomi Portable Compressor	00000000	

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.


1

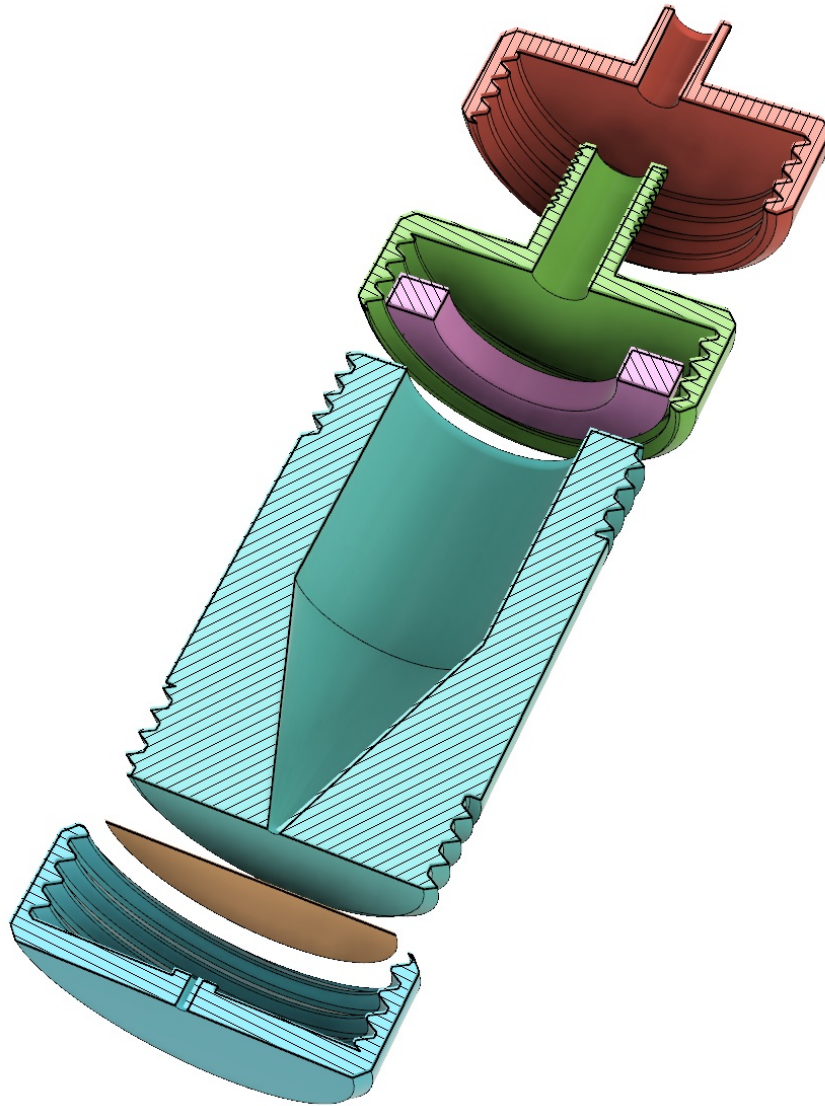


3d print all required parts.

2 insert any 100nm 24mm filter like for example


5m

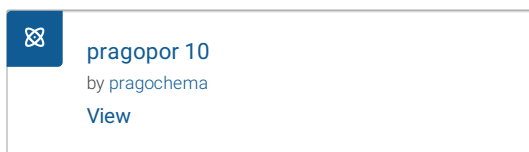

between blueparts and close with red cup

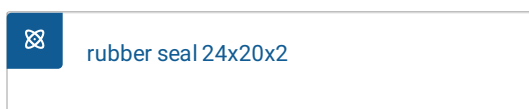


.Place any comon commercial 200nm filter between syringe and red syringe input port. colors are on attached image above



Target2™ Nylon Syringe Filter, 0.2µm, 30mm
by Thermo Fisher
Catalog #: F2500-2


- 3 pass sample thru both 200nm and 100nm filter by applying manual pressure on syringe. 5m

- 4 now that 1 mm² area of 100nm filter has enough pores clogged up with median 120nm sized virus particles. replace ^{5m} red syringe cap with green compressor input cap. Use any compressor which can turn itself off when reaching set pressure. I used



Xiaomi Portable Compressor
Portable Compressor
Xiaomi 00000000 [Link](#)

set compressor to turn off after reaching threshold amount of psi and turn the compressor on.



On sample without covid 19 virus threshold psi should not be reached

. If no compressor with proper psi sensitivity is available then monitor pressure in concentration chamber by any external sensitive manometer