



Upload image


Jul 31, 2020

SOLUTION- 11 - CD4+T cells isolation buffer

Marco Cosentino¹, Elisa Storelli¹, Alessandra Luini¹, Massimiliano LM Legnaro¹, Emanuela Rasini¹, Marco Ferrari¹, Franca Marino¹

¹Center for Research in Medical Pharmacology, University of Insubria (Varese, Italy)

1 Works for me dx.doi.org/10.17504/protocols.io.bi79khr6

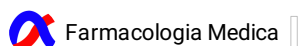

DOI

dx.doi.org/10.17504/protocols.io.bi79khr6

DOCUMENT CITATION

Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari, Franca Marino 2020. SOLUTION- 11 - CD4+T cells isolation buffer. **protocols.io**
dx.doi.org/10.17504/protocols.io.bi79khr6

LICENSE

This is an open access document distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 31, 2020

LAST MODIFIED

Jul 31, 2020

DOCUMENT INTEGER ID

39905

Composition g/ml:

BSA 0.250 g

EDTA 0.2235 g

Add PBS(Ca²⁺ and Mg²⁺ free) to 240 ml, adjust to pH7.4 with NaOH or HCl, bring to 250 ml with PBS.

The substances are in the closet 1 and in refrigerator 1 (Room TS08), weigh the substances with Gibertini balance (Room TS08).

Filter this buffer by using a sterile syringe and filter.

Sterile filters and syringes are in closet 1 (Room PS03).

Storage: 4 °C refrigerator 2 - (Room PS03).

BSA: catalog number A2153 , Sigma, Italy

EDTA: catalog number ED2SS, Sigma, Italy

Prymo siringe filter: catalog number EPSPE2230 Euroclone, Italy.