



FEB 27, 2023

OPEN ACCESS

DOI:
[dx.doi.org/10.17504/protocols.io.bupvvnv6](https://doi.org/10.17504/protocols.io.bupvvnv6)

External link:
<https://doi.org/https://doi.org/10.1101/2021.06.15.447579>

Protocol Citation: Tetsuya Komatsu, Kenji Ohya, Justice Opare Odoi, Shota Suganuma, Kotaro Sawai, Takayuki Wada, Tomotada Iwamoto, Fumito Maruyama 2023. Genomic DNA extraction from *Mycobacterium avium*. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.bupvvnv6>

MANUSCRIPT CITATION: Komatsu T, Ohya K, Ota A, Nishiuchi Y, Yano H, Matsuo K, Odoi JO, Suganuma S, Sawai K, Hasebe A, Asai T, Yanai T, Fukushima H, Wada T, Yoshida S, Ito T, Arikawa K, Kawai M, Ato M, Baughn AD, Iwamoto T, Maruyama F, Genomic features of subsp. isolated from pigs in Japan. GigaByte doi: <https://doi.org/10.1101/2021.06.15.447579>

🌐 Genomic DNA extraction from *Mycobacterium avium*

📁 In 1 collection

Tetsuya Komatsu¹, Kenji Ohya², Justice Opare Odoi³, Shota Suganuma², Kotaro Sawai², Takayuki Wada⁴, Tomotada Iwamoto⁵, Fumito Maruyama⁶

¹Aichi Prefectural Chuo Livestock Hygiene Service Center, Okazaki, Aichi, Japan;


²Faculty of Applied Biological Sciences, Gifu University, Gifu, Gifu, Japan;

³United Graduate School of Veterinary Sciences, Gifu University, Gifu, Gifu, Japan;

⁴Graduate School of Human Life Science, Osaka City University, Osaka, Osaka, Japan;

⁵Department of Infectious Diseases, Kobe Institute of Health, Kobe, Hyogo, Japan;

⁶Office of Academic Research and Industry-Government Collaboration, Hiroshima University, Higashi-Hiroshima, Hiroshima, Japan

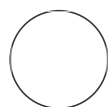

Tetsuya Komatsu

Aichi Prefectural Chuo Livestock Hygiene Service Center, Oka...

ABSTRACT

Mycobacterium avium subsp. *hominissuis* (MAH) is one of the most important agents causing non-tuberculosis mycobacterial infection in human and pigs. Genome analysis on MAH of human isolates has been proceeding, however, those of pigs are limited despite its potential source of infection to human. Here, we isolated MAH from pig lymph nodes or livers and obtained their genomes.

MATERIALS

TE: Tris-EDTA Buffer; 10 mM Tris-Cl (pH 8.0), 1 mM EDTA; sterilize by autoclaving. Phenol/chloroform/isoamyl alcohol; 25:24:1

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: May 01, 2021

Last Modified: Feb 27, 2023

PROTOCOL integer ID:
49621

Keywords: Mycobacterium, lymph node, pig, genome extraction

- 1 Suspend colonies on 7H11 agar into 300 µl of TE.
- 2 Centrifuge at 5,000 *g* for 10 min.
- 3 Discard supernatant, followed by suspension in 300 µl of acetone.
- 4 Centrifuge at 5,000 *g* for 10 min.
- 5 Discard acetone, followed by suspension in 200 µl of TE.

- 6** Add lysozyme at the final concentration of 100 µg/ml.
- 7** Incubate at 37°C overnight.
- 8** Add proteinase K at the final concentration of 0.2 mg/ml.
- 9** Incubate at 50°C for 3 hrs.
- 10** Add 1/10 volume of 3 M NaOAc.
- 11** Add equal volume of phenol/chloroform/isoamyl-alcohol and shake vigorously.
- 12** Centrifuge at 12,000 *g* for 10 min.
- 13** Transfer aqueous phase to a new tube.
- 14** Add equal volume of 2-propanol.

- 15 Centrifuge at 12,000 *g* for 10 min.
- 16 Remove supernatant and wash with 70% ethanol.
- 17 Suspended in 400 µl of TE, followed by adding RNase at the final concentration of 100 µg/ml.
- 18 Incubate at 37°C for 30 min.
- 19 Repeat steps 10 to 16.
- 20 Suspended in 100 µl of TE.