



MAY 22, 2023

OPEN ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.x54v9dojzg3e/v1

Collection Citation: Rolando DZ Lyles, M Julia Martinez, Benjamin Sherman, Kerry Burnstein 2023. Automation, live-cell imaging, and endpoint cell viability for 96-well plate drug screens.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.x54v9dojzg3e/v1>

License: This is an open access collection distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this collection and it's working

Created: Mar 01, 2023

Last Modified: May 22, 2023

COLLECTION integer ID:
77846

Automation, live-cell imaging, and endpoint cell viability for 96-well plate drug screens

Rolando DZ

Lyles^{1,2,3},

Kerry Burnstein^{1,4,3}

Benjamin


M Julia Martinez^{1,4}, Sherman^{1,4},

¹University of Miami Miller School of Medicine;

²Sheila and David Fuente Graduate Program in Cancer Biology;

³Sylvester Comprehensive Cancer Center;

⁴Department of Molecular & Cellular Pharmacology

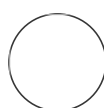

Rolando DZ Lyles

ABSTRACT

To streamline the identification of potentially active cancer therapeutics, here we describe a highly adaptable semi-automated approach to screen compounds simultaneously across a panel of cell lines. These protocols leverage automation to enhance robustness, reproducibility, and throughput while integrating the IncuCyte ZOOM live-cell imaging platform and the CellTiter-Glo endpoint viability assay to assess drug efficacy. The integration of both evaluative methods strategically bypasses the shortcomings of each approach individually allowing for more thorough and detailed analysis within a single drug screen. The expected output from protocol utilization includes traditional dose-response curves, IC50, area-under the curve (AUC), GR50, area-over the curve (AOC) as well as live-cell imaging to identify cell specific morphological changes, cytostatic, or cytotoxic effects of 72-hour drug treatments. This adaptable protocol can be employed across cancer model systems and represents a reproducible procedure to optimize 96-well plate cell growth conditions compatible with the integrated drug screen and simultaneously assess drug efficacy across multiple cell lines in future cancer research studies.

ATTACHMENTS

[S1.pdf](#)

Keywords: live-cell imaging, endpoint cell viability, 96-well plate drug screens

MATERIALS

For this protocol you will need:

CellTiter-glo (Promega, Cat# G7572) Luminescence Viability Assay or comparable endpoint cell viability assay

Luminometer compatible with the Promega CellTiter-Glo assay. This protocol features the Promega GloMax (Promega, Cat# GM3500) Explorer Multimode Microplate Reader.

A liquid handler system. This protocol features the Opentrons OT-2 Robot liquid handler system.

OT-2 Required attachments:

- P300 Single-Channel GEN2 Pipette
- P300 8-Channel GEN2 Pipette
- Opentrons 96 Tip Racks 300 µL
- 12-Channel Reservoirs for Automation (USA Scientific, Cat# 1061-8150)

White opaque welled, clear bottom 96-well tissue culture compatible plates. This protocol features Greiner Bio-one Cell culture microplate, 96 well, PS, F-Bottom, µClear (Sigma Aldrich, M0437-32EA)

Sartorius IncuCyte ZOOM Live-cell analysis system or any comparable live cell imaging platform.

- Compatible IncuCyte ZOOM Live cell analysis software.

Sterile tissue culture hood

Standard 8-channel p200 multi-channel pipette

White Opaque Lab Tape

Standard 25 mL reagent reservoirs (VWR, Cat# 89094-662)

Standard tissue culture reagents


Standard 96-well plate (clear)

Standard 24-well plate (clear)

ATTACHMENTS

[S1.pdf](#)

Protocol



NAME

96-well plate cell growth optimization for integrated live-cell and endpoint viability drug screening assay

VERSION 1


CREATED BY

Rolando DZ Lyles

OPEN

→

Protocol



NAME

96-well plate OT-2 liquid handler integrated live-cell and endpoint viability drug activity screen

VERSION 1

CREATED BY

Rolando DZ Lyles

OPEN

→