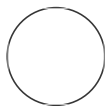




OCT 20, 2023

Proteolytic Peptide Desalting with C18 Hydrophilic-Lipophilic Balance (HLB) Cartridges

M A

J Bons¹, J P Rose¹, Watson¹, B Schilling¹¹Buck Institute for Research on Aging

M A Watson

ABSTRACT

Desalting of proteolytic peptide elution using C₁₈ Hydrophilic-Lipophilic Balance (HLB) Cartridges in preparation for downstream proteomic profiling.

OPEN  ACCESS**DOI:**

dx.doi.org/10.17504/protocols.io.3byl4q382vo5/v1

Protocol Citation: J Bons, J P Rose, M A Watson, B Schilling 2023. Proteolytic Peptide Desalting with C18 Hydrophilic-Lipophilic Balance (HLB) Cartridges.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.3byl4q382vo5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Oct 20, 2023

PROTOCOL integer ID:
89689

Keywords: Desalting, HLB, Peptides, Proteomics, Mass Spectrometry

- Oasis HLB 1cc cartridges, 10 mg sorbent (Waters, Cat. #WAT058882)
- Vacuum manifold
- Centrifuge
- Centrifugal vacuum concentrator
- 1.5-mL microcentrifuge tubes
- Condition buffer/Elution buffer: 50% acetonitrile (ACN), 0.2% formic acid (FA) in water
- Equilibration buffer/Wash buffer: 0.2% formic acid (FA) in water
- HPLC-grade water

- 1 Centrifuge the samples at 1,850 x *g* for 5 minutes at room temperature to pellet insoluble material.
- 2 Desalt the samples using Oasis HLB solid-phase extraction cartridges placed on top of a vacuum manifold as follows:
 - 2.1 Condition each cartridge two times with 800 µL of the condition buffer (50% ACN, 0.2% FA).
 - 2.2 Equilibrate each cartridge three times with 800 µL of the equilibration buffer (0.2% FA).

- 2.3 Load the peptide samples.
- 2.4 Wash each cartridge three times with 800 μL of the wash buffer (0.2% FA).
- 2.5 Elute peptides with 800 μL of the elution buffer (50% ACN, 0.2% FA) followed by a further 400 μL of the same elution buffer.
- 3 Vacuum dry the eluted peptide solution in a centrifugal vacuum concentrator.
- 4 Reconstitute the dried peptides in 0.2% FA for a final protein concentration of 1 $\mu\text{g}/\mu\text{L}$ and thoroughly mix the solution.
- 5 Centrifuge at 12,000 $\times g$ at room temperature for min, and store at -20°C .