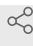




Jul 08, 2022

# Plasmid construction

 Forked from [Plasmid construction](#)Adam Yokom<sup>1</sup><sup>1</sup>University of California, Berkeley1 *Works for me* Share

This protocol is published without a DOI.



Adam Yokom

## DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](https://protocols.io) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](https://protocols.io), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

## ABSTRACT

plasmid construction

## PROTOCOL CITATION

Adam Yokom 2022. Plasmid construction. **protocols.io**  
<https://protocols.io/view/plasmid-construction-cczbsx2n>



## FORK NOTE

## FORK FROM

Forked from [Plasmid construction](#), Chunmei Chang

## KEYWORDS

ASAPCRN

#### LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

#### CREATED

Jul 08, 2022

#### LAST MODIFIED

Jul 08, 2022

#### PROTOCOL INTEGER ID

66307

#### MATERIALS TEXT

Q5® High-Fidelity DNA Polymerases (NEB)  
restriction enzymes (NEB)  
quick CIP (NEB)  
T4 ligase (NEB)  
Gel extraction kit (Bio Basic)

#### DISCLAIMER:

##### DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](https://protocols.io) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](https://protocols.io), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

- 1 Amplify the insert gene fragment by PCR with primers including 21 nt of overlapping sequence with the target gene.
- 2 Linearize the backbone with restriction enzymes (NEB).

- 3 Treat linearized backbone with quick CIP (NEB).
- 4 Run PCR products and linearized backbone in an agarose gel to confirm the size.
- 5 Purify the DNA from gel using a Gel extraction kit (Bio Basic).
- 6 Ligate the linearized backbone and the insert with the T4 ligase (NEB).
- 7 Transform the ligation products into home made competent cells.
- 8 Perform colony PCR to screen for colons that with inserted gene.
- 9 Sequence to verify that the inserted gene is correct.