



JUL 31, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.261gedpojv47/v1

Protocol Citation: Elisabeth Holzer 2023. EXPRESSION AND PURIFICATION OF HUMAN NEMO (GST-GFP-NEMO). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.261gedpojv47/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Jun 07, 2023


🌐 EXPRESSION AND PURIFICATION OF HUMAN NEMO (GST-GFP-NEMO)

Elisabeth Holzer^{1,2}

¹Max Perutz Labs, University of Vienna, Vienna Biocenter, Vienna, Austria ;

²Vienna Biocenter PhD Program, Doctoral School of the University of Vienna and Medical University of Vienna, Vienna, Austria

Liv

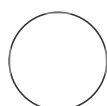

OLIVIA HARDING

ABSTRACT

This protocol describes how to express and purify human NEMO (IKK- γ) tagged N-terminally with GST and EGFP. The expression is performed with *E. coli* Rosetta pLysS cells. The protein is purified via a GST batch purification and gel filtration (SEC).

MATERIALS

pGEX-GST-Thrombin-TEV-GFP-NEMO (Addgene ID: 199781)

E. coli Rosetta (DE3) pLysS cells

LB medium with antibiotics: 50 μ l/ml ampicillin and 34 μ l/ml chloramphenicol

IPTG (Isopropyl- β -D-thiogalactopyranosid)

Thrombin (Serva, #36402.01)

Columns/Resin:

Glutathione Sepharose 4B (Cytiva, #17075605)

Superose 6 increase 10/300 column (Cytiva, #17517201)

Last Modified: Jul 31, 2023

BEFORE START INSTRUCTIONS

PROTOCOL integer ID:
83006

Keywords: ASAPCRN

Lysis Buffer:

50 mM HEPES pH 7.5

300 mM NaCl

2 mM MgCl₂

Freshly added: 2 mM β-Mercaptoethanol, Roche Protease Inhibitor (Merck, #5056489001), and DNase I (Sigma, #DN25-1G)

Wash Buffer 1:

50 mM HEPES pH 7.5

300 mM NaCl

Freshly added: 1 mM DTT

Wash Buffer 2:

50 mM HEPES pH 7.5

700 mM NaCl

Freshly added: 1 mM DTT

SEC Buffer:




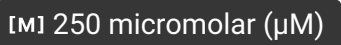


25 mM HEPES pH 7.5

150 mM NaCl






Freshly added: 1 mM DTT

Expression

16h 15m






- 1 Grow *E. coli* Rosetta (DE3) pLysS cells in  2 L LB medium at  37 °C until achieving an OD_{600 nm} of 0.4.
- 2 Reduce temperature to  18 °C and continue growth to OD_{600 nm} of 0.8.
- 3 Induce protein expression with  250 micromolar (μM) IPTG and continue growth for ~
 16:00:00 at  18 °C

16h

- 4 Centrifuge the cells at 3000 rcf ( 4 °C ,  00:15:00) 15m
- 5 Aspirate media and resuspend cell pellet in ~ 10 ml of 1X PBS, centrifuge the cells at 3000 rcf ( 4 °C ,  00:15:00) and take off the supernatant. 15m
- 6 Flash freeze the pellet in liquid nitrogen and store at  -80 °C until use.

Purification

46m 30s

- 7 Thaw pellet and resuspend in Lysis Buffer with freshly added  2 millimolar (mM) β -Mercaptoethanol, Roche Protease Inhibitor and DNase I
- 8 Sonicate sample  00:00:30 at 65% power for 5 cycles (Bandelin Sonopuls) 1m 30s
- 8.1 Repeat sonication for a total of 3X
- 9 Centrifuge at 48000 rcf at  4 °C for  00:45:00 45m
- 9.1 During this step, equilibrate  3 mL Glutathione Sepharose 4B beads by washing the slurry with water and Wash Buffer 1 (~ 6 ml).

10 Filter supernatant after centrifugation through a \rightarrow 0.45 μ m syringe filter

11 Incubate the sample on equilibrated beads on a tube roller at 4°C for 04:00:00 4h

12 Wash beads 5X with Wash Buffer 1

13 Wash beads 1X with Wash Buffer 2

Note

These wash steps remove nonspecific proteins



14 Wash beads 2X with Wash Buffer 1

Note

These washes remove high salt concentration

15 Cleave the resulting sample to produce GFP-NEMO by incubating with thrombin Overnight 4°C 1d

16 Centrifuge beads at 3000 rcf for 00:03:00 at 4°C 3m

- 17 Collect supernatant and filter through  .45 µm syringe filter
- 18 Concentrate product to a final volume of  500 µL using a 10 kDa MWCO concentrator
- 19 The concentrated and filtered protein can be applied onto a Superose 6 increase column (10/300 Cytiva) pre-equilibrated with SEC buffer
- 20 Fractions containing the purified proteins are pooled, concentrated, frozen in liquid nitrogen and stored at - 80°C

Note

This step can remove contamination and degradation products

Note

Protein purity can be determined by SDS PAGE analysis and measurement of protein concentration