



MAR 07, 2024

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.x54v9p9ypg3e/v1

Protocol Citation: Dominik Buchner 2024. Reconditioning PCR for removal of PCR bubbles in Illumina libraries. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.x54v9p9ypg3e/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Mar 07, 2024

🌐 Reconditioning PCR for removal of PCR bubbles in Illumina libraries

Dominik Buchner¹

¹University of Duisburg-Essen, Aquatic Ecosystem Research



Dominik Buchner

University of Duisburg-Essen, Aquatic Ecosystem Research

ABSTRACT

This protocol describes how to remove partly single-stranded PCR products ("PCR bubbles") from Illumina libraries. For more information about this phenomenon please see [this explanation from Illumina](#). For metabarcoding libraries, it can be hard to estimate the optimal input template or the perfect amount of PCR cycles and therefore overamplification happens frequently. PCR bubbles cannot be quantified reliably with fluorometric-based methods and may look different on different capillary electrophoresis devices.

PCR bubbles can lead to failed sequencing runs due to over- or underloading the flowcell.

GUIDELINES

Follow general lab etiquette. Wear gloves to prevent contamination of samples. Clean the workspace before starting and after finishing with 80% EtOH.

MATERIALS

Materials required:

Below all materials needed for the protocol are listed. Vendors and part numbers are listed but interchangeable depending on the supply situation.

Chemicals:

⊗ QIAGEN Multiplex PCR Plus Kit **Qiagen Catalog #206152**

Primers:

Illumina P5 5' - AATGATACGGCGACCACCGAGATCT - 3'

Illumina P7 5' - CAAGCAGAAGACGGCATACGAGAT - 3'




Buffers containing guanidine produce highly reactive compounds when mixed with bleach. Don't mix the extraction waste with bleach or solutions that contain bleach.


Reagents are potentially damaging to the environment. Dispose waste as mandated.

Quality control


- 1 Perform a quality control of your library. PCR bubbles may look different depending on the method used for quality control. See below and example of the Fragment Analyzer, Bioanalyzer and an agarose gel.

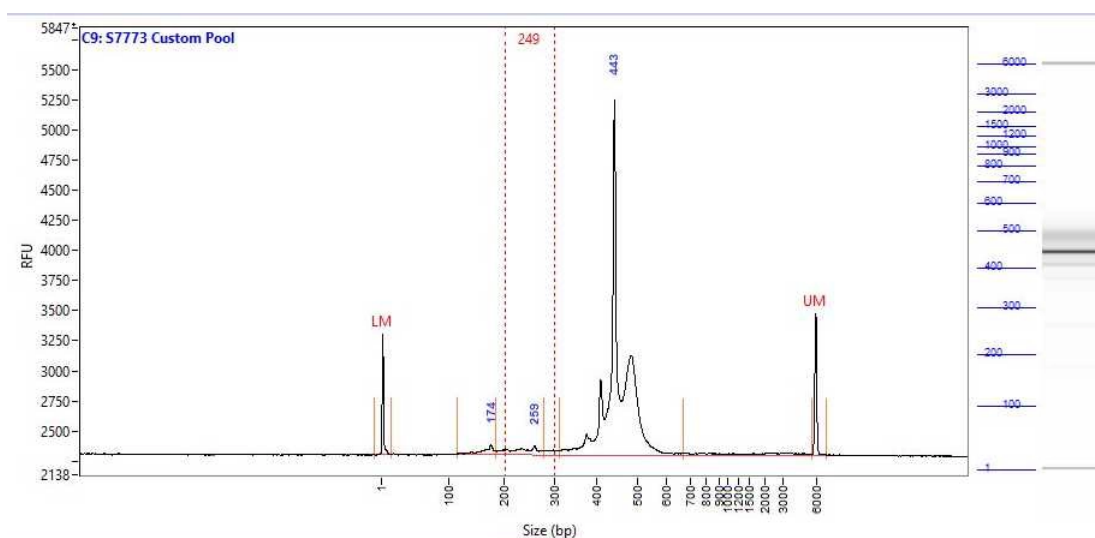

Example result of the Agilent fragment analyzer for a library that contains PCR bubbles. The desired peak is at 443 bp, the shoulder to the right of it is the PCR bubble.

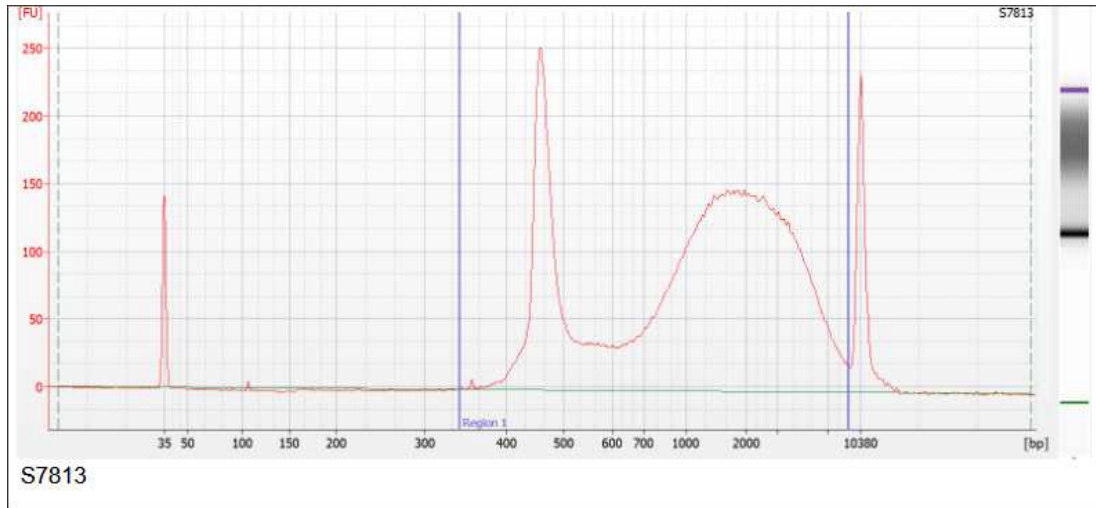


The same library on the Agilent bioanalyzer. The shoulder is moved further to the right, also the relationship of library to PCR bubble changed significantly.

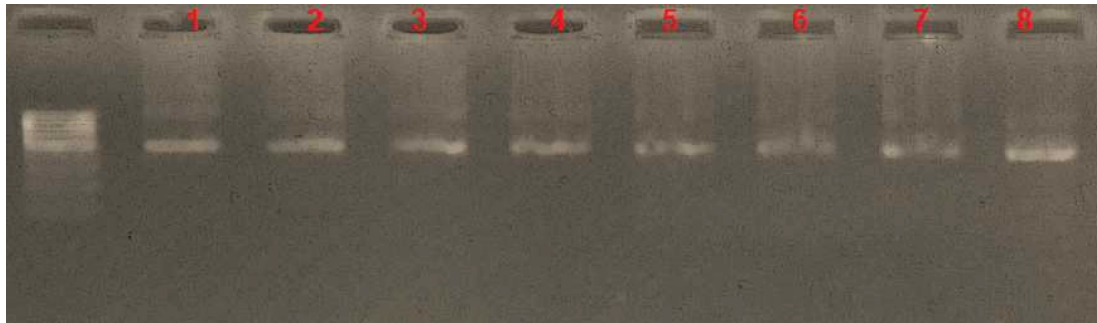

PCR reactions of the same library visualized on 1% agarose. Notice the faint band above the actual amplicon.

Library concentration (optional)

- 2 Concentrate your library by reducing the volume down to  100 μL . We usually do this with a spin-column based protocol, although this can be performed with magnetic beads as well.

Note

Please see:

Protocol



NAME

PCR cleanup and size selection with magnetic beads

CREATED BY

Dominik Buchner

PREVIEW

or

Protocol



NAME


Guanidine-based DNA extraction with silica-coated beads or silica spin columns

CREATED BY

Dominik Buchner

PREVIEW

Reconditioning PCR

- 3 Fill the concentration of your library and the project name into the Excel spreadsheet. The suggested master mix for the reconditioning PCR will be calculated accordingly. We usually go for  1250 ng of template input, however, this can be adjusted if necessary. master mix

Note

You can download the Excel spreadsheet here:



Mastermix calculator reconditioning PC...

4 Perform the PCR with 4 reactions of $50\ \mu\text{L}$.

PCR clean-up

5 Pool the 4 PCR reactions.

6 Perform a column-based PCR clean-up to exchange the buffer. This can also be done with magnetic beads. Elute the DNA in $100\ \mu\text{L}$.

Protocol



NAME


Guanidine-based DNA extraction with silica-coated beads or silica spin columns

CREATED BY

Dominik Buchner

PREVIEW

Size-selection

- 7 Perform a size selection with a ratio of 0.7x to remove residual primer dimers. Elute the final library in  50 µL .

Protocol



NAME

PCR cleanup and size selection with magnetic beads


CREATED BY

Dominik Buchner

PREVIEW


Perform final quality control

- 8 Perform a final quality control. Quantify the library concentration with a fluorometric-based method and perform quality control via electrophoresis. The shoulder should be gone and the library should be good for sequencing.

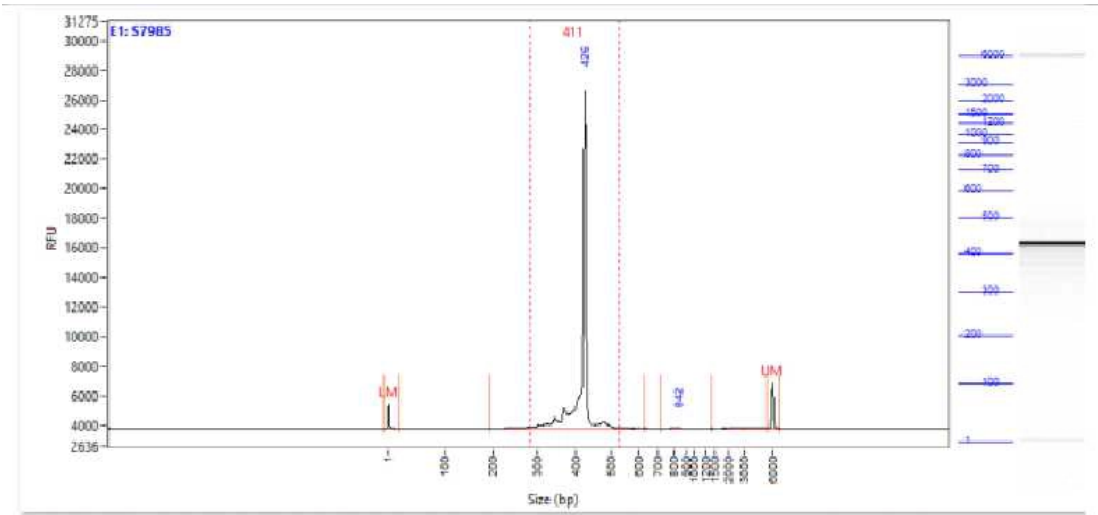

S7985 Fragment Analyzer


Example result of the Agilent fragment analyzer for a library after the removal of PCR bubbles.

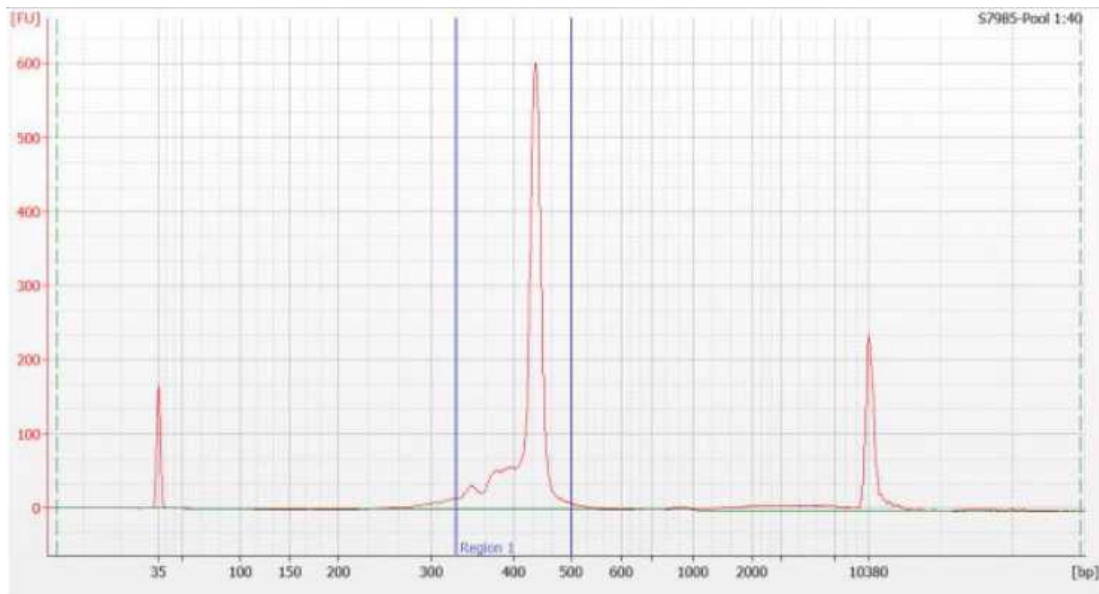

S7985 Bioanalyzer

The same library on the Agilent bioanalyzer after the removal of PCR bubbles.

Expected result

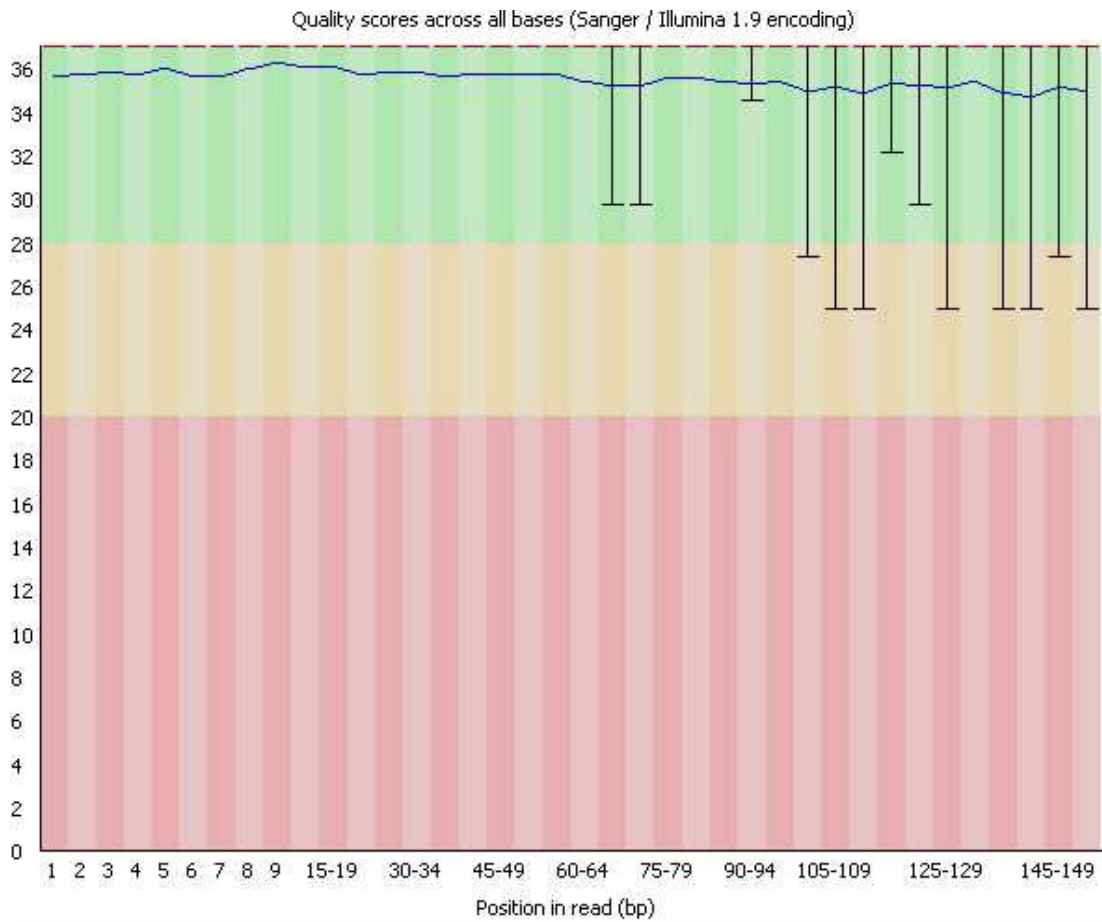

The quality of the sequencing run should be high after the reconditioning PCR.