

JAN 04, 2023

## OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.3byl4jzyjlo5/v1](https://dx.doi.org/10.17504/protocols.io.3byl4jzyjlo5/v1)

**Protocol Citation:** Javiera A Avilés, Tamara Matute, Isaac Núñez, Maira Rivera, Javiera Reyes, Jenny Molloy, Cesar A Ramirez-Sarmiento, Fernan Federici 2023. Recombinant expression and purification of HIV-1 RT. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.3byl4jzyjlo5/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jan 04, 2023

**Last Modified:** Jan 04, 2023

**PROTOCOL integer ID:**  
 74741

## Recombinant expression and purification of HIV-1 RT

Javiera A Avilés<sup>1</sup>, Tamara Matute<sup>2</sup>, Isaac Núñez<sup>2</sup>, Maira Rivera<sup>2</sup>, Javiera Reyes<sup>2</sup>, Jenny Molloy<sup>3</sup>, Cesar A Ramirez-Sarmiento<sup>2</sup>, Fernan Federici<sup>1</sup>

<sup>1</sup>Millennium Institute for Integrative Biology (iBio), Santiago, Chile;

<sup>2</sup>Institute for Biological and Medical Engineering, Pontificia Universidad Católica de Chile;

<sup>3</sup>Open Bioeconomy Lab, University of Cambridge

Reclone.org (The Reagent Collaboration Network)  
 Tech. support email: [protocols@reclone.org](mailto:protocols@reclone.org)  
[Click here to message tech. support](#)



Cesar A Ramirez-Sarmiento  
 Pontificia Universidad Catolica de Chile

## ABSTRACT

This protocol has been optimized for the recombinant expression of HIV-1 RT. The sequence of the plasmid encoding HIV-1 RT can be found on reclone.org.

The goal of this protocol was to eliminate the use of large volumes for dialysis and its fast buffer exchange into storage conditions.

## MATERIALS TEXT

### MATERIALS

- ☒ Sodium phosphate monobasic monohydrate **Sigma Aldrich Catalog #S9638**
- ☒ PMSF **Sigma Aldrich Catalog #P7626**
- ☒ Sodium phosphate dibasic **Sigma Aldrich Catalog #7558-79-4**
- ☒ DTT **Sigma Aldrich Catalog #D0632**
- ☒ Imidazole **Sigma Catalog #I5513**
- ☒ NaCl **Sigma Aldrich Catalog #53014**
- ☒ HisTrap FF Crude Column **Ge Healthcare Catalog #17528601**
- ☒ Glycerol **Sigma Aldrich Catalog #G5516**
- ☒ Dextrose **Sigma – Aldrich Catalog #D9434**
- ☒ Tween-20 **Sigma Aldrich Catalog #P9416**
- ☒ EDTA **Sigma Aldrich Catalog #ED2SS**
- ☒ Trizma Base **Contributed by users Catalog #93362**
- ☒ Lysozyme **Thermo Fisher Scientific Catalog #89833**
- ☒ Amicon Ultra-15 Centrifugal Filter Unit **Emd Millipore Catalog #UFC910024**

### Buffer A, pH 8.0

[M] 50 millimolar (mM) NaPO4, pH 8.0

[M] 50 millimolar (mM) dextrose  
 [M] 300 millimolar (mM) NaCl  
 [M] 1 millimolar (mM) EDTA  
 [M] 0.1 % volume Tween-20  
 [M] 40 millimolar (mM) Imidazole, pH 8.0

#### Buffer B, pH 8.0

[M] 50 millimolar (mM) Tris-HCl, pH 8.0  
 [M] 300 millimolar (mM) NaCl  
 [M] 1 millimolar (mM) EDTA  
 [M] 0.1 % volume Tween-20  
 [M] 10 % volume Glycerol  
 [M] 40 millimolar (mM) Imidazole, pH 8.0

#### Buffer C, pH 8.0

[M] 50 millimolar (mM) Tris-HCl, pH 8.0  
 [M] 300 millimolar (mM) NaCl  
 [M] 1 millimolar (mM) EDTA  
 [M] 0.1 % volume Tween-20  
 [M] 10 % volume Glycerol  
 [M] 150 millimolar (mM) Imidazole, pH 8.0

#### Buffer D, pH 8.0

[M] 50 millimolar (mM) Tris-HCl, pH 8.0  
 [M] 300 millimolar (mM) NaCl  
 [M] 1 millimolar (mM) EDTA  
 [M] 0.1 % volume Tween-20

#### Storage Conditions

[M] 25 millimolar (mM) Tris-HCl, pH 8.0  
 [M] 150 millimolar (mM) NaCl  
 [M] 0.5 millimolar (mM) EDTA  
 [M] 5 millimolar (mM) DTT  
 [M] 0.05 % volume Tween-20  
 [M] 50 % volume Glycerol

## DAY 1 - Plasmid transformation

1d

- 1 Transform  100 ng of plasmid containing HIV-1 RT into *E. coli*/BL21(DE3) competent cells using either heat shock or electroporation.

2h

- 2 Spread transformed cells in LB Agar plates supplemented with **1mM 0.05 mg/mL Kan**. Grow plate overnight at **37 °C**. 12h

## DAY 2 - Preinoculum

1d

- 3 Select a single colony from the LB agar plate to prepare a preinoculum in **10 mL** LB media supplemented with **1mM 0.05 mg/mL Kan**. Grow overnight at **250 rpm, 37°C**. 1d

## DAY 3 - Protein Overexpression











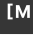

1d

- 4 Use the full volume of the preinoculum to inoculate **1 L** of LB media supplemented with **1mM 0.05 mg/mL Kan** (1% inoculation). Grow at **250 rpm, 37°C** until reaching an optical density at 600 nm ( $OD_{600}$ ) = 0.8. 4h
- 5 Upon reaching  $OD_{600}$  = 0.8, add **0.5 millimolar (mM) IPTG** and incubate for 2h at **220 rpm, 37°C**. 2h

## DAY 4 - Protein Purification by IMAC

6h

- 6 Centrifuge the cell culture **4000 x g, 4°C, 00:20:00**. Then, resuspend the cell pellet in **50 mL** of **Buffer A** freshly supplemented with **0.5 millimolar (mM) PMSF** and **0.2 mg/mL lysozyme**. 30m
- 7 Incubate the resuspended cells **80 rpm, Room temperature, 00:30:00**. 30m
- 8 Sonicate on ice for **00:04:00** using cycles of **00:00:06** ON and **00:00:06** OFF at 40% amplitude (Qsonica Q125, 125W). 30m
- 9 Centrifuge the unclarified lysate **20000 x g, 4°C, 00:20:00** and collect the supernatant. You might want to collect a small sample for SDS-PAGE afterwards. 30m

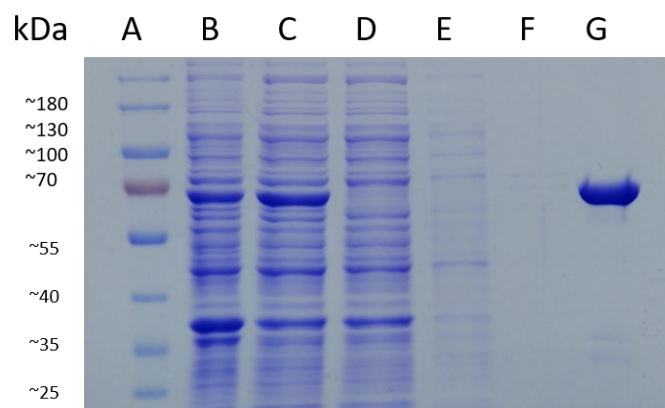
- 10 On a **1 mL HisTrap column (Ge Healthcare)** preequilibrated with 10 column volumes (c.v.) (here, 10 mL) of **Buffer A**, load the supernatant. Wash with 10 c.v. of **Buffer A**. Then, wash with 10 c.v. of **Buffer B**, and elute with 5 c.v. of **Buffer C**, collecting the eluted fractions every  0.5 mL in 1.5 ml tubes. 1h
- 11 To quickly pool the fractions containing the protein of interest, prepare a 96-well plate or 1.5 mL tubes with  40 µL of Bradford reagent and  160 µL of distilled water. Then, add  10 µL of each protein fraction and compare against a blank reference sample corresponding to  10 µL of **Buffer C**. You can determine your protein-containing fractions either by absorbance at 595 nm on a plate reader or visually by comparing the blue coloration of each fraction against the blank reference. Pool your fractions and collect a  10 µL sample for SDS-PAGE 5m
- 12 To decrease the imidazole concentration, perform a buffer exchange step with an Amicon Ultra-15 concentrator (Merck Millipore). Centrifuge  3000 x g, 10°C, 00:10:00, discard the flowthrough, add Buffer D to decrease the imidazole concentration and repeat this step, until the imidazole concentration reaches < 30 mM. 40m
- 13 Recover the concentrated protein and determine its concentration using the Bradford assay. Also, collect a  10 µL sample for SDS-PAGE. 5m
- 14 For storage, supplement your pooled fraction with  10 millimolar (mM) DTT. Then, add glycerol up to  50 % volume to reach **Storage Conditions**. Do consider that a final protein concentration of  1 mg/mL is appropriate for subsequent experiments. 5m
- 15 Generate  200 µL aliquots of the enzyme and store it at  -20 °C until required. 30m

## IMAC SDS-PAGE Result

10m

16

Expected result



SDS PAGE for the recombinant expression and purification of HIV-1 RT in *E. coli*/BL21 (DE3) using the open pTi vector (A) Protein marker; (B) Pellet (C) Clarified cell lysate; (D) Flowthrough; (E) First washing step; (F) Second washing step; (G) Pooled eluted fractions