

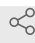


Sep 01, 2022

Overall protocol for MicroPOTS LCMS top down proteomics of kidney tissue sections

Mowei Zhou¹, James M Fulcher¹, Yen-Chen Liao¹, Ljiljana.PasaTolic¹¹Pacific Northwest National Laboratory

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.eq2lynm1qvx9/v1

Human BioMolecular Atlas Program (HuBMAP) Method Development Community

PNNL-TTD



Mowei Zhou
Pacific Northwest National Laboratory

ABSTRACT

This is the overall workflow for LCMS top down proteomics of kidney functional units from tissue sections using the MicroPOTS platform. The expected outcomes are proteoform identification and quantitation values from selected tissue functional units.

DOI

dx.doi.org/10.17504/protocols.io.eq2lynm1qvx9/v1

PROTOCOL CITATION

Mowei Zhou, James M Fulcher, Yen-Chen Liao, Ljiljana.PasaTolic 2022. Overall protocol for MicroPOTS LCMS top down proteomics of kidney tissue sections.
protocols.io
<https://protocols.io/view/overall-protocol-for-micropots-lcms-top-down-prote-cfhftj3n>



FUNDERS ACKNOWLEDGEMENT

National Institutes of Health (NIH) Common Fund, Human Biomolecular Atlas Program (HuBMAP)
Grant ID: UG3CA256959-01

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 19, 2022

LAST MODIFIED

Sep 01, 2022

PROTOCOL INTEGER ID

68871

Tissue collection

- 1 The tissue sections were prepared and shipped from Vanderbilt-TMC following the protocol below:



Cryostat Sectioning of Tissues for 3D Multimodal
Molecular Imaging
by **Jamie Allen,**
Vanderbilt University

PREVIEW

RUN



Sample preparation

- 2 Functional units (glomerulus, medullary, tubule) were dissected and collected into the microPOTS platform using the method below:



Laser Capture Microdissection of Tissue Functional Units
for microPOTS Top-Down Proteomics
by **James M Fulcher,**
Pacific Northwest National lab

PREVIEW

RUN



Data Acquisition

- 3 The samples were analyzed by LCMS top down proteomics as described below:



Top Down Proteomics Data Collection for Microdissected
Kidney Tissue Functional Units
by **Mowei Zhou,**
Pacific Northwest National Laboratory

PREVIEW

RUN



Data Analysis

- 4 LCMS datasets were analyzed for proteoform identification and quantitation. The final results are reported.



Proteoform Identification and Quantitation with TopPIC
and TDPortal for Human Tissues

by **Mowei Zhou,**

Pacific Northwest National Laboratory

PREVIEW

RUN

