



AUG 08, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.kxygx392dg8j/v1](https://dx.doi.org/10.17504/protocols.io.kxygx392dg8j/v1)

**Protocol Citation:** Erika Lara Flores, Joel Reyes, Mark Cookson, Michael Ward 2023. iNDI Papain Dissociation protocol for iNeurons Version 1. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.kxygx392dg8j/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Aug 08, 2023


# iNDI Papain Dissociation protocol for iNeurons Version 1

Erika Lara Flores<sup>1</sup>, Mark Cookson<sup>1</sup>, Joel Reyes<sup>2</sup>, Michael Ward<sup>1</sup>

<sup>1</sup>National Institutes of Health, NIH Center for Alzheimer's and Related Dementias (CARD);

<sup>2</sup>National Institutes of Health, National Institute of Neurological Disorders and Stroke (NINDS)

NIH Center for Alzheimer's and Related Dementias

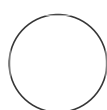



Erika Lara Flores

## ABSTRACT

Proteolytic enzymes are widely used in cell dissociation and papain has proved less damaging with some tissues and neuronal cultures and more effective than other proteases. Here we describe a papain dissociation protocol that has been used on our iPSC-derived neurons for scRNAseq and FACS assays readouts.

## MATERIALS

 Papain Dissociation System Worthington Biochemical Corporation Catalog  
 #LK003153






 TrypLE™ Express Enzyme (1X) no phenol red Thermo Fisher Scientific Catalog  
 #12604013

PBS

PBS-EDTA 0.05mM





PBS-0.04% BSA (For scRNAseq cell preparation, if applicable)

## Preparation of Reagents

- 1 **Equilibrate EBSS/phenol red** with 95%O<sub>2</sub>:5%CO<sub>2</sub> in incubator for several hours of overnight.
- 2 **Reconstitute the Ovomucoid mixture (Papain inhibitor-PI)** by adding  32 mL of EBSS and allow contents to dissolve. This will yield solution at an effective concentration of 10mg of ovomucoid inhibitor and 10mg of albumin per mL. Reconstitute for the first use, then store  1 mL at  4 °C .
- 3 **Papain solution:** Add  10 mL of TrypLE™ Express Enzyme (1X), no phenol red to the vial and place it in  37 °C bath for 10 minutes, this vial has 10 units/ml.




### Note

*When working with i3Neurons (WTC11 line), dissolve the papain with 20 mL of TrypLE™ Express Enzyme (1X).*

- 4 **DNase I vial D2:** Add 500ul of EBSS media and mix gently (2000 units/ml).
- 5 **Trituration Medium:**  
 30 mL of Neuronal Maturation medium (the same medium used for your differentiation)  
 30 µL of Rock inhibitor (1x)  
1 vial of DNase I
- 6 **Papain Inhibitor media:**  9 mL of Trituration medium +  1 mL of papain inhibitor.

## Dissociation protocol for iNeurons at day 28 of differentiat...

5m

- 7 Aspirate culture medium and wash gently with PBS (calcium/magnesium free).
- 8 Aspirate PBS and wash 2 times with PBS-0.5mM EDTA
- 9 Aspirate PBS-EDTA and add to one well of 6 well plate  1 mL of Papain solution (10 units/ml) to cells, and incubate at  37 °C for  00:05:00 .

5m

### Note


*Incubation time can vary and go up to 30 minutes and it depends on cell line, confluency and maturation day.*







*Incubation times for i3Neurons (WTC11):*

*Day 7: 3 minutes*

*Day 10: 5 minutes*

*Day 14: 10 minutes*

- 10 Aspirate papain solution when the cell bodies have an halo morphology like for and EDTA split. If the incubation is longer (up to 30 minutes) the cells will detach into papain solution, then continue with the next step.
- 11 Add  1 mL of trituration medium per well of a 6-well plate to dissociate cells and pipette medium over the plate up and down to detach cells and dissociate them as single-cell suspension.
- 12 Collect cells into a sterile conical 15 mL tube and adjust volume to 5 mL with PBS or trituration medium.

- 13 Centrifuge cells  00:05:00 at 200 - 300 x g at  Room temperature . 5m
- 14 Remove supernatant and resuspend the cell pellet gently with  1 mL of Papain inhibitor medium to wash off the papain.
- 15 Centrifuge cells  00:05:00 at 200 - 300 x g at  Room temperature 5m
- 16 Aspirate supernatant and resuspend cell pellet in  1 mL of Neuronal Maturation Medium.
- 17 After resuspending the cell pellet, count cells and use as necessary.

#### Note

For scRNAseq 10X genomics assay resuspend in cold PBS-0.04%BSA.  
For FACS assay resuspend in cold PBS-0.5mM EDTA.