



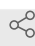
Version 3 ▼

Aug 24, 2022

# DNA extraction from fecal samples V.3

Yoshiyuki Matsuo<sup>1</sup><sup>1</sup>Kansai Medical University

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.3byl4k912vo5/v3](https://dx.doi.org/10.17504/protocols.io.3byl4k912vo5/v3) Yoshiyuki Matsuo

## ABSTRACT

DNA extraction method for metagenomic sequencing of the gut microbiota

## DOI

[dx.doi.org/10.17504/protocols.io.3byl4k912vo5/v3](https://dx.doi.org/10.17504/protocols.io.3byl4k912vo5/v3)

## PROTOCOL CITATION

Yoshiyuki Matsuo 2022. DNA extraction from fecal samples. **protocols.io**  
<https://protocols.io/view/dna-extraction-from-fecal-samples-b9qjr5un>  
Version created by Yoshiyuki Matsuo



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

May 20, 2022

## LAST MODIFIED

Aug 24, 2022

## PROTOCOL INTEGER ID

62955

## MATERIALS TEXT

### Reagents/Kits

- Phosphate-buffered saline (PBS)
- EZ-Beads (Promega/AMR, AMR76813M)
- Maxwell RSC Blood DNA Kit (Promega, AS1400)

### Equipment

- Vortex mixer
- High-speed microcentrifuge
- Block heater
- Micro Smash Beads Cell Disrupter (TOMY Digital Biology, MS-100)
- Maxwell RSC instrument (Promega, AS4500)

### Preparation of fecal samples

1 Place  **50-100 mg** of fecal sample into tube.

2 

Add  **1 mL** of PBS per  **100 mg** of feces.

3 


Mix thoroughly by vortexing.

4 

2m

Allow the sample to stand for  **00:02:00** to sediment large debris.

5 

Transfer  **300 µL** of the suspension to 1.5 mL tube.

6



3m

Centrifuge at **10000 x g, 00:03:00** .

7



Discard the supernatant.

8



Resuspend the pellet (~ **30 mg** of feces) in **300 µL** of PBS.

9



10m

Incubate at **70 °C** for **00:10:00** on the block heater.

10

Cool to **Room temperature** .

#### Mechanical cell disruption by bead beating

11




Transfer **300 µL** of the suspension to EZ-beads tube.

The EZ-Beads tube contains zirconium oxide beads of two different sizes (0.2 mm spheres and a large 5 mm bead) that can facilitate efficient cell lysis by bead beating.

12




Lyse cells either by using disruption device (12.1) or vortex mixer (12.2).

- 12.1 Place the EZ-beads tube in Micro Smash instrument and disrupt cells by shaking at  **2500 rpm, 00:02:00** .

Micro Smash Beads Cell Disrupter  
TOMY Digital Biology MS-100



Caution: Avoid using a disruption device with a high-speed linear reciprocating motion, as this may potentially result in breakage of the EZ-Beads tubes.

- 12.2 Place the EZ-beads tube on MN Bead Tube Holder attached to Vortex-Genie<sup>5m</sup> mixer and vortex for  **00:05:00** at maximum speed.

MN Bead Tube Holder  
Rubber-foam adapter for processing bead tubes  
with Vortex-Genie instrument  
**MACHEREY-NAGEL 740469**

13





Briefly spin the tube to collect contents.

Automated DNA extraction using Maxwell RSC Blood DNA Kit

23m

14



Add  **300 µL** of Lysis Buffer and  **30 µL** of Proteinase K Solution to the sample in EZ-beads tube.

15



Mix by inverting the tube.

16 

Briefly spin the tube.

17 

20m

Incubate at  **56 °C** for  **00:20:00** on the block heater.

18 

Briefly spin the tube.

19 

Transfer the supernatant (~  **500 µL** ) to 1.5 mL tube.

20 

3m

Centrifuge at  **18000 x g, 00:03:00** .

21 

Transfer the cleared lysate to Maxwell RSC Cartridge.

22 

Add  **50 µL** of Elution Buffer to elution tube.

23 Start the extraction run following the manufacturer's instructions.

Maxwell RSC instrument  
Automated nucleic acid purification platform  
Promega AS4500