

B



Sep 08, 2020

SPOT2 protocol

Huimin Zhao¹, stlane², Guanhua Xun¹

¹University of Illinois at Urbana-Champaign

1 Works for me

dx.doi.org/10.17504/protocols.io.bk28kyhw

stlane2

DOI

dx.doi.org/10.17504/protocols.io.bk28kyhw

PROTOCOL CITATION

Huimin Zhao, stlane2, Guanhua Xun 2020. SPOT2 protocol. **protocols.io** https://dx.doi.org/10.17504/protocols.io.bk28kyhw

LICENSE

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 08, 2020

LAST MODIFIED

Sep 08, 2020

PROTOCOL INTEGER ID

41792

STEPS MATERIALS

NAME	CATALOG #	VENDOR
Magnesium Sulfate (MgSO4) Solution - 6.0 ml	B1003S	New England Biolabs
E gene primer mix	N/A	
Non-primer oligos	N/A	
Bst 2.0 WarmStart DNA Polymerase - 8,000 units	M0538L	New England Biolabs
WarmStart RTx Reverse Transcriptase - 250 rxns	M0380L	New England Biolabs
Isothermal Amplification Buffer - 6.0 ml	B0537S	New England Biolabs
Nuclease-free Water - 25 ml	B1500S	New England Biolabs
Deoxynucleotide Solution Mix - 8 umol of each	N0447S	New England Biolabs
N gene primer mix	N/A	
Saliva sample	N/A	
Manganese(II) chloride tetrahydrate	M3634	Sigma Aldrich
Reporter probe 1	N/A	
Reporter Probe 2	N/A	
Non-CRISPR nuclease	N/A	

EQUIPMENT


NAME	CATALOG #	VENDOR
SPOT2 Device	N/A	
USB Laser Barcode Scanner	WN3300	Amazon
Thermal Immersion Cooker	P00252	

protocols.io

09/08/2020


 $\textbf{Citation:} \ \ \textbf{Huimin Zhao, stlane2} \ \textbf{,} \ \textbf{Guanhua Xun (09/08/2020)}. \ \textbf{SPOT2 protocol.} \ \underline{\textbf{https://dx.doi.org/10.17504/protocols.io.bk28kyhw}}$

- 1 Label test kit bag with patient's name and a unique sample number for the test run. Provide a 20 uL microcap and PCR tube and have patient collect a saliva sample and place it into then seal the PCR tube. Label PCR tube with sample number matching the test kit bag.
- Place PCR tubes into floating rack and insert into sous vide water bath at 95C. Incubate for 5 minutes for sample pretreatment.

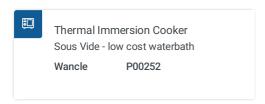

8 95 °C © 00:05:00

2.1 During incubation, scan barcodes on test kits and enter patient identifying information into the SPOT2 device.

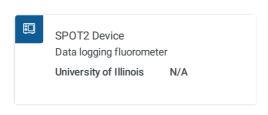



3 After completion of pretreatment, use the second provided 20 uL microcap to transfer a small volume of saliva from the PCR tube to the capillary tube, containing SPOT assay mastermix. Place capillary tubes into floating rack and insert into sous vide water bath at 63C for 30 minutes.

8 63 °C © 00:30:00

SPOT assay mastermix:

	Initial concentration	Final concentration	Amount (µL)
Upper compartment			
WarmStart® Bst 2.0	8000 units/mL	160 units/mL	2
WarmStart® RTx	15,000 units/mL	150 units/mL	1
Isothermal amplification buffer	10X	0.5X	8
dNTPs	10 mM	0.7 mM	5.6
MgSO4	100 mm	4 mM	3.2

N gene primer mix	10X	0.25X	2
E gene primer mix	10X	0.25X	2
Saliva samples			5
Non-CRISPR nuclease	5 mg/mL or 55 μM	1.375 uM	2
MnCl2	50 mM	0.5 mM	0.8
Non-primer oligos (total 6 oligos)	100 μΜ	625 nM	3
Reporter probe 1	100 μΜ	156.25 nM	0.125
Reporter probe 2	100 μΜ	312.5 uM	0.25
Nuclease-free water			44.025
Total			80

83

Bst 2.0 WarmStart DNA Polymerase -8,000 units

by New England Biolabs Catalog #: M0538L

⊒2 µl

WarmStart RTx Reverse Transcriptase -250 rxns

by New England Biolabs Catalog #: M0380L

□1 μl

88

Isothermal Amplification Buffer - 6.0 ml

by New England Biolabs Catalog #: B0537S

⊒8 μl

Deoxynucleotide Solution Mix - 8 umol of each

by New England Biolabs Catalog #: N0447S

■5.6 μl [M]0.7 Milimolar (mM)

Magnesium Sulfate (MgSO4) Solution -6.0 ml

by New England Biolabs Catalog #: B1003S

■3.2 µl [M]4 Milimolar (mM) 88 N gene primer mix Catalog #: N/A E gene primer mix Catalog #: N/A Saliva sample Catalog #: N/A Nuclease-free Water - 25 ml by New England Biolabs Catalog #: B1500S Non-CRISPR nuclease Catalog #: N/A Manganese(II) chloride tetrahydrate by Sigma Aldrich Catalog #: M3634 ■0.8 µl [M]0.5 Milimolar (mM) Non-primer oligos Catalog #: N/A 83

Reporter probe 1 Catalog #: N/A



4 Transfer floating capillary rack to sous vide at 95C for 5 minutes.

8 95 °C © 00:05:00

5 Individually transfer capillary tubes to SPOT2 device and press "Run" on the SPOT software. An integrated RFID chip on each capillary will automatically link test result to the patient information entered previously.

