



OCT 31, 2023

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.x54v9p3qzg3e/v1

Collection Citation: Ester Kalef-Ezra, Ben Harvey, Katherine Roper, Christos Proukakis 2023. SureSelect XT HS2 DNA to prepare libraries for single-cell Whole Genome Sequencing (scWGS) after single-cell Whole Genome Amplification (scWGA). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.x54v9p3qzg3e/v1>

License: This is an open access collection distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this collection and it's working.

🌐 SureSelect XT HS2 DNA to prepare libraries for single-cell Whole Genome Sequencing (scWGS) after single-cell Whole Genome Amplification (scWGA)

Ester Kalef-Ezra^{1,2}, Ben Harvey³, Katherine Roper³, Christos Proukakis^{1,2}

¹Department of Clinical and Movement Neurosciences, UCL Queen Square Institute of Neurology, London, UK;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815;

³Diagnostics and Genomics Group, Agilent Technologies LDA UK Ltd

ASAP Collaborative Research Network

University College London



Ester Kalef-Ezra
University College London

Created: Oct 09, 2023

DISCLAIMER

Last Modified: Oct 31, 2023

COLLECTION integer ID:
89822

Keywords: ASAPCRN

Funders

Acknowledgement:

Aligning Science Across Parkinson's through the Michael J. Fox Foundation for Parkinson's Research (MJFF)
Grant ID: 000430

Acknowledgements and Funders:

We thank ICH.ZCR UCL Genomics Sequencing and Agilent for technical support.

This research was funded in part by Aligning Science Across Parkinson's [Grant ID: 000430] through the Michael J. Fox Foundation for Parkinson's Research (MJFF).

Note

ALL IMAGES CAN BE ENLARGED BY CLICKING AT THEM

Critical notes!

Please follow Good Laboratory Practices.

- To prevent DNA contamination, clean all surfaces and equipment before use with DNA AWAY Surface Decontaminant, separate the lab in pre- and post-PCR areas, and use filtered sterile pipette tips.
- For each protocol step that requires the removal of tube cap strips, reseal the tubes with a fresh strip of caps.
- In all master mixes to be used manually, include 5-10% excess volume in the calculations.

ABSTRACT

We adapted the SureSelect XT HS2 DNA protocol to prepare libraries using as input material samples after single-cell Whole Genome Amplification (scWGA), instead of genomic DNA samples, for library preparation for Illumina single-cell Whole Genome Sequencing (scWGS). This protocol can be employed either manually (Section 2: Option A) or in combination with automation using the Bravo Automated Liquid Handling Platform (Section 2: Option B).

ATTACHMENTS

[861-2221.pdf](#) [2022_SureSelect_Automation_G9985-90020.pdf](#) [SureSelect Enzymatic Fragmentation kit.pdf](#) [SureSelect XT HS2 DNA.pdf](#)


GUIDELINES

We have used this protocol for DNA products after single-cell Whole Genome Amplification (scWGA) by either PicoPLEX (Takara), ResolveDNA (BioSkryB) and droplet MDA (Samplix). We recommend to bead purify the PicoPLEX and ResolveDNA products prior library preparation. In all cases, the amplicon quantity and size should be assessed using Qubit and TapeStation or other methods prior to SureSelect XT HS2 DNA library preparation.


Bravo error/notification messages

Two error/notification messages are encountered when VWorks Bravo software is initializing:

1. If you encounter the G-axis error message shown below, select Ignore and Continue, leaving device in current state.

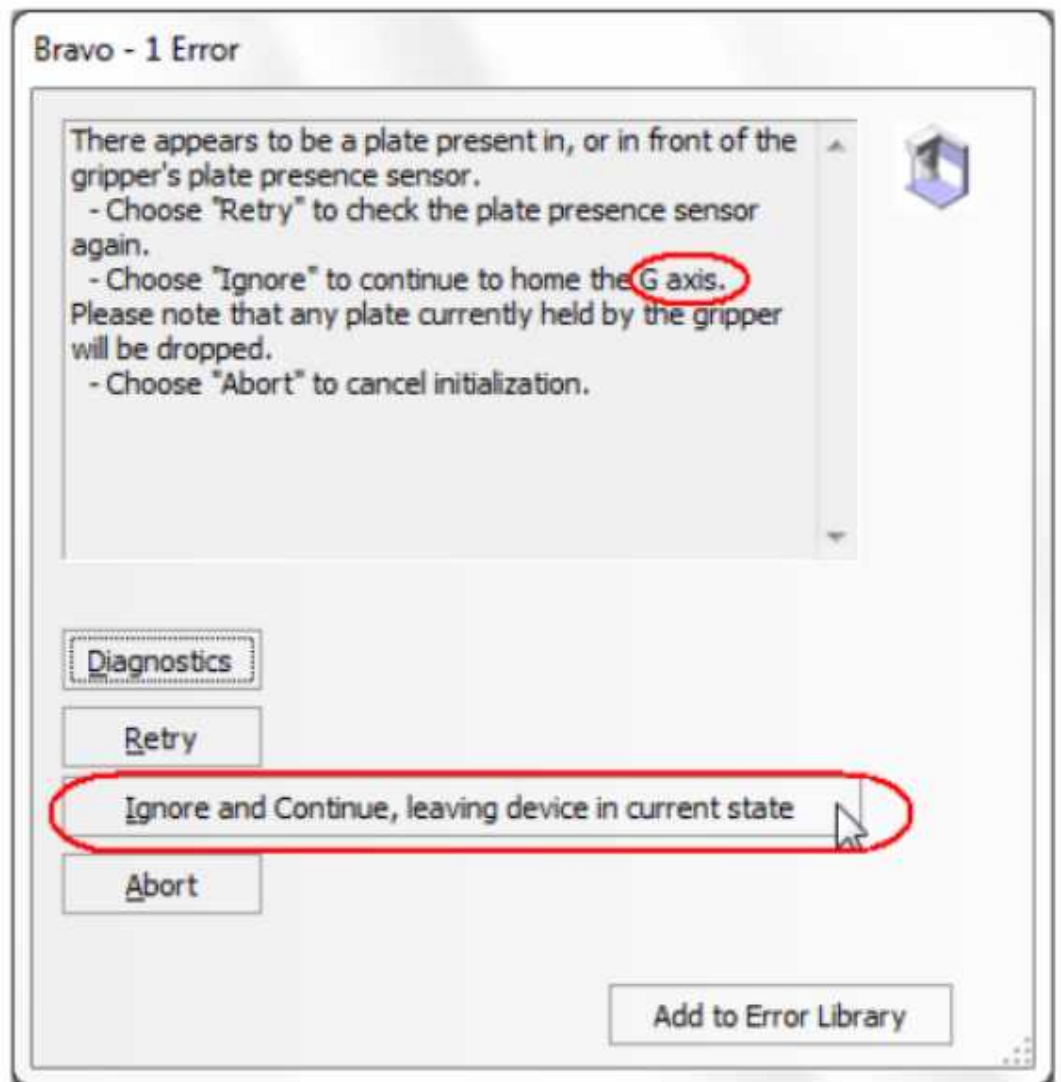

2. If you encounter the W-axis error message shown below, select Retry.

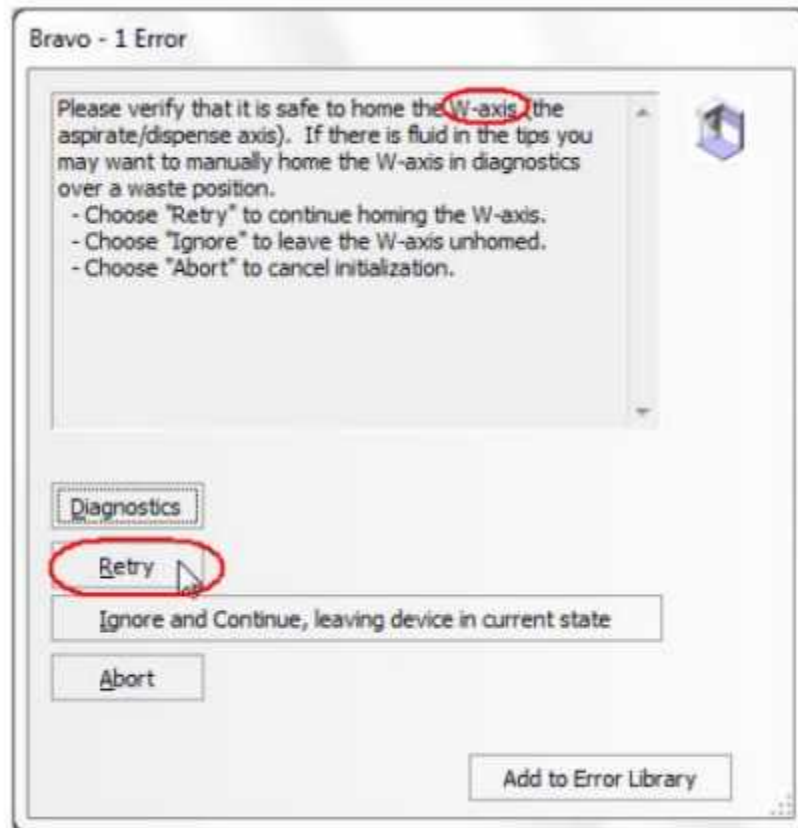




MATERIALS

Commercial Reagents:

1. Equipment and consumables for DNA quality control:







DNA quantification. Here we used Qubit (Thermo Fisher Scientific):

Equipment	
Qubit Fluorometer	NAME
Fluorometer	TYPE
Invitrogen	BRAND
Q33238	SKU
https://www.thermofisher.com/order/catalog/product/Q33238#/Q33238	LINK



-  Qubit™ dsDNA BR Assay Kit Thermo Fisher Scientific Catalog #Q32853
-  Qubit® dsDNA HS Assay Kit Thermo Fisher Scientific Catalog #Q32854
-  Qubit™ Assay Tubes Invitrogen - Thermo Fisher Catalog #Q32856

2. Nucleic acid analysis platforms. Here we used TapeStation (Agilent):

Equipment	
4150	NAME
TapeStation System	TYPE
Agilent	BRAND
G2992AA	SKU
https://www.agilent.com/store/en_US/Prod-G2992AA/G2992AA	LINK


-  Optical tube strips (8x Strip) Agilent Technologies Catalog #401428
-  Optical tube strip caps (8x strip) Agilent Technologies Catalog #401425
-  D1000 ScreenTape Agilent Technologies Catalog #5067-5582
-  D1000 Reagents Agilent Technologies Catalog #5067-5583
-  High Sensitivity D1000 ScreenTape Agilent Technologies Catalog #5067-5584
-  High Sensitivity D1000 Reagents Agilent Technologies Catalog #5067-5585

SureSelect XT HS2 performed manually:







-  SureSelect XT HS enzymatic fragmentation kit, 16 reactions Agilent Technologies Catalog #5191-4079
- or
-  SureSelect Enzymatic Fragmentation Kit, 96 reactions Agilent Technologies Catalog #5191-4080
 - SureSelect XT HS2 DNA Library Preparation Kit with Index Primer Pairs (Choose

one of the following: Agilent G9981A (index 1-16), G9985A (index 1-96), G9985B (index 97-192), G9985C (index 193-288) or G9985D (index 289-384))




SureSelect XT HS2 with automation (Bravo):


-  SureSelect Enzymatic Fragmentation Kit, 96 Reactions, Auto Agilent Technologies Catalog #5191-6764
- SureSelect XT HS2 DNA Library Preparation Kit with Index Primer Pairs (Choose one of the following: Agilent G9985A (index 1-96), G9985B (index 97-192), G9985C (index 193-288) or G9985D (index 289-384))


Consumables:


-  UltraPure[®] DNase/RNase-Free Distilled Water Thermo Fisher Catalog #10977049
-  Low TE Buffer Invitrogen - Thermo Fisher Catalog #12090-015
-  Ethanol (molecular biology grade, ≥99.8%) Merck MilliporeSigma (Sigma-Aldrich) Catalog #51976-500ML-F
- Plasticware compatible with the selected thermal cycler. Here we used PCR 8-tube strips w/o caps
 -  8-tube strips for qPCR/PCR, without caps VWR International Catalog #732-1517
- Low binding filtered tips (sterile)
- Low-binding 1.5ml DNase/RNase-free (sterile). Here we used DNA LoBind Tubes,
 -  DNA LoBind[®] Tubes Eppendorf Catalog #0030108418
- Powder-free gloves
-  DNA AWAY[®] Surface Decontaminant, Surface decontaminant; 8.5 oz. (250mL) Thermo Fisher Catalog #7010PK
- Cleaning wipes for surface cleaning. Here we used Conti Washcloth Dry Brosch Direct PH5959


Consumables for SureSelect XT HS2 with automation:

- Processing plate:
 -  Armadillo PCR Plate, 96-well, clear, clear wells, barcoded Thermo Fisher Scientific Catalog #BC2396
- PCR plate:
 -  PCR plates, 96-well, Armadillo[™] Standard Semi-skirted VWR International Catalog #731-0649
- Deep well plate for AMPure beads and master mixes:
 -  Nunc[®] DeepWell[®] Plates with Shared-Wall Technology, 96U, 1mL, sterile Thermo Fisher Catalog #260251
- Waste Plate:

 Axygen™ Storage Microplates Fisher Scientific Catalog #P-2ML-SQ-C


▪  VersiCap Mat, 96-well, domed cap strips Thermo Fisher Scientific Catalog #AB1810

▪  Agilent Bravo 250uL Pipette Tips, Sterile, Filtered, for 96LT head Contributed by users Catalog #19477-022

▪  Reservoir, single cavity, polypropylene, 300 mL, 96 pyramids base geometry, 4 mm height, 25/pk Agilent Technologies Catalog #201244-100

Equipment:

- Thermal Cycler with 96-well, 0.2 ml block. Here we used Multigene Optimax Gradient Thermal cycler with 96 well block (Labnet MB0520)
- Manual version: Magnetic separator compatible with 0.2 ml tubes, strips of plates. Here we used either Magnetic Separator - PCR Strip (Takara 635011) or ResolveDNA Dual Volume Strip Tube Magnet (BioSkrayB PN100226)
- Plate or strip tube centrifuge. Here we used SciSpin Mini microfuge, blue, 7000rpm (SciQuip SS-6050) and plate centrifuge (Starlab N2631-0008)
- 1.5 ml tube centrifuge. Here we used Mini-Centrifuge (Fisherbrand 16617645)
- PCR cooler. Here we used PCR-Cooler, Blue, Capacity: 0.2 mL (Eppendorf 1019228)
- General lab pipettes (single and multi-channel)
- Nuclease-free filtered pipette tips sterile, if possible, use low binding tips
- 96-well plate mixer. Here we used Vortex mixer. Here we used WhirliMixer (Fisons Scientific 1993-520) and iSwix VT Digital Vortex Mixer (Appleton Woods ST6000)
- Ice bucket
- For automated version: Bravo NGS Option A robot (Agilent G5573A)

 SMARTer-Seq™ Magnetic Separator - PCR Strip Takara Bio Inc. Catalog #635011

Equipment

SciSpin MINI Microfuge BLUE, inc 2 rotors of SS-6058 & SS-6059	NAME
SciQuip	BRAND
SS-6050	SKU
https://sciquip.co.uk/scispin-mini-microfuge.html	LINK

Equipment	
Plate Centrifuge	NAME
StarLab	BRAND
N2631-0008	SKU
https://www.starlabgroup.com/en/product/plate-centrifuge-n2631-0008.html	LINK

Equipment	
Mini-Centrifuge 100-240V, 50/60Hz Universal Plug, Grey	NAME
minicentrifuge	TYPE
Fisherbrand™	BRAND
16617645	SKU
https://www.fishersci.co.uk/shop/products/fisherbrand-standard-mini-centrifuge/16617645	LINK

Equipment	
iSwix VT Digital Vortex Mixer	NAME
Appleton	BRAND
ST6000	SKU
https://www.appletonwoods.co.uk/product/iswix-vt-digital-vortex-mixer-parent/	LINK

SAFETY WARNINGS



Safety Warnings:

Please follow the Safety Data Sheets (SDS) for all reagents for safe handling and safety hazards.

DISCLAIMER

Acknowledgements and Funders:

We thank ICH.ZCR UCL Genomics Sequencing and Agilent for technical support.

This research was funded in part by Aligning Science Across Parkinson's [Grant ID: 000430] through the Michael J. Fox Foundation for Parkinson's Research (MJFF).

Note

ALL IMAGES CAN BE ENLARGED BY CLICKING AT THEM

Critical notes!

Please follow Good Laboratory Practices.

- To prevent DNA contamination, clean all surfaces and equipment before use with DNA AWAY Surface Decontaminant, separate the lab in pre- and post-PCR areas, and use filtered sterile pipette tips.
- For each protocol step that requires the removal of tube cap strips, reseal the tubes with a fresh strip of caps.
- In all master mixes to be used manually, include 5-10% excess volume in the calculations.

ATTACHMENTS


861-
2221.pdf

2022_SureSe
lect_Auto...

SureSelect
Enzymatic...

SureSelect
XT HS2 DN...

Protocol




NAME

Section 1: Enzymatic DNA Fragmentation (Manually)

VERSION 1

CREATED BY




Ester Kalef-EzraUniversity College London

OPEN

→

Protocol




NAME

Section 2: NGS library preparation for sequencing

VERSION 1

CREATED BY




Ester Kalef-EzraUniversity College London

OPEN

→

Protocol




NAME

Section 3: Libraries quality control (QC)

VERSION 1

CREATED BY



Ester Kalef-EzraUniversity College London

OPEN

→