



SEP 01, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.14egn36mml5d/v1

Protocol Citation: Daniel Manrique-Castano 2023. Staining of GFAP and IBA1 in PDGFR-B/Td-tomato brain sections. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.14egn36mml5d/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working

Created: Sep 01, 2023

Staining of GFAP and IBA1 in PDGFR-B/Td-tomato brain sections

Daniel Manrique-Castano¹

¹Université Laval

Protocols



Daniel Manrique-Castano
 Université Laval

ABSTRACT

This protocol is suitable for the staining of GFAP and IBA1 in PDGFR-B/Td-Tomato fixed mouse brain sections.

GUIDELINES

Read the entire protocol before starting the procedure.

Note that this protocol uses 3 hours (room temperature) incubation for primary antibodies.

Do not submit the tissue to antigen retrieval or other acidic solutions to preserve the Td-tomato signal.

Works with the listed antibodies and dilutions. For other references, preliminary tests are highly recommended.

MATERIALS


Last Modified: Sep 01, 2023

PROTOCOL integer ID: 87251


Keywords: Staining, GFAP, IBA1, Immunofluorescence

Equipment	
ImmEdge® Hydrophobic Barrier PAP Pen	NAME
Vector	BRAND
H-4000	SKU
https://vectorlabs.com/products/histology/immedge-hydrophobic-barrier-pen	LINK
Immunohistochemistry / Immunocytochemistry, Immunofluorescence, In situ hybridization	SPECIFICATIONS


Permeabilization solution = [IM] 0.3 % (v/v) Triton


 Triton X-100 Sigma
Aldrich Catalog #T8787-50ML

/ = [IM] 0.1 % (v/v)


 TWEEN 20 Sigma
Aldrich Catalog #P7949

/


[IM] 0.3 Molarity (m)  Glycine Contributed by users in 1x PBS

Antibody buffer = [IM] 0.1 % (v/v)  TWEEN 20 Sigma
Aldrich Catalog #P7949


/ [IM] 1 % (v/v)


 Normal Donkey Serum Merck MilliporeSigma (Sigma-Aldrich) Catalog #D9663-10ml

[IM] 1 Mass / % volume


 Bovine Serum Albumin (BSA) Sigma
Aldrich Catalog #A7906


Primary antibodies:


 GFAP Monoclonal Antibody (2.2B10) Thermo Fisher Scientific Catalog #13-0300

 Anti Iba1 Rabbit antibody FUJIFILM Wako Pure Chemical Corporation Catalog #019-19741


Secondary antibodies:

 Alexa Fluor® 647 AffiniPure Donkey Anti-Rat IgG (H L) Jackson ImmunoResearch Laboratories, Inc. Catalog #712-605-153


 Anti-Rabbit Invitrogen - Thermo Fisher Catalog #A21206

 DAPI (46-Diamidino-2-Phenylindole Dilactate) Invitrogen - Thermo Fisher Catalog #D3571

Mounting media:


 Fluoromount-G Electron Microscopy Sciences Catalog #17984-25

PROTOCOL MATERIALS


 Alexa Fluor® 647 AffiniPure Donkey Anti-Rat IgG (H L) Jackson ImmunoResearch Laboratories, Inc. Catalog #712-605-153

Materials

 Glycine **Contributed by users** Materials


 TWEEN 20 Merck MilliporeSigma (Sigma-Aldrich) Catalog #P7949

In Materials, Materials and [3 steps](#)


 Fluoromount-G Electron Microscopy Sciences Catalog #17984-25

Materials, Step 9


 Anti-Rabbit Invitrogen - Thermo Fisher Catalog #A21206 Materials

 Triton X-100 Merck MilliporeSigma (Sigma-Aldrich) Catalog #T8787-50ML

Materials

 Bovine Serum Albumin (BSA) Merck MilliporeSigma (Sigma-Aldrich) Catalog #A7906

Materials

 Normal Donkey Serum Merck MilliporeSigma (Sigma-Aldrich) Catalog #D9663-10ml

Materials



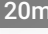
SAFETY WARNINGS


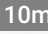


DAPI is highly toxic. Handle it with care.

Tissue preparation and blocking


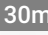
20m

- 1 Take out sections from -80 and place them in an incubator/plate for  00:20:00 at  37 °C  20m
. This step is performed to ensure tissue attachment to the crystal slides.

- 2 Draw a hydrophobic barrier on each slide using **ImmEdge® Hydrophobic Barrier PAP Pen** and let dry for about  00:10:00 minutes. This will prevent buffers from escaping the tissue area.  10m

Note


There are other hydrophobic pens on the market. However, this one is recommended for its quality, consistency, and durability.

- 3 To initially rehydrate and permeabilize the tissue, place the slides in a slide jar containing the **permeabilization solution** with shaking for  00:30:00  30m


Note


The proposed permeabilization solution contains **glycine**, which is suitable for reducing PFA-related autofluorescence, especially in the 488 channel.

Antigen retrieval or any treatment with acid is not recommended given that it destroys the endogenous Td-tomato signal from PDGFR-B+ cells.


- 4 To prevent unspecific binding, incubate brain sections in **permeabilization solution** containing  1h



 5 % (v/v)

 Normal Donkey Serum Merck MilliporeSigma (Sigma-Aldrich) Catalog #D9663-10ml and

 1 Mass / % volume

 Bovine Serum Albumin (BSA) Sigma Aldrich Catalog #A7906 for

 01:00:00 at  Room temperature .



Note

Please note that blocking serum must be chosen according to **secondary antibody species**. For the procedure depicted in this protocol, all secondary antibodies are raised in donkey. However, goat secondary antibodies will also be suitable.

Additionally, consider **avoiding BSA** when **primary antibodies** come from goat or sheep. BSA can generate unspecific background.

Antibody incubation

3h

- 5 When blocking is finished, decant the buffer (no washing is required) and incubate **primary antibodies** in **antibody buffer** for  03:00:00 at  Room temperature according to **table 1**.

3h



A	B	C	D	E
Antibody	Company	Reference	Specie	Dilution
Gfap	Invitrogen	13-0300	Rat	1:500
Iba1	Wako	019-19741	Rabbit	1:500

Table 1. Primary antibodies

Note

Please note that antibody buffers **do not contain Triton X-100**. This detergent tends to break the hydrophobic barrier and the sections may not be adequately incubated. In addition, given that the sections are already permeabilized, we do not contemplate the use of this detergent.

Please note that the proposed antibody incubation lasts only **3 hours**. Although in most cases staining protocols recommend overnight incubation for primary antibodies, the tests performed in our lab showed better efficiency at room temperature.

The indicated **Iba1 antibody** is superior to others available in the market.


6

When primary antibody incubation is finished, wash the sections  00:05:00 5x with

5m



[M] 0.05 % (v/v)

 TWEEN 20 Sigma
Aldrich Catalog #P7949


in PBS.

7

Incubate **secondary antibodies** in [M] 0.1 % (v/v)

1h



 TWEEN 20 Sigma
Aldrich Catalog #P7949

for  01:00:00 at  Room temperature




according to **table 2**.

A	B	C	D	E
Antibody	Company	Reference	Channel	Dilution
Donkey anti-rat	Jackson	712-605-153	647	1:500
Donkey anti-rabbit	Invitrogen	A21206	488	1:500
Dapi	Invitrogen	D3571	405 (Dapi)	1:5000

Table 2. Secondary antibodies


8

When secondary antibody incubation is finished, wash the sections  00:05:00 3x with


10m



[M] 0.05 % (v/v)


 TWEEN 20 Sigma
Aldrich Catalog #P7949

in PBS. Follow this with

 00:05:00 2x

washes with **PBS** to remove all detergent traces.

9

Clean the remaining buffer on the slides using absorbent tissue and mount the sections with a drop of  Fluoromount-G Electron Microscopy
Sciences Catalog #17984-25