

Aug 30, 2024

## JAX-Sen: Mouse pancreas dissociation for single-cell RNA sequencing

DOI

[dx.doi.org/10.17504/protocols.io.ewov19zzylr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov19zzylr2/v1)

Ramalakshmi Ramasamy<sup>1</sup>, Juliana Alcoforado Diniz<sup>1</sup>, Patrick Fleming<sup>1</sup>, Jessica Garofalo<sup>1</sup>, Dylan Baker<sup>1</sup>, Paul Robson<sup>1,2</sup>

<sup>1</sup>The Jackson Laboratory for Genomic Medicine, Farmington, CT, USA;

<sup>2</sup>Department of Genetics and Genome Sciences, University of Connecticut School of Medicine, Farmington, CT, USA

Cellular Senescence Net...



Harshpreet Chandok

The Jackson Laboratory

OPEN  ACCESS



DOI: [dx.doi.org/10.17504/protocols.io.ewov19zzylr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov19zzylr2/v1)

**Protocol Citation:** Ramalakshmi Ramasamy, Juliana Alcoforado Diniz, Patrick Fleming, Jessica Garofalo, Dylan Baker, Paul Robson 2024. JAX-Sen: Mouse pancreas dissociation for single-cell RNA sequencing. **protocols.io**

<https://dx.doi.org/10.17504/protocols.io.ewov19zzylr2/v1>

**License:** This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**We use this protocol and it's working**

**Created:** July 02, 2024

**Last Modified:** August 30, 2024

**Protocol Integer ID:** 102769

**Keywords:** JAX-Sen



**Funders Acknowledgement:**

**National Institute on Aging  
(NIA) JAX-Sen Senescence  
Tissue Mapping Center Grant  
Grant ID: ID: U54 AG079753**

## Abstract

These samples are part of the JAX-Sen project in the SenNeT Consortium. We aim to study and characterize senescence in the C57Bl/6 mouse pancreas. We compare young (6 months old) and old (24 months old) mouse pancreas using scRNAsequencing. This protocol describes the single-cell dissociation of the pancreas before library preparation and sequencing.



## Reagents and Materials

- 1
  - 15 ml conical tubes – 4 per sample
  - 1.5 or 2 ml Protein lobind tubes
  - Micro scissors/ scalpel
  - Ice-cold PBS
  - Histopaque-1077 – Sigma H8889
  - Wide-bore pipette tips
  - Collagenase solution
  - Wash buffer
  - 100 µm cell strainer (Corning, 431752)
  - ACK lysis buffer – 2 ml/sample
  - SS\_buffer – 5 ml/sample
  - 40 um cell strainer
  - 70 uM strainer
  - 10 ml syringe and 18G needle if available

### 2 **SS\_buffer:**

| A    | B               | C                | D     | E      | F      | G      | H     | I  |
|------|-----------------|------------------|-------|--------|--------|--------|-------|----|
|      | Stock co<br>nc. | Working<br>conc. | 50 ml | 100 ml | 250 ml | 500 ml | 20 ml |    |
| EDTA | 0.5 M           | 2 mM             | 0.2   | 0.4    | 1      | 2      | 0.08  | ml |
| BSA  | 10%             | 2%               | 10    | 20     | 50     | 100    | 4     | ml |
| FBS  | 100%            | 15%              | 7.5   | 15     | 37.5   | 75     | 3     | ml |
| DPBS | 1X              | make up          | 32.2  | 64.4   | 161    | 322    | 12.88 | ml |

#### 2.1 **Collagenase solution:**

For each mouse, dissolve 4mg collagenase powder in 10mL serum-free Medium 199 (or MEM). Prepare fresh daily. Fill 5mL syringe with 3mL collagenase P (Roche Applied Science 11 249 002 001) and 15mL tube with 5mL collagenase; label tubes and store on ice until ready to use. Note: each lot of collagenase must be tested experimentally to determine optimal working concentration. Use 0.4mg/mL as starting point.

- **Wash buffer:**

DMEM containing 20% FBS.

- **Stopping media:** DMEM + 10% FBS + 1:100 dilution Glutamate CMRL media




## Procedure

- 2.2 Collect shipment from the warehouse.
- 2.3 Fill up ice buckets and prep them with reagents:  
Bucket 1 (small): PBS, HBSS, DMEM+FBS, SS\_buffer  
Bucket 2 (large): For mincing tissue
- 2.4 Confirm if the sample is cold. Take pictures.
- 2.5 Place tubes on ice and mince the tissue (<1mm) with sterile micro scissors.
- 2.6 After mincing, transfer the sample into labeled 15 ml tubes containing **5 ml of 0.4mg/ml collagenase solution**. Incubate at 37°C for 10 mins.
- 2.7 Remove and shake tubes vigorously until pancreatic tissue is completely dissociated with no large chunks visible, 10-20 seconds. Alternatively, triturate using 10 mL pipet until well-dispersed.
- 2.8 **Immediately** dilute digest with cold wash buffer to 15mL total volume. Invert several times to mix. **Keep all the tubes cold.**
- 2.9 Centrifuge at **300 x g, 4°C, 1 min** to pellet pancreatic tissue.
- 2.10 Decant supernatant into waste beaker. Add **12mL wash buffer** and resuspend pellet by shaking or with 10mL pipette.
- 2.11 Repeat wash (steps #2.9-2.10) **twice**. After the final spin, decant supernatant, then invert open tubes onto a paper towel and tap gently to remove all residual buffer.
- 2.12 Add **5mL Histopaque-1077** and gently resuspend pellet using 10mL pipette.
- 2.13 Hold tube at 45° angle and *slowly* layer **5mL wash buffer** on top of gradient layer (**10mL syringe** is helpful for slowly dispensing buffer). Do not allow buffer to mix with Histopaque! Very

carefully transfer tubes to centrifuge.

2.14 **Turn off centrifuge brake.** Centrifuge at 300xg/4°C for 25 minutes.

2.15 *Islets should be visible floating at interface between Histopaque and wash buffer, or slightly below.* Remove islets from interface using sterile transfer pipette and transfer to new 15mL conical tube. Centrifuge at 300xg/4°C for 25 minutes. *Other cells (mostly acinar): in the pellet.*

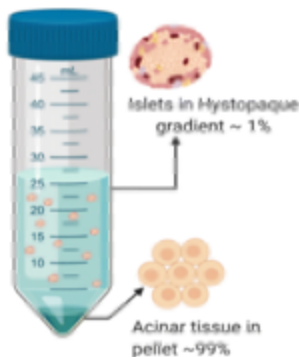

Split the 2 layers and treat them as separate samples going forward.

***Proceed to islets protocol and acinar+ducts protocols that we use for human pancreatic cells.***

## Acinars and Ducts:

2.16 Wash pellets from histopaq gradient **3x** with 10 ml cold PBS

2.17 Add 3 ml TrypLE E to the pellet and incubate at 37C (waterbath), mixing with pipette approximately every 5 mins. Check at 10 mins and every 5 mins following for single cells. (Ducts should take 10-15 mins, and acinar should take slightly longer ~20 mins. But we have a mixture here, so check under the microscope and stop reaction accordingly)

2.18 Add **3X stopping media (9 ml).**

2.19 Centrifuge at 1300 rpm/4°C for 2 minutes. Discard supernatant.



2.20 Wash again with **3 ml stopping media** and pellet cells. Resuspend in 1 ml of stop media.

2.21 Filter through 40  $\mu$ M Flowmi filter.

2.22 Count cells and proceed to Flex protocol.

### Islets:

2.23 Centrifuge the islets suspension at 130g/RT for 5 minutes.

2.24 Add 1.5 ml accutase, mix with pipette and incubate at 37C. Mix every 2 mins, checking at 6 mins.

2.25 Add CMRL media to approximately 15 ml total volume.

2.26 Centrifuge at 500xg/4°C for 5 mins.

2.27 Resuspend in 1 ml of CMRI media.

2.28 Filter through 40  $\mu$ M filter.

2.29 Count cells and proceed to Flex protocol.