

Apr 28, 2024 Version 4

SARS-CoV-2 Mpro fluorescence dose response V.4

DOI

dx.doi.org/10.17504/protocols.io.81wgbye9nvpk/v4

Haim Barr^{1,2}, Noa Lahav^{1,2}

¹Weizmann Institute of Science; ²ASAP Drug Discovery Consortium

Haim Barr: General acknowledgement: The Wohl Drug Discovery institute, The Nancy and Stephen Grand Israel National Center for Personalized Medicine.;

ASAP Discovery



Mary-Ann Xavier

Diamond Light Source

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.81wgbye9nvpk/v4

Protocol Citation: Haim Barr, Noa Lahav 2024. SARS-CoV-2 Mpro fluorescence dose response. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.81wgbye9nvpk/v4> Version created by **[Mary-Ann Xavier](#)**

License: This is an open access protocol distributed under the terms of the **[Creative Commons Attribution License](#)**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working

Created: April 12, 2024

Last Modified: April 28, 2024

Protocol Integer ID: 98134

Keywords: SARS-CoV-2, Main Protease, dose response, Coronavirus

Funders Acknowledgement:

**National Institutes of
Health/National Institute Of
Allergy and Infectious
Diseases (NIH/NIAID)**
Grant ID: U19AI171399

Disclaimer

The content is solely the responsibility of the authors and does not necessarily represent the official views of the National Institutes of Health.

Abstract

This is a **functional, biochemical assay** used to identify treatments for viral infectious disease that target SARS-COV-2 Main Protease (MPro).

Utilizing a **direct enzyme activity measurement method**, the experiment was performed in a 384-well plate reading the fluorescence intensity. This assay tested the mode of action of inhibition.

It was developed at the Weizmann Institute of Science, as a part of the ASAP Drug Discovery Consortium.

Experiment Assay Concentrations

A	B	C
Reagent	Final Assay Concentration	Units
SARS Mpro	5	nM
SARS Substrate	375	nM
HEPES (pH 7.3)	20	mM
NaCl	50	mM
Glycerol	10	% by volume
TWEEN 20	0.01	% by volume
TCEP	1	mM

For more information, please check out the "Materials" Section



Guidelines

Plate Information:

Total Assay Volume: 20 μ L

Compounds Top Assay Concentration: 100 μ M

Dilution Factor: 2

Dose Response Points: 12

Number of Replicates: 2

Backfill with DMSO: Yes

Materials

Assay Buffer Reagents (Concentration listed are Stock Solution Concentrations)

1. [M] 40 millimolar (mM) ☒ HEPES 1M Solution pH 7.3 Fisher Scientific Catalog #AAJ16924K2 (or similar)
2. [M] 100 millimolar (mM) ☒ Sodium chloride Merck MilliporeSigma (Sigma-Aldrich) Catalog #S9888-25G (or similar)
3. [M] 50 % volume ☒ Glycerol Merck MilliporeSigma (Sigma-Aldrich) Catalog #G5516 (or similar)
4. [M] 10 % volume ☒ TWEEN® 20 Merck MilliporeSigma (Sigma-Aldrich) Catalog #P9416 (or similar)
5. [M] 1000 millimolar (mM) ☒ Tris(2-carboxyethyl)phosphine hydrochloride Merck MilliporeSigma (Sigma-Aldrich) Catalog #75259 (TCEP) (or similar)

***Note:** all components are added fresh to the assay buffer before each experiment

Additional Reagents:

[M] 710 micromolar (μ M) **SARS MPro Enzyme**

- The Enzyme original stock was originally [M] 750 micromolar (μ M) and was diluted to create aliquots of [M] 20000 nanomolar (nM) using a **storage buffer** (50 mM Tris pH 7.5, 1 mM DTT, 50 mM NaCl, 1 mM EDTA, 50% Glycerol).
- Before an experiment, the 20000 nM aliquots were **diluted with Assay Buffer** to create a [M] 10 nanomolar (nM) solution to be loaded into the Combi.

[M] 20000 micromolar (μ M) **SARS MPro Substrate**

- SARS MPro Substrate Stock ([5-FAM]-AVLQSGFR-[Lys(DabcyI)-K-amide) was purchased and dissolved in **DMSO** and yielded a concentration of [M] 20000 micromolar (μ M)
- Before an experiment, the SARS MPro Substrate Stock had an *intermediate dilution step* with **DMSO** to yield a [M] 100 micromolar (μ M) SARS MPro Substrate Solution. Then, before adding the SARS MPro Substrate to the Combi, it was diluted again with **Assay Buffer** to yield a concentration of [M] 750 nanomolar (nM). The final concentration of **SARS MPro Substrate** for the assay was [M] 375 nanomolar (nM)



Safety warnings


⚠ Please be sure to wear proper Personal Protective Equipment (PPE) while performing this experiment.

Before start

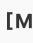

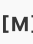
Note: Inhibitor compounds stock concentration is **20 mM**. Compounds are pre-dispensed into 384 plates and stored at -200°C until use.



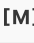

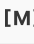



Prepare 384 Well Plate

- 1 **PRIME** Multi-Drop Combi Tube Dispensing Cassette with **Assay Buffer** by selecting the **PRIME** button on the Combi Dispenser until the tubes are filled completely.
- 1.1 **DISPENSE**  10 µL Assay Buffer to Columns **1** and **23** of assay plate
 - **Note:** These will represent the ***inhibitor control columns*** (Contain: Substrate, Assay Buffer, DMSO; **no experimental compounds**)
- 1.2 **EMPTY** Multi-Drop Combi Tube Dispensing Cassette (by selecting the **EMPTY** button on the Combi Dispenser until the tubes of the cassette are emptied). Discard the Assay Buffer discharged from the cassette.

Prepare Reagents

- 2 **PRIME** Multi-Drop Combi Tube Dispensing Cassette with  10 nanomolar (nM) SARS MPro by selecting the **PRIME** button on the Combi Dispenser until the tubes were filled completely.
 - **Note:** Be sure to cycle dispensing several times on a clean plate lid (This confirms there are no bubbles in the Dispensing Cassette).
- 2.1 **DISPENSE**  10 µL  10 nanomolar (nM) SARS MPro to Columns **2** through **22** and also Column **24**.

Note:


 -  10 nanomolar (nM) SARS Mpro is two times the final concentration for the assay. It is diluted to be a final concentration of  5 nanomolar (nM) SARS MPro .
 - Column 2 and Column 24 are ***neutral control columns*** (Contain: Enzyme, Substrate, DMSO; **no experimental compounds**)
- 2.2 **EMPTY** Multi-Drop Combi Tube Dispensing Cassette (by selecting the **EMPTY** button on the Combi Dispenser until the tubes of the cassette are emptied). Discard the  10 nanomolar (nM) SARS MPro discharged from the cassette.
- 3 **CENTRIFUGE**  15000 rpm, Room temperature, 00:01:00 plate to remove bubbles 1m
- 4 **INCUBATE** plate for  00:15:00 at  Room temperature 15m



5 **PRIME** Multi-Drop Combi Tube Dispensing Cassette with **Assay Buffer** by selecting the **PRIME** button on the Combi Dispenser until the tubes are filled completely. Then, **EMPTY** the Multi-Drop Combi Tube Dispensing Cassette (by selecting the **EMPTY** button on the Combi Dispenser until the tubes of the cassette are emptied). **Discard the Assay Buffer discharged from the cassette.**


6 **PRIME** Multi-Drop Combi Tube Dispensing Cassette with [M] 750 nanomolar (nM) SARS Substrate by selecting the **PRIME** button on the Combi Dispenser until the tubes were filled completely.

- **Note:** Be sure to cycle dispensing several times on a clean plate lid (This confirms there are no bubbles in the Dispensing Cassette).



6.1 **DISPENSE**  10 µL [M] 750 nanomolar (nM) SARS Substrate into Columns **1 through 24** (the full plate)

Note:

- [M] 750 nanomolar (nM) SARS Substrate is two times the final concentration for the assay. It is diluted to be a final concentration of [M] 375 nanomolar (nM) SARS Substrate

7 **CENTRIFUGE** plate at  15000 rpm, Room temperature, 00:01:00 in plate centrifuge to remove bubbles

1m

8 **INCUBATE** plate at  Room temperature for  00:30:00

30m

⚠ **Make sure the plate is protected from light!**

Recommended: Clean/Empty the Multi-Drop Combi Reagent Dispenser and Dispensing Cassette during this incubation step

Read Plate Fluorescence

9 **READ** and **RECORD** the plate Relative fluorescence units (RFU) via the **"SARS Endpoint protocol"** on the **PERAstar FS Control Software**.

Expected result

Gain 300 should yield ~10,000 RFU in full reaction and ~6,000 RFU in Buffer Control