




May 28, 2021

Glyoxal fixation of mammalian cells for immunofluorescence

Karla LH LH Feijs¹¹Institute of Biochemistry and Molecular Biology, RWTH Aachen University Hospital

1 Works for me

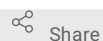

Share

dx.doi.org/10.17504/protocols.io.bvcyn2xw

FeZa lab

Karla LH Feijs
Institute of Biochemistry and Molecular Biology, RWTH Aachen...

ABSTRACT

This is our standard protocol for glyoxal fixation of mammalian cells grown on glass cover slips for immunofluorescence microscopy.

DOI

dx.doi.org/10.17504/protocols.io.bvcyn2xw

PROTOCOL CITATION

Karla LH LH Feijs 2021. Glyoxal fixation of mammalian cells for immunofluorescence . **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bvcyn2xw>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 28, 2021

LAST MODIFIED

May 28, 2021

PROTOCOL INTEGER ID

50296

MATERIALS TEXT

Standard materials to make the required buffers:

ddH₂O
PBS
Ethanol
Glyoxal, 40%
Glacial acetic acid
BSA
Ammonium chloride

Mountant containing DAPI

All solutions can be made up to one day in advance and stored in the fridge.

Glyoxal solution

- 2.835 ml ddH₂O
- 0.789 ml ethanol (absolute)
- 0.313 ml glyoxal (40% stock solution)
- 0.03 ml glacial acetic acid (100%)

Adjust the pH to 5-6 with 1M NaOH using pH paper

Quenching solution

- 0.535 mg NH₄Cl

Adjust to 100ml with PBS

Blocking solution (PBS with 2.5% BSA and 0.1% Triton X-100)

- 2.5g BSA
- 1 ml Triton X-100 (10% stock solution)

Adjust to 100 ml with PBS

Antibody buffer (PBS with 1% BSA and 0.1% Triton X-100)

- 1g BSA
- 1 ml Triton X-100 (10% stock solution)

Adjust to 100 ml with PBS

Mountant with DAPI

Fix the cells 5m

- 1 Wash the cells with warm DMEM without FCS. 5m
- 2 Add 300µl **glyoxal solution** per well (12-well plate) and incubate on ice for 30 minutes ⚡ **On ice** ⌚ **00:30:00** ,
then at room temperature for 20 minutes ⚡ **Room temperature** ⌚ **00:20:00** . 40m
- 3 Wash 2 times with PBS 10m

4 Add 300µl **quenching solution** for 20 minutes 📍 **Room temperature** ⌚ 00:20:00 20m

5 Wash once with PBS 5m

At this stage, fixed cells can be stored in the fridge for at least a week.

Staining the cells 2h 35m

6 Add 500µl **blocking solution** and incubate at room temperature for 30 minutes 📍 **Room temperature** ⌚ 00:30:00 30m

In meantime, prepare the primary antibody. The dilution is dependent on the antibody, but can range from undiluted for non-purified hybridoma supernatants to 1:1000 for affinity purified antibodies. The antibodies are diluted in **antibody buffer**.

7 There are two ways to incubate the antibodies: 45m
1. Prepare enough so that the slide can be covered in the well with the antibody solution, ±250-300µl.
2. To save antibody, prepare only 50µl and pipette this onto parafilm. Place the coverslip on top of this drop *with the cells facing the drop* and incubate in a humid chamber.

Incubate at least 45 minutes at room temperature or alternatively overnight at 4°C 📍 **Room temperature** ⌚ 00:45:00

8 Place the coverslips back in the 12-well plate with *the cells facing upwards*. 15m
Wash 3 times with PBS.

9 Dilute the secondary antibody in **antibody buffer**, generally in a range of 1:1000-1:2000 for most antibodies. Add 250-300µl per well and incubate for at least 45 minutes at room temperature or alternatively overnight at 4°C 45m
📍 **Room temperature** ⌚ 00:45:00 .

IMPORTANT: this incubation step should be in the dark, as the conjugated secondary antibodies are light-sensitive.

10 Wash 2 times in PBS, followed by 2 times in ddH2O. 15m

11 Place a drop of mountant containing DAPI on a microscope slide and place the coverslip on this drop with the *cells facing the drop*. 5m

IMPORTANT: warm up the mountant to 37°C to avoid unequal DAPI-staining. Tip the coverslip sideways on a tissue to remove excess water before mounting.

12 Dry the slides overnight in a dark place at room temperature. Store at 4°C afterwards until imaging.