




Sep 21, 2021

Cloning into pSL2680 CRISPR Plasmid - iGEM IISER Pune 2021

misaal.bedi¹¹Indian Institute of Science Education and Research, Pune

In Development

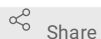

Share

dx.doi.org/10.17504/protocols.io.bw9kph4w

iGEM IISER Pune India 2021

misaal.bedi

ABSTRACT

This protocol can be used to clone the guide RNA and Homologous Repair Template into the pSL2680 plasmid which was a gift from Himadri Pakrasi (Addgene plasmid # 85581 ; <http://n2t.net/addgene:85581> ; RRID:Addgene_85581).

DOI

dx.doi.org/10.17504/protocols.io.bw9kph4w

PROTOCOL CITATION

misaal.bedi 2021. Cloning into pSL2680 CRISPR Plasmid - iGEM IISER Pune 2021. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bw9kph4w>

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Cpf1 Is A Versatile Tool for CRISPR Genome Editing Across Diverse Species of Cyanobacteria. Ungerer J, Pakrasi HB. Sci Rep. 2016 Dec 21;6:39681. doi: 10.1038/srep39681. 10.1038/srep39681 PubMed 28000776

LICENSE

— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 10, 2021

LAST MODIFIED

Sep 21, 2021

PROTOCOL INTEGER ID

52236

MATERIALS TEXT

pSL2680 was a gift from Himadri Pakrasi (Addgene plasmid # 85581 ; <http://n2t.net/addgene:85581> ; RRID:Addgene_85581)

Restriction enzymes: AarI, KpnI + their buffers

Antibiotics: Kanamycin

LB media

Distilled Water

Plasmid Miniprep Kit

Gel Purification Kit

Gel electrophoresis set up

Ligase buffer

PNK Enzyme

NEB Gibson Assembly Kit

Cloning the gRNA

- 1 The gRNA is to be cloned before the cassette into the plasmid.

Day 1

- 2 Take **30 mL** of LB media with **0.05 µg/µl** Kanamycin and scratch off some solid culture from the vial and add to the media in a flask. Grow this **Overnight** at **37 °C**

Day 2

- 3 Perform miniprep on the entire **30 mL** culture (since the pSL2680 plasmid is of a low copy number). Spin it down and resuspend it in **500 µl** resuspension buffer in the fridge. Divide it into two sets and follow the rest of the miniprep protocol. In the end, elute from the column with **50 µl** of distilled water.

- 3.1 Digest the **50 µl** of miniprep with AarI using the following reaction mix:

AarI buffer - **10 µl**

plasmid prep - **50 µl**

50X AarI oligo - **2 µl**

AarI enzyme - **4 µl**

Distilled water - **34 µl**

Incubate at **37 °C** for 4 hours.

- 3.2 Gel extract the digest. Run the entire digest in a double or triple lane on a 0.7-1% agarose gel. Excise the largest band and place it in a **2 mL** Eppendorf tube. Melt at **55 °C** (Tip: Inverting the tube every five minutes helps in faster melting)

3.3 Perform gel purification using a gel purification kit and follow its protocol.

3.4 In the end, elute with $\square 20 \mu\text{l}$ of water (ensure that a small amount of water is placed in the center of the membrane).

Day 3

4 Anneal the gRNA oligos with the following reaction mixture:

$\square 100 \text{ Micromolar } (\mu\text{M})$ stock of the gRNA left primer - $\square 10 \mu\text{l}$

$\square 100 \text{ Micromolar } (\mu\text{M})$ stock of the gRNA right primer - $\square 10 \mu\text{l}$

Ligase buffer - $\square 5 \mu\text{l}$

Distilled water - $\square 25 \mu\text{l}$

PNK enzyme - $\square 1 \mu\text{l}$

Use the following thermocycler program:

Heat to $\delta 95^\circ\text{C}$ for 5 minutes, and then ramp to $\delta 4^\circ\text{C}$ at $\delta 0.1^\circ\text{C}$ /second.

4.1 Dilute the annealed oligos 1:50. Not a lot is required so $\square 2 \mu\text{l}$ with $\square 98 \mu\text{l}$ of water can be used.

Ligate the oligos into pSL2680 which was digested the previous day with the following reaction mix:

pSL2680 gel extract - $\square 8 \mu\text{l}$

Annealed oligos that were dilute 1:50 - $\square 0.5 \mu\text{l}$

Ligase buffer - $\square 1 \mu\text{l}$

Ligase - $\square 0.5 \mu\text{l}$

Incubate \odot **Overnight** at $\delta 16^\circ\text{C}$

Day 4

5 Transform the entire $\square 10 \mu\text{l}$ ligation reaction into the E. coli XL1 - blue or Top10 strains as follows:

Thaw cells δ **On ice** for 10 minutes.

Add $\square 10 \mu\text{l}$ ligation reaction to cells and stir. Don't pipette up and down.

Incubate δ **On ice** for another 10 minutes.

Heat shock for 45 seconds at $\delta 42^\circ\text{C}$.

Place δ **On ice** for 2 minutes.

Add $\square 500 \mu\text{l}$ of liquid LB media.

Incubate at 37°C for 1 hour.

- 5.1 Centrifuge the reaction mixture after incubation and remove the supernatant. Add the pellet to $50\ \mu\text{l}$ of the media.
- Plate transformants on LB media with $0.05\ \mu\text{g}/\mu\text{l}$ Kanamycin and add IPTG+X-gal.
- Spread $40\ \mu\text{l}$ of $20\ \text{mM}$ IPTG and $20\ \text{mg}/\text{ml}$ X-gal on a plate with LB media with $0.05\ \mu\text{g}/\mu\text{l}$ Kanamycin and let it soak in for at least an hour.
- Prepare the plates before starting the transformation.

Day 5

- 6 Pick 8 white colonies from the transformants and patch them on LB media with $0.05\ \mu\text{g}/\mu\text{l}$ Kanamycin.
- Start $10\ \text{mL}$ cultures of 2 colonies off the plate in LB media with $0.05\ \mu\text{g}/\mu\text{l}$ Kanamycin.

Day 6

- 7 Isolate the plasmid from the 2 cultures incubated the previous day. Determine the concentration and send for sequencing. Once the correct sequence is confirmed, prepare stocks.

Cloning of cassette

- 8 After the gRNA has been cloned into the plasmid, we can begin cloning the cassette to make modifications in the *S. elongatus* UTEX 2973 genome.

Day 1

- 9 Synthesize the repair template using high-fidelity PCR.

Gel extract the PCR reactions and purify them, and elute in $25\ \mu\text{l}$ of distilled water.

Digest the plasmid carrying the gRNA using KpnI and make sure to dephosphorylate it, to prevent recircularization of the backbone.

Use the following reaction mixture:

Buffer - $15\ \mu\text{l}$

Plasmid Prep - $75\ \mu\text{l}$

KpnI restriction enzyme - $3\ \mu\text{l}$

FastAP (dephosphorylates vector) -

Distilled Water - $53\ \mu\text{l}$

Incubate for 4 hours at δ 37 °C

Gel extract the digest and elute in \square 25 μ l of distilled water.

Day 2

- 10 Concentrate the gel extractions and perform Gibson Assembly and concentrate the DNA. Allow samples to dry for 2.5-3 hours.

For Gibson Assembly, exactly \square 5 μ l of correctly mixed and highly concentrated DNA is required. Use \square 200 ng of the vector and a 2X molar ratio of each fragment to be assembled by Gibson Assembly.

Note: Use molar ratios of the fragments, not weight.

The vector and PCR fragments should be mixed using the

Prepare the mixture in a PCR tube with a total volume of \square 5 μ l .

Add \square 15 μ l of the Gibson master mix to the \square 5 μ l of DNA.

Mix them together and PCR at δ 50 °C for an hour and then hold at δ 4 °C

Day 3

- 11 Transform at least \square 10 μ l of the Gibson reaction in XL1-blue (or Top10) using the same protocol as described in step 5 of Day 4 in the 'Cloning the gRNA' section of this protocol.
Plate on LB media with \square 0.05 μ g/ μ l Kanamycin

Day 4

- 12 Check the transformed colonies by PCR.
Run the reactions on a gel and select two positive colonies.

Start with \square 10 mL of each colony in LB media with \square 0.05 μ g/ μ l Kanamycin for plasmid isolation.

Day 5

- 13 Isolate plasmids from the two colonies isolated the previous day and elute in \square 50 μ l .
Sequence the insert in the plasmid. If the sequencing is correct, this protocol has been successfully completed and all that remains is to transform the plasmid into cyanobacteria.

Note: Don't forget to freeze the recombinant E.coli in permanent stock.