


Jul 17, 2022

ONT DNA Barcoding Fungal Amplicons w/ MinION & Flongle

Stephen Douglas Russell¹¹The Hoosier Mushroom Society

1 Works for me

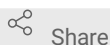

Share

dx.doi.org/10.17504/protocols.io.36wgq7qykvk5/v1

The Hoosier Mushroom Society



Stephen Douglas Russell

ABSTRACT

This collection of protocols outline a working methodology to DNA barcode fungal specimens. This process will work for dried tissue, fresh tissue, or with DNA template that has already gone through an alternative extraction protocol. It takes you from computer/software setup through DNA extraction, amplification, library preparation, sequencing, and primary data analysis.

This particular workflow is designed for use with **Flongle 9.4.1 flowcells** and **V12 Ligation** chemistry.

DOI

dx.doi.org/10.17504/protocols.io.36wgq7qykvk5/v1

COLLECTION CITATION

Stephen Douglas Russell 2022. ONT DNA Barcoding Fungal Amplicons w/ MinION & Flongle. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.36wgq7qykvk5/v1>



LICENSE

This is an open access collection distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jun 09, 2022

LAST MODIFIED

Jul 17, 2022

COLLECTION INTEGER ID

64290

MATERIALS TEXT

Total equipment needs and cost:

PCR tube rack x10 ([Amazon](#)) : \$33.78 (would be good to have at least 10)
Tweezers (flat, non-serrated tips - Ebay): \$10.00
Teasing needle ([Amazon](#)): \$12.35
Mini centrifuge ([Ebay](#)): \$70.00
PCR tube rack 1.5mL ([Amazon](#)): \$26.42
0.5 -10uL multichannel pipette ([Amazon](#)): \$155.51 (may be able to find used cheaper on Ebay)
50 - 300uL multichannel pipette ([Amazon](#)): \$172.89 (may be able to find used cheaper on Ebay)
Summit Professional Freezer -20C (Facebook Marketplace / [1](#) / [2](#)): \$150 (gotta get lucky)
Magnetic bead separator for 1.5mL eppi tubes ([Ebay](#)): \$59.00
Tip disposal bucket
Gel electrophoresis system ([miniPCR](#)): \$300
Heat block ([Amazon](#)): \$179.99
Hula mixer (Ebay): \$200.00 (optional)
Flongle Starter Pack: \$1,460.00
MinION Mk1B Starter Pack: \$1000.00
Thermocycler: \$100 - \$1600

Total upfront equipment cost: ~\$3,200 - \$3,700, excluding thermocycler and computer
Eppendorf Mastercycler (used) - \$1600 (can find less expensive models used, reasonable units in the \$200-\$400 range)
System 76 laptop - \$4,000 (can probably find something that works for much cheaper...depending on your use case)

Consumables:

0.2mL PCR tubes ([Amazon](#)): \$12.83
0.2 non-skirted 96-well PCR plates ([USA Scientific/Amazon](#)): \$22.00/10 x3 = \$66.00
8-strip PCR caps ([USA Scientific](#)): \$11.00/125 strips x3 = \$33.00
Eppendorf DNA LoBind 1.5mL tubes ([USA Scientific](#)): \$31.75 per 250 tubes
15mL tubes ([Amazon](#)): \$17.99
10uL-100uL Pipette ([Amazon](#)): \$32.39
1000uL pipette ([Amazon](#)): \$32.39
10uL filtered pipette tips ([Amazon](#)): \$63.99
200uL filtered pipette tips ([Amazon](#)): \$57.77
1000uL pipette tips ([Amazon](#)): \$13.28
Alcohol swabs ([Amazon](#)): \$5.49
Kimwipes ([Amazon](#)): \$24.30
Fine-tip Sharpies ([Amazon](#)): \$1.84
Eliminase ([Ebay](#)): \$60.00

Total upfront consumable purchases: \$400 - \$500

Primers

ONT-tagged Forward Primers \$87.50
ONT-tagged Reverse Primers \$396.16

Total upfront primer cost: ~\$500

Reagents:

X-Amp DNA Reagent ([IBI Scientific](#)): 50mL kit is \$128.47. Cost per sample between \$0.039 (15uL) -

\$0.064 (25uL)

PCR Master Mix ([Empirical Bioscience](#)): \$206.90 shipped (\$10.34 per plate; \$0.108 per reaction/sample)

Ethanol: \$56.18 per 1L

[NEBNext Ultra II End Repair/dA-Tailing Module - 24 rxns](#) **New England**

Biolabs Catalog #E7546S

\$283.00 per 24 reactions

[Molecular Water](#) **IBI**

Scientific Catalog #IB42130

[HighPrep™ PCR Clean-up System](#) **MagBio Genomics**

Inc. Catalog #AC-60005

: \$117.88 per

50 mL. \$0.047 per rxn.

[Ligation Sequencing Kit \(Q20\)](#) **Oxford Nanopore**

Technologies Catalog #SQK-LSK112

:

\$694.43 per 6 reactions

[NEBNext Quick Ligation Module](#) **New England**

Biolabs Catalog #E6056S

: \$361.00 per 20

reactions

Total upfront reagent costs: ~\$1,600

Extraction ongoing costs:

Total per sample: \$0.093 - \$0.118

Total per 96 samples: \$8.93 - \$11.33

Total per Flongle run (480 samples) : \$44.65 - \$56.65

Total per Flongle run (672 samples) : \$62.50 - \$79.30

PCR ongoing costs:

Total per sample: ~\$0.225

Total per 96 samples: ~\$22.50

Total per Flongle run (480 samples) : ~\$110.00

Total per Flongle run (672 samples) : ~\$150.00

dA-tailing ongoing costs:

Total per Flongle run (1/2 rxns): \$5.95

Total per MinION run: \$11.85

Total per 96 samples: \$0.061

Total per sample (Flongle: 480 samples): \$0.012

Total per sample (Flongle: 672 samples): \$0.0089

Ligation ongoing costs:

Total per Flongle run (1/2 rxns): \$66.95

Total per MinION run: \$133.89

Total per 96 samples: \$13.38

Total per sample (Flongle: 480 samples): \$0.139

Total per sample (Flongle: 672 samples): \$0.0996

Sequencing ongoing costs:

Flongle cells: \$810 per 12 cells + \$50 per shipment (likely want 3-4 shipments) = ~\$100 per cell

Total upfront costs:

~\$5700 - \$6200 + thermocycler + compute

Total ongoing costs:

Per Flongle run: \$333 per 480 samples = \$0.69 per specimen

















Per Flongle run: \$402 per 672 samples = \$0.59 per specimen

ABSTRACT

This collection of protocols outline a working methodology to DNA barcode fungal specimens. This process will work for dried tissue, fresh tissue, or with DNA template that has already gone through an alternative extraction protocol. It takes you from computer/software setup through DNA extraction, amplification, library preparation, sequencing, and primary data analysis.

This particular workflow is designed for use with **Flongle 9.4.1 flowcells** and **V12 Ligation** chemistry.

FILES

-   **ONT Sequencing IT/Compute Pop!_OS 22.04 Setup
Version 1
by Stephen Douglas Russell**
-   **Preparing ONT-tagged Primers and Master Mix for Fungal DNA Barcoding
Version 2
by Stephen Douglas Russell**
-   **ONT DNA Extraction for Fungal Barcoding
Version 1
by Stephen Douglas Russell**
-   **ONT Post-PCR Pooling & Purification for Fungal Barcoding
Version 2
by Stephen Douglas Russell**
-   **ONT dA-tailing for Fungal Barcoding
Version 2
by Stephen Douglas Russell**
-   **ONT Q20+ (V12) Adapter Ligation for Fungal DNA Barcoding
Version 2
by Stephen Douglas Russell**
-   **ONT Flongle Flowcell Loading with Q20+ (V12) Chemistry
Version 2
by Stephen Douglas Russell**
-   **ONT Basecalling, Demultiplexing, and Analysis for Fungal Barcodes
Version 1
by Stephen Douglas Russell**