



DEC 08, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.e6nvwdq47lmk/v1](https://dx.doi.org/10.17504/protocols.io.e6nvwdq47lmk/v1)

**Protocol Citation:** Martin T. Henrich, Fanni F. Geibl 2023. DAB staining protocol for subsequent stereological cell counting. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.e6nvwdq47lmk/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

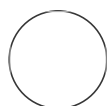
**Created:** Dec 04, 2023

## DAB staining protocol for subsequent stereological cell counting

Martin T. Henrich<sup>1,2</sup>, Fanni F. Geibl<sup>1,2</sup>

<sup>1</sup>Northwestern University, Feinberg School of Medicine;

<sup>2</sup>Philipps University Marburg, Clinic for Psychiatry and Psychotherapy



enrico.zampese

### ABSTRACT

This protocol describes the steps for performing a double chromogen staining using DAB and SK- 4700. Stained sections can subsequently be imaged and used for stereological cell quantification with Stereo Investigator software from MBF Bioscience or other Stereology software.

### ATTACHMENTS

[n75mb93sf.pdf](#)

### MATERIALS

#### Materials:

- Pre-cut brain sections (30-50 µm thick)
- 0.1 M PB
- Triton-X100 (detergent)
- Methanol (e.g. from Sigma-Aldrich)
- Ethanol absolute (e.g. from Sigma-Aldrich)
- 30% H<sub>2</sub>O<sub>2</sub> (e.g. from Merck)
- Xylene (e.g. from J.T.Baker)
- Na-hypochlorite (e.g. Sigma-Aldrich)
- Normal donkey serum (S30-100ML, Sigma-Aldrich)
- Primary antibodies (e.g. Anti-ChAT AB144P, Merck; anti-TH ab113, abcam; anti-p62, ab109012, abcam; anti-alpha-synuclein(phosphoS129), ab51253, abcam)
- Biotinylated secondary antibody
- Vectastain® Elite® ABC HRP Kit (Vector Laboratories)
- 3,3'-diaminobenzidine (DAB) (e.g. from Serva)
- SK-4700 SG Peroxidase Kit (Vector Laboratories)
- 12-well-plates with netwell inserts (e.g. Corning Costar Netwell)
- Quick-hardening Mounting Medium (e.g. Eukis, Sigma Aldrich)
- Microscopy slides

**Last Modified:** Dec 08, 2023

**PROTOCOL integer ID:** 91986

- Glass coverslips
- Designated container for storage of microscopy slides
- Orbital shaker (e.g. Heidolph Duomax 1030)
- Plajorm shaker (e.g. Heidolph Unimax 1010)
- Microscope system configured for brightfield (and multi-channel fluorescent work) with Stereo Investigator software from MicroBrightField (MBF) (e.g. Zeiss Axio Imager.M2 with MBF extension modules needed for Stereo Investigator software)

#### Recommended PPE:

- Lab coat/disposable gown
- Safety goggles
- Examination gloves
- FFP2 mask/respirator

#### Primary antibodies:



Goat anti-ChAT antibody Merck Millipore (EMD Millipore) Catalog #AB144P



Anti-Tyrosine Hydroxylase antibody Abcam Catalog #ab113



Recombinant Anti-SQSTM1 / p62 antibody [EPR4844] - Autophagosome Marker Abcam Catalog #ab109012



Recombinant Anti-Alpha-synuclein (phospho S129) antibody Abcam Catalog #ab51253

## DAB staining - Before the procedure:

12h 47m



- 1 Prepare [M] 0.1 Molarity (M) PB (PB).
- 2 Prepare PBT ([M] 0.1 Molarity (M) PB with 0.3% Triton X-100).
- 3 Prepare a quenching solution (for 40 mL solution: 32 mL [M] 0.1 Molarity (M) PB, 4 mL


methanol (100%),  4 mL  $\text{H}_2\text{O}_2$  (30%).



- 4 Prepare NDS solution (5% normal donkey serum diluted in PBT).

## DAB Staining Procedure



12h 47m

- 5 Place brain sections ( 30  $\mu\text{m}$   50  $\mu\text{m}$ ) in 12-well-plates with netwell inserts (up to 6-10 sections per netwell insert, depending on the size of the sections).



- 6 Place 12-well-plates on platform shaker and wash sections for 3x 5 min in  0.1 Molarity (M) PB. Exchange PB solution in between washing steps.

- 6.1 Place 12-well-plates on platform shaker and wash sections for  00:05:00 in  0.1 Molarity (M) PB. Exchange PB solution in between washing steps (1/3).



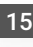
5m

- 6.2 Place 12-well-plates on platform shaker and wash sections for  00:05:00 in  0.1 Molarity (M) PB. Exchange PB solution in between washing steps (2/3).

5m

- 6.3 Place 12-well-plates on platform shaker and wash sections for  00:05:00 in  0.1 Molarity (M) PB. Exchange PB solution in between washing steps (3/3).

5m

- 7 Quench sections for  00:15:00 at  Room temperature in quenching solution in 12-well-plate  15m

### Note

We recommend  4 mL solution per well for good results.

8 Thereafter, wash sections for 4x 5 min in **0.1 Molarity (M)** PB on platform shaker. Exchange PB solution in between washing steps.



8.1 Thereafter, wash sections for **00:05:00** in **0.1 Molarity (M)** PB on platform shaker. Exchange PB solution in between washing steps (1/).



8.2 Thereafter, wash sections for **00:05:00** in **0.1 Molarity (M)** PB on platform shaker. Exchange PB solution in between washing steps (2/4).



8.3 Thereafter, wash sections for **00:05:00** in **0.1 Molarity (M)** PB on platform shaker. Exchange PB solution in between washing steps (3/4).



8.4 Thereafter, wash sections for **00:05:00** in **0.1 Molarity (M)** PB on platform shaker. Exchange PB solution in between washing steps (4/4).

9 Block sections for **01:00:00** at **Room temperature** in NDS solution in 12-well-plates.



1h

#### Note

We recommend **4 mL** solution per well for good results.

10 For incubation with primary antibody (e.g. anti neuronal nuclei (NeuN), Merck Millipore, MAB377, 1:100) transfer sections in a new 12-well-plate but without netwell inserts. This allows better shaking **Overnight**. Solution for incubation with primary antibody should contain NDS solution and the respective primary antibody diluted according to manufacturer recommendation. Incubate sections on







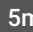

orbital shaker with gentle shaking at  4 °C  Overnight .


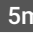

#### Note




We recommend at least 1 ml solution for each well.





**11**  On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections 4x 5 min in PB on a platform shaker at  Room temperature . Change washing solution after each washing step.

**11.1** On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  00:05:00  in PB on platform shaker at  Room temperature . Change washing solution after each washing step (1/4).

**11.2** On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  00:05:00  in PB on platform shaker at  Room temperature . Change washing solution after each washing step (2/4).




**11.3** On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  00:05:00  in PB on platform shaker at  Room temperature . Change washing solution after each washing step (3/4).




**11.4** On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  00:05:00  in PB on platform shaker at  Room temperature . Change washing solution after each washing step (4/4).

**12**  Next, transfer sections back in 12-well-plate without netwell inserts for incubation with biotinylated secondary antibody (e.g., biotinylated donkey anti-mouse, Jackson ImmunoResearch, 715-065-151, 1:500) diluted in NDS solution. Incubate for  01:00:00  on orbital shaker at  Room temperature .


#### Note

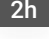



We recommend at least 1 ml solution for each well.

- 13**  Next, without washing, incubate sections in 12-well-plate without netwell inserts on orbital shaker at **Room temperature** for **01:00:00** in avidin-biotin-peroxidase solution using the Vectastain ABC Kit, diluted according to manufacturer's recommendation. **1h**
- 14** Prepare 5% DAB working solution fresh before color reaction (for 40 ml solution: **35.2 mL** **0.1 Molarity (M)** PB, **4 mL** DAB solution (**5 undetermined**), **0.8 mL** H<sub>2</sub>O<sub>2</sub> (1%)). Wear an FFP2/N95 mask or a respirator for this step and all further steps in which you work directly with DAB. Only carry out this work under a fume hood.
- 15**  Without washing, transfer brain sections back in 12-well-plates with netwell inserts and start DAB color reaction by incubating sections in DAB working solution for **00:06:00** at **Room temperature** under fume hood, depending on desired color intensity. Gently shake sections during incubation to avoid sections sticking together. **6m**
- 16**  Stop color reaction by washing sections 4x 5 min. in **0.1 Molarity (M)** PB on plajorm shaker. **5m**  
Exchange PB solution in between washing steps. Use a fresh 12-well-plate which was not previously in contact with DAB solution. While washing use diluted Na-hypochlorite solution to neutralize DAB solution on used 12-well-plate and inserts. Let plastics neutralize in diluted Nahypochlorite solution **Overnight**.
- 16.1** Stop color reaction by washing sections for **00:05:00** in **0.1 Molarity (M)** PB on plajorm shaker (1/4). **5m**
- 16.2** Stop color reaction by washing sections for **00:05:00** in **0.1 Molarity (M)** PB on plajorm shaker (2/4). **5m**
- 16.3** Stop color reaction by washing sections for **00:05:00** in **0.1 Molarity (M)** PB on plajorm shaker (3/4). **5m**
- 16.4** Stop color reaction by washing sections for **00:05:00** in **0.1 Molarity (M)** PB on plajorm shaker (4/4). **5m**

- 17 Block sections again for  01:00:00 at  Room temperature in NDS solution in 12-well-plates or  1h orbital shaker.

Note

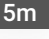


We recommend  4 mL solution per well for good results.

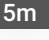


- 18 For incubation with primary antibody (e.g. anti-tyrosine hydroxylase (TH), Merck Millipore, AB152,  2h 1:1000), transfer sections in a new 12-well-plate but without netwell inserts. This allows better shaking  Overnight . solution for incubation with primary antibody should contain NDS solution and the respective primary antibody diluted according to manufacturer recommendation. Incubate sections on orbital shaker with gentle shaking at  4 °C  Overnight .

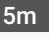


Note

We recommend at least  1 mL solution for each well.

- 19 On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections 4x 5 min in PB on platform shaker at  Room temperature . Change washing solution after each washing step.

- 19.1 On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  5m  00:05:00 in PB on platform shaker at  Room temperature . Change washing solution after each washing step (1/4).

- 19.2 On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  5m  00:05:00 in PB on platform shaker at  Room temperature . Change washing solution after each washing step (2/4).

- 19.3 On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for  5m  00:05:00 in PB on platform shaker at  Room temperature . Change washing solution after each washing step (3/4).

**19.4** On the next day, transfer sections into 12-well-plate with netwell inserts and wash sections for **00:05:00** in PB on platform shaker at **Room temperature**. Change washing solution after each washing step (4/4). **5m**

**20** Next, transfer sections back in 12-well-plate without netwell inserts for incubation with biotinylated secondary antibody (e.g., biotinylated donkey anti-rabbit, Jackson ImmunoResearch, 711-065-152, 1:500) diluted in NDS solution. Incubate for **01:00:00** on orbital shaker at **Room temperature**. **1h**

**Note**

We recommend at least **1 mL** solution for each well.

**21** Next, without washing, incubate sections in 12-well-plate without netwell inserts on orbital shaker at **Room temperature** for **01:00:00** in avidin-biotin-peroxidase solution using the Vectastain ABC Kit, diluted according to manufacturer's recommendation. **1h**

**22** Prepare Peroxidase Substrate Kit (SG Peroxidase Substrate Kit, Vector Labs, SK-4700) as recommended by the manufacturer. Wear an FFP2/N95 mask or a respirator for this step and all further steps in which you work directly with SK-4700. Only carry out this work under a fume hood.

**23** Without washing, transfer brain sections back in 12-well-plates with netwell inserts and start SK-4700 color reaction by incubating sections in prepared peroxidase substrate solution for 3- **00:06:00** at **Room temperature** under fume hood, depending on desired color intensity. Gently shake sections during incubation to avoid sections sticking together. **6m**

**24** Stop color reaction by washing sections 4x 5 min. in **0.1 Molarity (M)** PB on platform shaker. **5m**  
Exchange PB solution in between washing steps. Use a fresh 12-well-plate which was not previously in contact with SK-4700 solution. While washing use diluted Na-hypochlorite solution to neutralize SK-4700 solution on used 12-well-plate and inserts. Let plastics neutralize in diluted Na-hypochlorite solution **Overnight**.



- 24.1** Stop color reaction by washing sections for  00:05:00 in  0.1 Molarity (M) PB on platform  5m shaker (1/4).
- 24.2** Stop color reaction by washing sections for  00:05:00 in  0.1 Molarity (M) PB on platform  5m shaker (2/4).
- 24.3** Stop color reaction by washing sections for  00:05:00 in  0.1 Molarity (M) PB on platform  5m shaker (3/4).
- 24.4** Stop color reaction by washing sections for  00:05:00 in  0.1 Molarity (M) PB on platform  5m shaker (4/4).
- 25** After washing, mount sections on microscopy slides using a fine brush.
- 26** Let sections dry  Overnight  5m .
- 
- 27** Immerse dried slides in a series of ethanol (70%, 96%, 100%) 30 sec. each for stepwise dehydration. Work under a fume hood.
- 28** Let slides rest for  00:10:00 in 100% xylene solution. Work under a fume hood.  10m
- 29** Next, quickly apply a small amount of hard-drying mounting medium (e.g., Eukis) sufficient to cover the

sections. Carefully avoid the formation of air bubbles. Gently apply a coverslip over the sections and the mounting medium. Work under a fume hood.

### **DAB staining - After the procedure:**

- 30 Let slides cure for 2 days under a fume hood.
- 31 Dispose of waste and excess reagents/solution according to institutional guidelines.
- 32 Clean tools/working station.
- 33 Microscopy slides should be stored in a designated container at room temperature until time of observation.

### **Stereological cell counting - Before the procedure:**

- 34 Turn on microscope and computer according to specific manuals/instructions.
- 35 Turn on MBF Stereo Investigator software.

### **Stereological cell counting - Procedure:**

- 36 Load microscopy slides in the designated stage.

- 37 Use brightfield settings for imaging.
- 38 Carefully follow all steps of the Stereo Investigator Workflow according to the manual. The workflow is step-by-step and very intuitive.
- 39 We typically use a 60x oil immersion lens for imaging.
- 40 After cell counts have been collected export counting data into an excel sheet for further data analysis.

### **Stereological cell counting - After the procedure:**

- 41 Clean immersion objectives with a lens wipe and the appropriate cleaning solution.
- 42 Turn off software/microscope/laser according to specific instructions.