


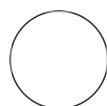

JUL 11, 2023

# A new metabarcoding approach to survey diversity at the species level of Arcellinida (Amoebozoa: Tubulinea)

R. González-Miguéns<sup>1</sup><sup>1</sup>Real\_Jardín\_Botánico\_CSIC

R. González-Miguéns

## ABSTRACT

PCR protocol for the paper: "A needle in a haystack: a new metabarcoding approach to survey diversity at the species level of Arcellinida (Amoebozoa: Tubulinea)"

OPEN  ACCESS

### DOI:

[dx.doi.org/10.17504/protocols.io.yxmvm2389g3p/v1](https://dx.doi.org/10.17504/protocols.io.yxmvm2389g3p/v1)

### Protocol Citation:

R. González-Miguéns 2023. A new metabarcoding approach to survey diversity at the species level of Arcellinida (Amoebozoa: Tubulinea).

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.yxmvm2389g3p/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

### Protocol status:

Working  
We use this protocol and it's working

**Created:** Nov 22, 2022

**Last Modified:** Jul 11, 2023

### PROTOCOL integer ID:

73112

## Polymerase Chain Reaction protocol

**Primers:**

Primer forward: LCO 1490 (5' GGTCAACAAATCATAAAGATATTGG 3')

Primer reverse: HCO 2198 (5' TAAACTTCAGGGTGACCAAAAAATCA 3')

**Reagents for the PCR:**

A	B
Polymerase	6 µL
BSA	1 µL
Primer forward	1 µL (10 µmol)
Primer reverse	1 µL (10 µmol)
Destilated water	1 µL
Sample (eDNA)	1 µL

Total reaction mix per sample= 10 µL


**Thermocycler program:**

A	B	C	D
	Temperature	Time	cycles
Initial denaturation	96 °C	5 min	1
Denaturation	94 °C	15 s	40
Annealing	40 °C	15 s	
Extension	72 °C	90 s	
Final extension	72 °C	10 min	1

Denaturation temperatures and times may vary depending on the polymerase used. This protocol was tested with the MyTaq™ Red Mix™ polymerase.

**Gel electrophoresis:**

Analyze the results of your PCR reaction via gel electrophoresis. The resultant fragment was 692 bp long:

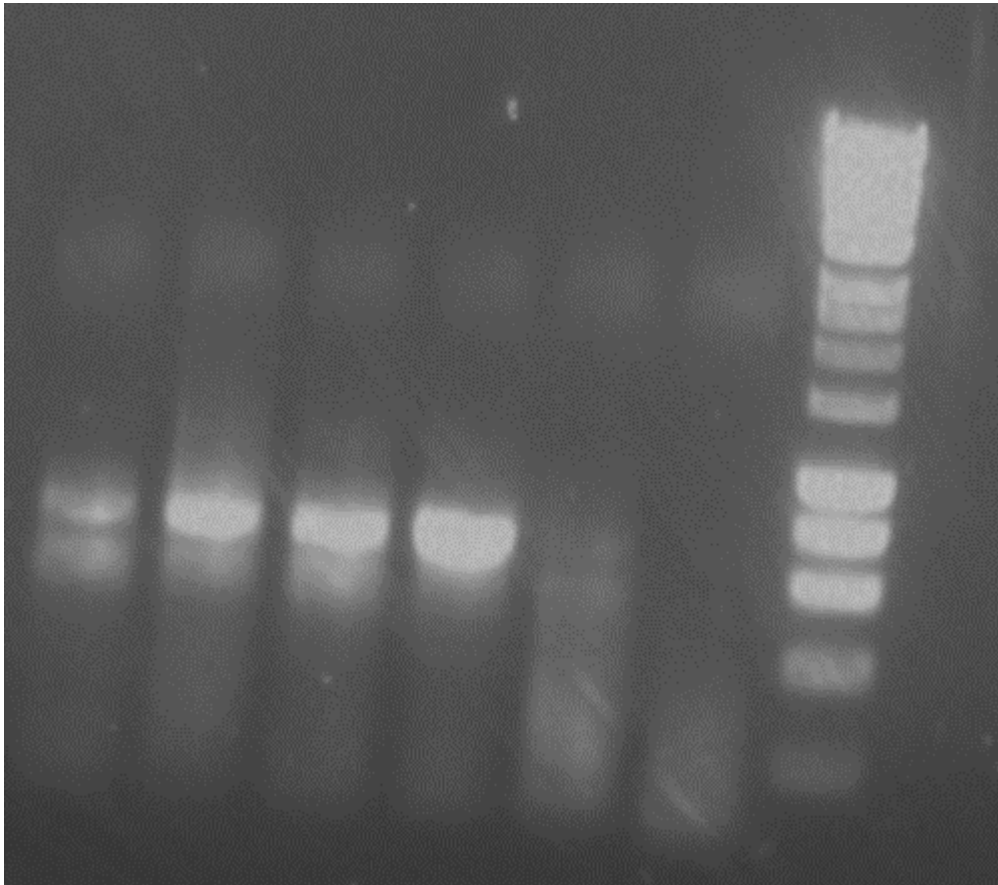

example of four PCR products with one negative control on the right. HyperLadder™ 1kb as the molecular weight marker.

## 2

### Second PCR

#### Primers:

Primer forward: LCO 1490 (5' GGTCAACAAATCATAAAGATATTGG 3')

Primer reverse: ArCOIR (5' CCACYNGAATGWGCTARAATACC 3')

#### Reagents for the PCR:

A	B
Polymerase	12 µL
Primer forward	1 µL (10 µmol)
Primer reverse	1 µL (10 µmol)
Destilated water	6 µL

A	B
Sample (first PCR product)	1 $\mu$ L

Total reaction mix per sample= 20  $\mu$ L

#### Thermocycler program:

A	B	C	D
	Temperature	Time	cycles
Initial denaturation	96 °C	5 min	1
Denaturation	94 °C	15 s	40
Annealing	55 °C	15 s	
Extension	72 °C	90 s	
Final extension	72 °C	10 min	1

Denaturation temperatures and times may vary depending on the polymerase used. This protocol was tested with the MyTaq™ Red Mix™ polymerase.

#### Gel electrophoresis:

Analyze the results of your PCR reaction via gel electrophoresis. The resultant fragment was 407 bp long