



JUL 18, 2023

## OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.ewov1q4bkgr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1q4bkgr2/v1)

**Protocol Citation:** Julian Liber 2023. 20 minute PCR Enzymatic Cleanup.  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.ewov1q4bkgr2/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

**Created:** Jul 17, 2023


**Last Modified:** Jul 18, 2023

**PROTOCOL integer ID:**  
 85121

## 20 minute PCR Enzymatic Cleanup

Julian Liber<sup>1</sup>

<sup>1</sup>Department of Biology - Duke University

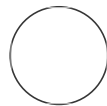

Julian Liber

### ABSTRACT

A rapid microplate-compatible protocol for cleanup of PCR products prior to Sanger sequencing. Verify that the PCR produces a single band, and retain some PCR product for cleanup and sequencing. Set up at room temperature, and react on the thermocycler.

**Mechanism:** PCR products contain 1) primers and 2) dNTPs which will interfere with downstream Sanger sequencing reactions. Exonuclease I degrades ssDNA (primers) and calf intestinal phosphatase degrades dNTPs. Functions like [NEB #E1050](#).

### MATERIALS

- Quick CIP (5U/uL) [NEB #M0525](#)
- Thermolabile Exonuclease I (20U/uL) [NEB #M0568](#)
- rCutSmart buffer [NEB #B6004](#) (or CutSmart buffer)
- NEBuffer r3.1 [NEB #B6003](#) (or NEBuffer 3.1)
- PCR-grade water

## Making the enzyme master mix

- 1 Add to a 1.5 mL sterile conical tube

- 760 uL water
- 94 uL 10x rCutSmart buffer
- 96 uL 10x Exonuclease I buffer (NEBuffer r3.1)
- 4 uL Quick CIP (5U/uL) [NEB #M0525](#)
- 10 uL Thermolabile Exonuclease I (20U/uL) [NEB #M0568](#)

2 Mix, then aliquot 200 uL into separate tubes. Store at -20C. Good for at least 2 years.

## PCR Cleanup

3 Mix 1.2 uL of master mix with 1.5 uL PCR product in PCR strip tubes or plates.

### Note

A 1:1 ratio of enzyme mix:PCR product can also be used. Pipetting 2 uL may be easier, and you can do equal parts in such a case. If doing a large number of reactions, I generally use 2 uL because it is faster to pipette that volume and ensure it stays in the tube.

h

4 Seal tubes or plates and place on a thermocycler. React at 37°C for 15 minutes, then 80°C for 1 minute.

5 Add the product to Sanger sequencing reactions. 1 uL of product, 1 uL of 10 uM primer, and 10 uL of water per tube. If the band is faint, you can use 2 uL or more of product.

### Note

Follow the sequencing facility's requirements for template DNA and primer concentration. The above amounts are recommendations based on what has worked for us.

