

JUL 16, 2023

Preparation of parts HygR-LHRZ and ZeoR-LHRZ

Dariusz Abramczyk¹

¹University of Edinburgh



Dariusz Abramczyk

OPEN ACCESS

dx.doi.org/10.17504/protocol s.io.4r3l271o3g1y/v1

Protocol Citation: Dariusz Abramczyk 2023. Preparation of parts HygR-LHRZ and ZeoR-LHRZ. protocols.io https://dx.doi.org/10.17504/p rotocols.io.4r3l271o3g1y/v1

License: This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: In development

We are still developing and optimizing this protocol. This is a protocol for preparation of a double-feature part (antibiotic selection marker-LHRZ)

Created: Mar 09, 2023

Last Modified: Jul 16, 2023

PROTOCOL integer ID:

78443


ABSTRACT

This cloning protocol is for preparation tandem pair HygR-LHRZ and ZeoR-LHRZ as a part of combinatory assembly system for preparation of integration/insertion arrays. https://www.protocols.io/edit/nanochromosome-arrays-combinatory-assembly-cprkvm4w

Either the Hyg^R-LHR^Z or the Zeo^R-LHR^Z tandem pairs are designed to be use for all arrays and are displayed on 3'-end each of arrays.


During inchworming/swapping HR double cross-over recombination process, the antibiotic selection marker is replaced in every HR event while the LHRZ works as a homology recombination region used during all inchworming/swapping rounds. (see https://www.protocols.io/edit/nanochromosome-arrays-combinatory-assembly-cprkvm4w)

Both parts are AsiSI free (this is for an optional excise site the cassette from the plasmid to isolate an insertion array before Pichia transformation)

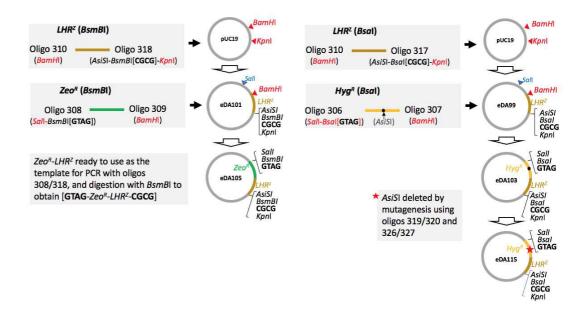

General scheme for preparation library of AntR-LHRZ

, inte..

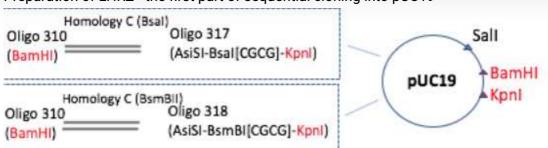

0.2 PCR preparation part LHRZ for cloning with Kpnl/BamHI (3'-end of LHRZ with Bsal or BsmBI RE site delivered by an oligo)

02:00:00

PrimeSTAR GXL Premix Contributed by users Catalog #R051B

A	expected size 998bp	expected size 998bp
components	LHRZ for (BsmBI)	LHRZ for (Bsal)
eDA8 (125pg/uL)	1	1
oligo 310	1.25	1.25
oligo317		1.25
oligo 318	1.25	
PrimeStar Premix 2x	25	
water	up to 50	


LHR(non-coding DNA) delivered from the DNA library design by dr Valentin Zulkover and synthesised by prof Adele Marston Adele Marston lab

A	В	С
310 (F)	LHRZ (BamHI	CCGGATCCC CGTGTAGGT CTACAAACT GG
318 (R)	LHRZ (Kpnl- BsmBl-AsiSI)	GCGGTACCC GTCTCAcgcg GCGATCGCC CTGCAGGTT GAGTTGGCG AAGGTGCG

3

A	В	С
317 (R)	LHRZ (Kpnl- Bsal-AsiSI)	GGGGTACCG GTCTCAcgcg GCGATCGCC CTGCAGGTT GAGTTGGCG AAGGTGCG

Program	time	number of cycles
98oC	10 sec	33 x
55oC	10sec	
68oC	1min 30 sec	
4oC	hold	

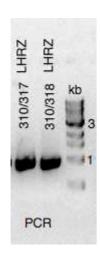

0.3

PCR clean-up QIAquick PCR Purification Kit (Qiagen)QIAquick PCR Purification Kit Elution with water 50uL

0.4

KpnI and BamHI digestion of PCR product and pUC19 (all enzymes NEB) BamHI/ KpnI

Kpnl-HF - 20,000 units **New England Biolabs Catalog** #R3142L

BamHI-HF - 50,000 units **New England Biolabs Catalog** #R3136L

A	В	С	D

30m

A	В	С	D
10 x CutSmart	8	8	8
BamHI-HF (20u/uL)	1.5	1.5	1
Kpnl-HF (20u/uL)	1.5	1.5	1
PCR LHRZ 310/317	50		
PCR LHRZ 310/318		50	
pUC19 298 ng/uL			10
water	up to 80 ->		

RE digestion reactions



0.5 🕙 00:30:00

PCR clean-up QIAquick PCR Purification Kit (Qiagen)QIAquick PCR Purification Kit

Elution with water 50uL

Checking DNA concentration

LHRZ 310/317 (BamHI/KpnI) - ~100ng/ul

LHRZ 310/318 (BamHI/KpnI) - ~100ng/ul

0.6 Ligation reactions T4 DNA ligase (NEB)

T4 DNA Ligase - 20,000 units **New England Biolabs Catalog**

A	1	2
pUC B/K ~100ng/uL) BamHI/KpnI	1	1
T4 DNA ligase	1	1
LHRZ (100ng/ul) B/K(310/317)	5.5	
10 X buffer	1	1
water	up to 10 ->	
LHRZ (100ng/ul) (310/318) B/K		5.5

Overnight 16oC

30m

0.7 👏 02:00:00

E.coli transformation (DH5alpha chemically LiAc competent cells)<u>E.coli chemical competent cells</u>




0.8

65 01:00:00

Several colonies obtained on both agar plates. Single colonies transferred on new LB +100ug/mL ampicillin plate, growth overnight and a single colonies submitted for bacterial colony PCR.


see protocol colony PCR

Material loaded on agarose gel for verification

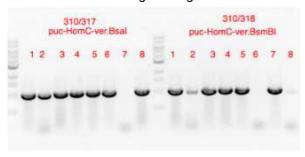

E.coli colony PCR


Positively verified clones re-cultured in LB + 100 ug/ml

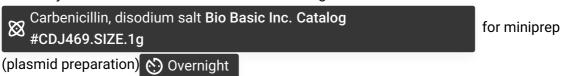

0.9 Plasmid isolation (miniprep) Qiagen miniprep 👏 01:00:00

PCR verified and Sanger sequencing intermediate plasmid LHRZ(Bsal)-pUC19 called eDA99 eDA99.dna

1 Preparation of ZeoR-LHRZ and HygR-LHRZ

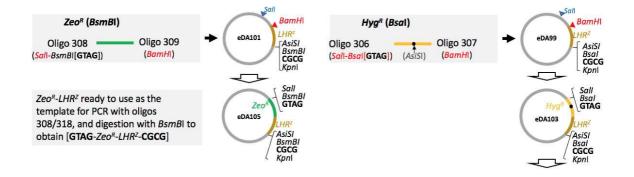

The final step - insertion of ZeoR or HygR into plasmid eDA101 and eDA99 correspondingly

1.1 PCR amplification of Zeocin expression cassette (Bsal-free) for cloning with BsmBl Zeocin and hygromycin cassettes originally derived from plasmids pGS-GnTI and pGS-GnTII (Glycoswitch) glycoswitch plasmids

PrimeSTAR GXL Premix Contributed by users Catalog #R051B

components	ZeoR	HygR
eDA48 (175 (ng/uL)		1
eDA78 (125pg/uL)	1	
oligo 308	1.25	
oligo309	1.25	
oligo 306		1.25
oligo 307		1.25
PrimeStar Premix 2x	25	25
water	up to 50	up to 50uL

A	В	С
308 (F)	Zeocin cassette (Sall-BsmBl) overhang GTAG	CTGTCGACC GTCTCAgtag CCCACACAC CATAGCTTC AAAATG
309(R)	Zeocin cassette (BamHI)	GGGATCCGC AAATTAAAG CCTTCGAGC G

A	В	С
306 (F)	HygR (Sall- Bsal)	CCTGTCGAC GGTCTCAgta gGACATGGA GGCCCAGAA TACCC
307 (R)	HygR (BamHI)	GCGGATCCC AGTATAGCG ACCAGCATT CACATAC

1.2 🕙 00:30:00

PCR clean-up QIAquick PCR Purification Kit (Qiagen)QIAquick PCR Purification Kit Elution with water 50uL

1.3 Sall and BamHI digestion of PCR product and pUC19 (all enzymes NEB)BamHI

BamHI-HF - 10,000 units **New England Biolabs Catalog** #R3136S

Sall-HF - 2,000 units **New England Biolabs Catalog** #R3138S

components	В	С	D	E
10x cutsmart buffer NEB	8	8	8	8
BamHI-HF (20u/uL)	1.5	1.5	1.5	1.5
eDA99 70 ng/uL			30	
eDA101 81 ng/uL				30
PCR hygR 306/307		50		
PCR zeo 308/309	50			
Sall-HF (20u/uL)	1.5	1.5	1.5	1.5
water	up to 80 ->			


1.4 🕙 00:30:00

PCR clean-up QIAquick PCR Purification Kit (Qiagen)QIAquick PCR Purification Kit Elution with water 50uL

30m

4h

30m

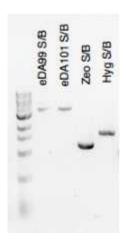

loaded 2 ul of each fraction

1.5 Ligation reactions T4 DNA ligase (NEB)

T4 DNA Ligase - 20,000 units **New England Biolabs Catalog** #M0202S

T4 DNA ligase	1	1
10 X buffer	1	1
eDA101 (S/B)		3
zeo 308/309 (S/B)		5
eDA99 (S/B)	3	
hygR 306/307 (S/B)	5	
water	up to 10 -	

👏 Overnight 16oC

1.6 👏 02:00:00

E.coli transformation (DH5alpha chemically LiAc competent cells) <u>E.coli chemical competent</u> cells

Selection on LB + 100ug/mL carbenicillin and growth at Overnight 37 °C

Carbenicillin, disodium salt Bio Basic Inc. Catalog #CDJ469.SIZE.1g


2h

1.7 🕙 01:00:00

.Several colonies obtained on both agar plates. Single colonies transferred on new LB +100ug/mL ampicillin plate, growth overnight and a single colonies submitted for bacterial colony PCR

see protocol colony PCR

(oligo 131/40 for verification ZeoR-LHRC-puC see plasmid map eda105.dna (oligo 306/42 for verification HygR-LHRC-puC see plasmid map eda103.dna Material loaded on agarose gel for verification

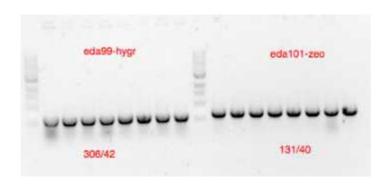

All colonies verified positively

Positively verified clones re-cultured in LB ampicillin for miniprep (plasmid preparation)

🖒 Overnight

1.8 Plasmid isolation (miniprep) Qiagen miniprep 01:00:00

PCR verified and Sanger sequencing intermediate plasmid ZeoR-LHRC-puC see plasmid

🛭 eda105.dna


Tandem paired ZeoR-LHRZ part (BsmBI-free, AsiSI-free,) ready for cloning with BsmBI as a PCR template for a use in the protocol

https://www.protocols.io/edit/nanochromosome-arrays-combinatory-assembly-cprkvm4w

PCR verified and Sanger sequencing intermediate plasmid HygR-LHRC-puC see plasmid map

eda103.dna
The part is Bsal-free but contains internal AsiSI in hph gene, requires
mutagenesis

2 Mutagenesis of AsiSI in eDA103 plasmid by Gibson assembly

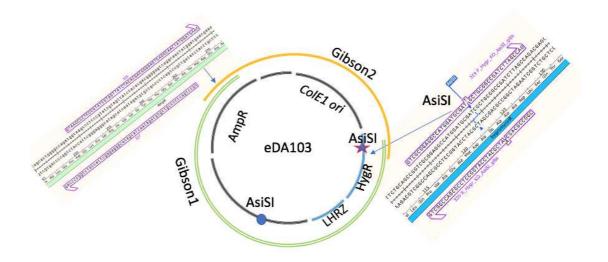

Schematics of mutagenesis of AsiSI site in HygR by Gibson assembly.

Note

eDA103 contains two AsiSI sites, but only AsiSI site in HygR needs to be removed as a part of integration/insertion array. Another AsiSi site (in pUC19 vector) is not taking under a consideration.

2.1 © 01:30:00 PCR amplification for further Gibson assembly reaction

PrimeSTAR GXL Premix Contributed by users Catalog #R051B

PrimeStarHS 2x premix (Takara) High Fidelity

components	for Gibson 1	for Gibson 2
water	up to 50 uL	for Gibson 1
Primestar2x	25	25
oligos 327/320		1.25/1.25
oligos 319/326	1.25/1.25	
eDA103 (74pg/uL)	1.5	1.5

PROGRAM TAK-2.30	time reaction	cycles
98C	10sec	
55C	10 sec	

PROGRAM TAK-2.30 time reaction		⊗y⊗ykes
68C	2min 30 sec	
4C	hold	

PCR reaction treatment with DpnI for 00:30:00 at 3 37 °C

2.2 **(:)** 00:30:00

PCR clean-up QIAquick PCR Purification Kit (Qiagen)QIAquick PCR **Purification Kit**

Elution with water 50uL.


DNA Molar concentration of purified PCR products was calculated using bioline conventor

DNA concentration convertor

Gibson 1 fragment 86 ng/uL

Gibson 2 fragment 46 ng/ul

PCR products verified on 1% agarose gel (4 uL loaded)

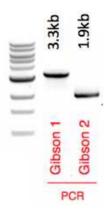

© 00:30:00

2.3 Gibson assembly of 3.3kb fragment 1 and 1.9kb fragment 2

37 °C 37 °C

30m

30m

Gibson Assembly Cloning Kit - 10 rxns New England Biolabs Catalog #E5510S

A	В	
Gibson1 319/326 (3.3kb) 40nM 3.3 kb	1	
Gibson2 327/320 (2.05 kb) 40nM	1	
Gibson 2x mix	2.5	
water	0.5	

Gibson assembly protocol protocol NEB Gibson Assembly Master Mix

Immediately after Gibson reaction is E.coli transformation (DH5alpha chemically LiAc competent cells) E.coli chemical competent cells

Selection on LB +100ug/mL ampicillin and growth at 🚫 Overnight



2.4 Single colonies verified by colony PCR




Several colonies obtained on both agar plates. Single colonies transferred on new LB +100ug/mL ampicillin plate, growth overnight and a single colonies submitted for bacterial colony PCR

see protocol colony PCR


Oligos 41 and 307 used for primary verification.

Oligo sequences available in the plasmid eDA115 map 0 eda115.dna



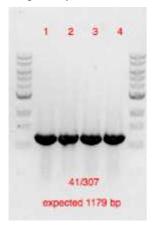

2.5 Positively verified clones re-cultured in LB + 100 ug/m

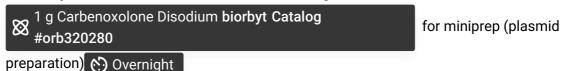



2.6 Plasmid isolation (miniprep) Qiagen miniprep (5) 01:00:00

2.7 Restriction digestion verification. Plasmid eDA115 eda115.dna digested with AsiSI (NEB)

AsiSI - 2,500 units **New England Biolabs Catalog** #R0630L

Α	В	С
AsiSI	1	1
eDA115 326 ng/uL	3	
eDA103 103 ng/uL		8
10 x CutSmart	3	3
water	up to 30 ->	

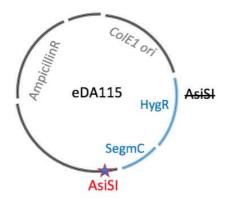

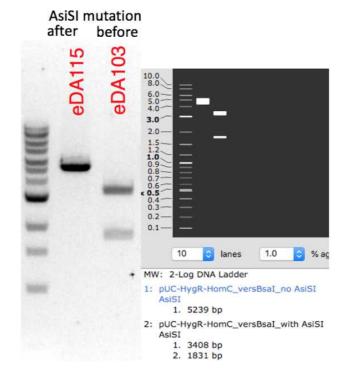

Tandem paired Hyg-LHRZ part (Bsal-free, AsiSI-free,) ready for cloning with Bsal as a PCR template for a use in the protocol

https://www.protocols.io/edit/nanochromosome-arrays-combinatory-assembly-cprkvm4w