



NOV 02, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.e6nvwo5zlmk/v1](https://dx.doi.org/10.17504/protocols.io.e6nvwo5zlmk/v1)

**Protocol Citation:** Vera A. Tang, Sean M Cook, Joanne Lannigan, Jennifer Jones, Joshua A Welsh 2023. Resource 7: rEV immunophenotyping. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.e6nvwo5zlmk/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working

**Created:** Jan 30, 2023

## Resource 7: rEV immunophenotyping

In 1 collection

Sean M

Vera A. Tang<sup>1</sup>, Cook<sup>2</sup>,  
 Joshua A Welsh<sup>4</sup>

Joanne Lannigan<sup>3</sup>, Jennifer Jones<sup>2</sup>,


<sup>1</sup>University of Ottawa;

<sup>2</sup>Laboratory of Pathology, Translational Nanobiology Section, Centre for Cancer Research, National Institute of Health, National Institutes of Health;

<sup>3</sup>Flow Cytometry Support Services;

<sup>4</sup>Translational Nanobiology Section, Laboratory of Pathology, Center for Cancer Research, National Cancer Institute, National Institutes of Health

Translational Nanobiology Section

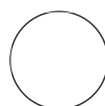

Jennifer Jones

### DISCLAIMER

This protocol summarizes key steps for a specific type of method, which is one of a collection of methods and assays used for EV analysis in the NCI Translational Nanobiology Section at the time of submission of this protocol. Appropriate use of this protocol requires careful, cohesive integration with other methods for EV production, isolation, and characterization.

**Last Modified:** Nov 02, 2023

**PROTOCOL integer ID:** 76114

**Funders**

**Acknowledgement:**

NIH

Grant ID: ZIA BC011502

NIH

Grant ID: ZIA BC011503

NIH

Grant ID: 4UH3TR002881-03

## ABSTRACT

Flow cytometry (FCM) is a common extracellular particles (EPs), including viruses and extracellular vesicles (EVs), characterization method. Frameworks such as MIFlowCyt-EV exist to provide reporting guidelines for metadata, controls, and data reporting. However, tools to optimize FCM for EP analysis in a systematic and quantitative way are lacking. Here, we demonstrate a cohesive set of methods and software tools that optimize FCM settings and facilitate cross-platform comparisons for EP studies. We introduce an automated small particle optimization (SPOT) pipeline to optimize FCM fluorescence and light scatter detector settings for EP analysis and leverage quantitative FCM (qFCM) as a tool to further enable FCM optimization of fluorophore panel selection, laser power, pulse statistics, and window extensions. Finally, we demonstrate the value of qFCM to facilitate standardized cross-platform comparisons, irrespective of instrument configuration, settings, and sensitivity in a cross-platform standardization study utilizing a commercially available EV reference material.

## MATERIALS

Exosome standards fluorescent Merck MilliporeSigma (Sigma-Aldrich) Catalog #SAE0193

, referred to as 'rEVs'

DPBS Thermo Fisher Scientific Catalog #14190144 , referred to as 'DPBS'

5 mL Round-bottom tube Corning Catalog #352052 , referred to as 'FACS tube'

Protein LoBind Tubes (5 mL) Thermo Fisher Scientific Catalog #0030122356

referred to as 'low-binding Eppendorf tubes'

Low Protein Binding Microcentrifuge Thermo Fisher Scientific Catalog #88379

96-well V bottom plate Contributed by users

PE anti-human CD81 (TAPA-1) Antibody BioLegend Catalog #349505 , referred to as 'anti-CD81-PE'

APC anti-human CD81 (TAPA-1) Antibody BioLegend Catalog #349509 , referred to as 'anti-CD81-APC'

Pacific Blue™ anti-human CD81 (TAPA-1) Antibody BioLegend Catalog #349515

, referred to as anti-CD81-PB'

QbSure Cytex Catalog #B7-10005 , referred to as 'QbSure'

100 nm polystyrene NIST bead Thermo Fisher Scientific Catalog #3100A



referred to as '100 nm NIST bead'

## Sample Preparation

35m



- 1 Briefly centrifuge the rEVs  100 x g, 4°C, 00:05:00 before opening

5m

- 2 Add  100 µL of  4 °C deionized water. Pipette up and down to mix.

### Note

Do not vortex

- 3 Dilute the reconstituted rEVs 1 in 5 in PBS for staining\* Take  50 µL of the rEV stock and add  200 µL of PBS.

#### Note

*The sample concentration chosen for staining should take into consideration the need to dilute the sample further after staining to reduce the concentration of free antibodies, since the samples will be acquired without a washing step to remove these antibodies as in the case of cells.*

- 4 Prepare the antibody dilutions for a 2x staining concentration of each of the antibodies to be used in the antibody titration. See the example below for sample calculations prepared for a staining concentration of 2 µg/ml.














A	B	C	D	E	F
	AB stock concentration (µg/ml)	2x staining conc (µg/ml)	Stock volume for 2x staining conc (µL)	Total volume for 2x staining conc (µL)	Volume of PBS to be added (uL)
anti-Human CD81-PE	120	4	2	60	58
anti-Human CD81-APC	200	4	2	100	98
anti-Human CD81-PB	300	4	2	150	148


Sample calculations prepared for a staining concentration of 2 µg/ml. All Antibodies clones are 5A6.

- 5 Samples will now be prepared using the following plate map.



A	B	C	D	E	F	G	H
	<b>CD81 PE</b>	<b>CD81 APC</b>	<b>CD81 Pacific Blue</b>	<b>CD81 PE</b>	<b>CD81 APC</b>	<b>CD81 Pacific Blue</b>	
	rEV	rEV	rEV	PBS	PBS	PBS	
<b>0 µg/mL</b>	A1	A2	A3	A4	A5	A6	
<b>0.0625 µg/mL</b>	B1	B2	B3	B4	B5	B6	

	A	B	C	D	E	F	G	H
	0.125 µg/mL	C1	C2	C3	C4	C5	C6	
	0.25 µg/mL	D1	D2	D3	D4	D5	D6	
	0.5 µg/mL	E1	E2	E3	E4	E5	E6	
	1 µg/mL	F1	F2	F3	F4	F5	F6	
	2 µg/mL	G1	G2	G3	G4	G5	G6	QbSure Beads
		PBS	PBS	PBS	PBS	PBS	PBS	100 nm PS NIST-Traceable Beads

- 6 In a 96-well V-bottom plate, add  10 µL 1:5 rEV solution to wells A1-A3
  
- 7 Add  10 µL DPBS to wells A1-A3 and reverse pipet to mix. These wells will serve as the rEV controls
  
- 8 Add  20 µL DPBS to wells A4-A6. These wells will serve as buffer only controls
  
- 9 In the same 96-well V-bottom plate from  go to step #6, add  10 µL 1:5 rEV solution to wells B1-G3. Add  10 µL DPBS to wells A1-G1. Add  10 µL from the PE working solution tubes to wells B1-G1, add  10 µL from the APC working solution tubes to wells B2-G2, add  10 µL from the PB working solution tubes to wells B3-G3.
  
- 10 In the same 96-well V-bottom plate from  go to step #6, add  10 µL from the corresponding PE working solutions tube to wells B4-G4. Add  10 µL from the corresponding APC working solutions tube to wells B5-G5. Add  10 µL from the corresponding PB working solutions tube to wells B6-G6. These wells will serve as the DBPS+AB control wells.


11 Cover and incubate this plate for  00:30:00 at RT.

30m


12 Using a separate 96-well V-bottom plate, add  199  $\mu\text{L}$  DPBS into wells A1-G6. Pipet  1  $\mu\text{L}$  from the incubation plate into the new plate in the same wells and reverse pipette to mix.

#### Note

Steps 11 and 12 will be repeated for the ON incubations as well to bring rEVs down to running concentration that doesn't swarm the instrument detectors.

13 Add  200  $\mu\text{L}$  DPBS into wells A7-G7. These wells will help reduce sample carryover after each rEV+AB combination is acquired.

14 Re-cover incubation plate and let incubate ON at RT to repeat measurements next day as directed in the following steps.

15 Label as FACS tube as 'QbSure', and add  500  $\mu\text{L}$  DPBS. Vortex the QbSure beads for 5 sec and add 3 drops QbSure beads into FACS tube.

## Cytometer Setup


16 Setup cytometer to acquire at FCM<sub>PASS</sub> output gains for both light scatter and fluorescent detectors as determined by the FCM<sub>PASS</sub> detector optimization module outputs. For the Aurora, an optimal gain template is returned to import into the Spectraflo software.



17 On the Cytex Aurora, set window extension to 0. On the CytoFLEX platform turn on 'High Acquisition Mode'.

## Note

To set the window extension to 0 on the Cytex Aurora, navigate to the 'Lasers' tab under the sample acquisition settings and set the window extension to 0.

To change the acquisition mode on the CytoFLEX S, click on the 'Advanced' menu on the top-of-the-screen ribbon. Next click on 'Event Rate Settings' menu, change the acquisition mode to 'High', and press 'Ok'.

- 18 Set the cytometer triggering threshold to the violet SSC parameter and run a DPBS control first to ensure the cytometer is clean and thresholds/event rates have remained unchanged. The background event rate should be ~1000 events/sec. **All samples should be acquired with the lowest flow rate, typically ~10-15  $\mu\text{L min}^{-1}$ .**
- 19 Acquire all wells for at least 60 sec.
- 20 In an open well, run the 100 nm polystyrene NIST beads at the same settings as the rEVs (diluted in DPBS to a concentration of 5E6 p/mL) until at least 10000 bead events have been acquired.
- 21 In an open well, add  200  $\mu\text{L}$  from the QbSure FACS tube and collect 10,000 bead events at the same settings as the rEVs.

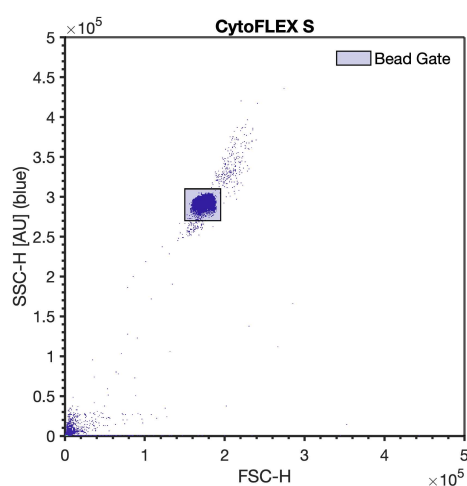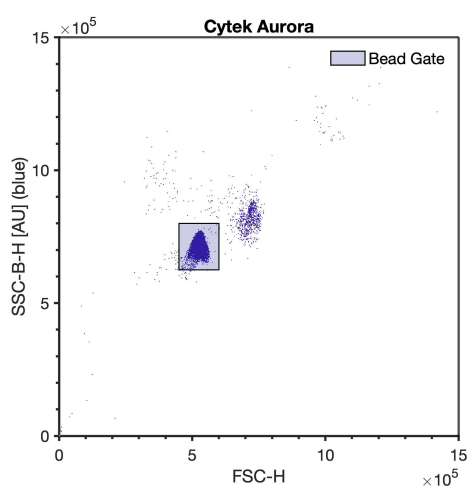

QbSure beads gated on blue SSC-H and FSC-H on the Cytek Aurora and CytoFLEX S. The violet SSC settings remain unchanged. In order to limit .fcs file size, a FSC scatter trigger may be

used until 10000 QbSure bead events are acquired.

#### Note

You may need to adjust the B-SSC and FSC detector settings for the bead acquisition. Do not adjust fluorescence detector settings or high sensitivity SSC detector settings. Assuming the QbSure beads were cross calibrated in the previous protocols, this allows for the calibration of the rEV fluorescent data into calibrated units.

- 22 If calibration of data into standard units is desired, the 100 nm polystyrene NIST bead and QbSure beads can be used to calibrate data by following the FCM<sub>PASS</sub> experiment calibration protocols.

#### Protocol



NAME

**FCMPASS Protocol Collection**

CREATED BY

**Joshua A Welsh**

**PREVIEW**

- 23 Don't forget to repeat the acquisition for the samples incubating ON.