



Jul 28, 2021

Analysis of glycosphingolipids from human plasma

David A Priestman^{1,2}, Danielle te Vrugte¹, Kerri-Lee Wallom¹, María E Fernández-Suárez¹, Maria Leondaraki^{1,3}, Carissa Drake¹, Frances M Platt^{1,2,4,5}

¹Department of Pharmacology, University of Oxford, Mansfield Road, Oxford, OX1 3QT, United Kingdom;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815;

³Faculty of Medicine, Goethe University Frankfurt am Main, Theodor-Stern-Kai, 60590 Frankfurt am Main, Germany;

⁴Royal Society Wolfson Research Merit Award Holder; ⁵Wellcome Trust Investigator in Science.

2 Works for me

Share

dx.doi.org/10.17504/protocols.io.busvwnw6

Oke Avwenagha

ABSTRACT

Interest in the role of cellular glycosphingolipids (GSLs) in health and disease led to us developing a sensitive method to analyse the full complement of GSL structures present in mammalian cells, fluids and tissues. The original qualitative method we developed was published in 2004 and measured the oligosaccharides selectively released from glycosphingolipids using a ceramide glycanase enzyme derived from the medicinal leech. We have now updated and refined this protocol with the focus on achieving sensitive and reproducible quantitation of GSLs in control and patient plasma samples. The method uses the fluorescent compound anthranilic acid (2-AA) to label oligosaccharides prior to analysis using normal-phase high-performance liquid chromatography. The labelling procedure is rapid, selective, and easy to perform. With the inclusion of a 2AA-labelled chitotriose calibration standard, it is possible to obtain accurate and reproducible molar quantities of individual GSL species.

DOI

dx.doi.org/10.17504/protocols.io.busvwnw6

PROTOCOL CITATION

David A Priestman, Danielle te Vrugte, Kerri-Lee Wallom, María E Fernández-Suárez, Maria Leondaraki, Carissa Drake, Frances M Platt 2021. Analysis of glycosphingolipids from human plasma. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.busvwnw6>

FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's through the Michael J. Fox Foundation for Parkinson's Research (MJFF).
Grant ID: ASAP-0478

MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Analysis of fluorescently labeled glycosphingolipid-derived oligosaccharides following ceramide glycanase digestion and anthranilic acid labeling DCA Neville, V Coquard, DA Priestman, DJM te Vrugte, DJ Sillence, Raymond A Dwek, Frances M Platt, Terry D Butters Analytical biochemistry 331 (2), 275-282, (2004). doi: 10.1016/j.ab.2004.03.051.

KEYWORDS

Glycosphingolipids, Glucosyl Ceramide, HPLC, Oligosaccharide Analysis

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 06, 2021

LAST MODIFIED

Jul 28, 2021

OWNERSHIP HISTORY

May 06, 2021 Urmilas

May 20, 2021 Oke Avwenagha

PROTOCOL INTEGER ID

49717

GUIDELINES

This protocol requires the use of some hazardous materials. As such, users must be appropriately trained and hazardous materials stored, used, and disposed of in accordance with your institution's health and safety policies, and approved laboratory policies, risk assessments and codes of practice.

MATERIALS TEXT

Acetonitrile 1.00030

Gradient grade for liquid chromatography LiChrosolv® Reag. Ph Eur

https://www.sigmaaldrich.com/catalog/product/mm/100030?lang=en®ion=GB&clid=EA1aIQobChMI1JyZu_Xa8AIvCwUGAB10IqVLEAAAYASAAEgJHmfD_BwE

Anthranilic acid

A89855 Sigma-Aldrich

Reagent grade, ≥98%

<https://www.sigmaaldrich.com/catalog/product/sial/a89855?lang=en®ion=GB>

Boric acid

ReagentPlus®, ≥99.5%

B0252 Sigma-Aldrich

https://www.sigmaaldrich.com/GB/en/product/sigald/b0252?cm_sp=Insite-_-caSrpResults_srpRecs_srpModel_boric%20acid%20reagent%20plus-_-srpRecs3-1

Chloroform

Suitable for HPLC, ≥99.8%, amylene stabilized

34854-M Sigma-Aldrich

<https://www.sigmaaldrich.com/catalog/product/sial/34854m?lang=en®ion=GB>

Discovery® DPA-6S SPE Tube

https://www.sigmaaldrich.com/GB/en/search/dpa-6s?focus=products&page=1&perPage=30&sort=relevance&term=DPA-6S&type=product_name

Gibco™ DPBS, no calcium, no magnesium

14190094

<https://www.thermofisher.com/order/catalog/product/14190094#/14190094>

Glycosphingolipid standards

<https://www.matreya.com/Departments/Product-Categories/Glycosphingolipids.aspx>

Kinesis SPE Columns: TELOS® C18(EC) 100mg/1ml SPE Columns

<https://kinesis.co.uk/kinesis-spe-columns-telosr-c18-ec-100mg-1ml-spe-columns-210-100m-001t.html>

Ludger - BioQuant 2AA Labelled Chitotriose Standard

Cat. #: BQ-CAA-CHI-01 Batch: B37I-02

<https://www.ludger.com/docs/products/bq/bq-cao-chi/bq-cao-chi-01-b37i-02-cofa.pdf>

Ludger - 2-AA Labeled Glucose Homopolymer Ladder

Cat. # CAA-GHP-30

<https://www.glycoprofiling.com/docs/products/cao/ludger-cao-ghp-30-guide.pdf>

Methanol

34860 Sigma-Aldrich

Suitable for HPLC, ≥99.9%

https://www.sigmaaldrich.com/catalog/product/sigald/34860?lang=en®ion=GB&cm_sp=Insite-_-caSrpResults_srpRecs_srpModel_methanol%2034860-_-srpRecs3-1

rEGCase I: recombinant Endoglycoceramidase custom-synthesized by GenScript

<https://www.genscript.com/>

Sarstedt 1.5 ml Micro tubes with screw cap and seal

Product reference number 72.692

<https://www.sarstedt.com/en/products/laboratory/screw-cap-micro-tubes-reaction-tubes/screw-cap-micro-tubes/product/dmsarproducts/specificationPdf/Product/72.692.005/>

Sodium acetate trihydrate

S1304 Sigma-Aldrich

Meets USP testing specifications

https://www.sigmaaldrich.com/catalog/product/sial/s1304?lang=en®ion=GB&clid=CjwKCAjwJ2FBhAuEiwAlKu19jp2DL5CW1Ke85Y0Ksb5boyV-5zKyUkzvg3uVVuZXV3abjGZbVruJhoCescQAvD_BwE

Sodium cyanoborohydride

Reagent grade, 95%

156159 Sigma-Aldrich

<https://www.sigmaaldrich.com/catalog/product/aldrich/156159?lang=en®ion=GB>

TSKgel® Amide-80 HPLC Column

<https://www.sigmaaldrich.com/GB/en/product/supelco/813071>

SAFETY WARNINGS

This protocol requires the use of some hazardous solvents, reagents and chemicals. Refer to the Safety Data Sheets (SDS) provided by supplier and applicable Control of Substances Harmful to Health (COSHH). The correct personal protective equipment must be worn, and incidents reported in line with your institution's policy and procedures.

BEFORE STARTING

Check that you have the required reagents, solvents, chemicals, equipment and PPE.

GSL preparation from plasma

1 Use  50 µl plasma for GSL quantification.

2 

Add **150 µl** de-ionised water to make the volume up to **200 µl** in a 1.5 ml screw-cap tube.

3 

Add **0.8 mL** of chloroform/methanol (1:2, v/v) to give (C/M/W 4:8:3 final).



Chloroform and methanol are very toxic, and methanol is flammable. Refer to the Safety Data Sheets.

4 

Leave **Overnight** at **4 °C**.

5 Vortex.

6 

10m

Centrifuge at 16,000 x g for **00:10:00** at **Room temperature**.

7 

Transfer supernatant (about **1 mL**) to new tube and separate into two phases: by adding **0.2 mL** PBS and then **0.2 mL** chloroform.

8 Vortex.




9 

10m

Centrifuge at 16,000 x g for **00:10:00** at room temperature.

10 Remove very carefully the lower phase to a new tube and retain the upper phase.

11 Dry down the lower phase under a stream of (oxygen-free) nitrogen in heating block (**42 °C**).

- 12 When dry, re-suspend the lower phase in  **20 µl** chloroform/methanol (1:3).
- 13 Add upper phase to lower phase and vortex.
- 14 Pre-equilibrate C18 columns (telos, Kinesis, UK) with 4 x  **1 mL** methanol and 3 x  **1 mL** deionised water.
- 15 Load lower/upper phase mix onto column, let drip through gravity flow.



16 

Rinse sample tube with  **1 mL** water, apply to column to wash.

17 

Wash column with 3 x  **1 mL** water.

18 Elute GSLs into a new tube with:

18.1  **2 mL** chloroform/methanol (98:2). Push through first  **0.5 mL** . You can use syringe with adapter that fits into top of column.

18.2  **2 mL** chloroform/methanol (1:3).


18.3  **1 mL** methanol.

19 

Vortex and leave  **Overnight** at  **4 °C** or carry on to enzymatic digestion.

GSL digestion with EGC'ase I

16h

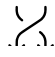20 



Vortex (5 mL , C18) and dry down samples under a stream of nitrogen in heating block (42 °C).

21 When about 150 µl sample remaining, transfer to 1.5 mL screw-cap tube.

22

Rinse sample tube with 200 µl C:M 2:1, vortex and combine with the rest of the sample in the screw-cap tube.

23 Dry down, under a slow stream of nitrogen in heating block (42 °C).

24 Re-suspend in 50 µl C:M 2:1, vortex, dry down under a very slow stream of nitrogen.

NB It is essential that ALL the sample is dried in the bottom of the tube.

25

Add 90 µl enzyme/buffer to each sample and vortex:

25.1 rEGCase: 4 µl enzyme (stored in freezer) plus 86 µl buffer per sample.

25.2 rEGCase buffer (stored at 4 °C): 0.6 % Triton in 50 Milimolar (mM) sodium acetate pH 5.2.

26

16h

Incubate at 37 °C for 16:00:00 .

rEGCaseI: recombinant Endoglycoceramidase I was custom synthesised by GenScript.

2AA labelling of glycans released from GSLs

1h 30m

27

Labelling mix: 30 mg/mL 2AA in labelling buffer and 45 mg/mL sodium cyanoborohydride
Labelling buffer: 4% sodium acetate and 2% boric acid in methanol

Dissolve 2AA in labelling buffer first with vortexing.

28  

Add 2AA in buffer to sodium cyanoborohydride and vortex.



Sodium cyanoborohydride is very toxic. Refer to the Safety Data Sheet.

29  


Add **310 µl** labelling mix to the **90 µl** sample (digest) in **1.5 mL** screw-cap tube.

30 

1h 30m

Incubate in oven at **80 °C** for **01:00:00**, vortexing at **00:30:00**.

31 Allow to cool to **Room temperature**.

32 

Add **1 mL** acetonitrile:water (97:3).

33 Transfer from screw-cap tube to **15 mL** tubes.

34 Rinse screw-cap tube with 2 x **1 mL** acetonitrile:water (97:3) and add to **15 mL** tube.

35 Pre-equilibrate 50 mg Discovery® DPA-6S SPE Tube (supplied by Sigma-Aldrich) with:




1 mL acetonitrile

2 x **1 mL** water

2 x **1 mL** acetonitrile

36 Apply samples ( 3.4 mL) to equilibrated DPA-6S columns, let drip through gravity flow.


37  

Wash columns with acetonitrile:water (95:5) → add  1 mL into the  15 mL tubes to wash and add to columns, and then add 3 x  1 mL acetonitrile:water (95:5) directly onto the columns.

38 

Elute with  600 µl water into new tubes.

39 

Add  60 µl :  140 µl (sample:acetonitrile) to HPLC vials and vortex.

40 For HPLC, inject  5 µl -  50 µl .

41 DPA-6S eluate is stable stored at  -20 °C .

HPLC protocol

1h

42 Purified 2AA-labeled oligosaccharides are separated and quantified by normal-phase high-performance liquid chromatography (NP-HPLC) as previously described (Neville *et al.*, 2004, see the reference below).

[Analysis of fluorescently labeled glycosphingolipid-derived oligosaccharides following ceramide glycanase digestion and anthranilic acid labeling](#)

DCA Neville, V Coquard, DA Priestman, DJM te Vrugte, DJ Sillence, Raymond A Dwek, Frances M Platt, Terry D Butters

Analytical biochemistry 331 (2), 275-282, (2004).

doi: 10.1016/j.ab.2004.03.051.

43 The HPLC system consists of a Waters Alliance 2695 separations module and an in-line Waters 2475 multi λ- fluorescence detector set at Ex λ360 nm and Em λ425 nm.

44 The solid phase used is a 4.6 × 250 mm TSK gel-Amide 80 column maintained at  30 °C (Anachem, Luton, UK).

45 The chromatographic flow rate is 0.8 mL/min, and run time was 60 min. The total run time is 01:00:00 .

1h

Time (min)	0	6	35	37	39	41	42	54	60
Flow (mL/min)	0.8	0.8	0.8	0.8	0.8	0.8	1.2	1.2	0.8
% Solvent A	71.6	71.6	52.8	23.0	23.0	71.6	71.6	71.6	71.6
% Solvent B	8.4	8.4	27.2	57.0	57.0	8.4	8.4	8.4	8.4
% Solvent C	20.0	20.0	20.0	20.0	20.0	20.0	20.0	20.0	20.0

Gradient conditions for Normal Phase HPLC. All chromatography was controlled and data were collected and processed using Waters Empower software.

46 The mobile phases are acetonitrile (solvent A), de-ionised water (solvent B) and 100 Milimolar (mM) ammonium acetate, 3.85 (solvent C).

47 Starting conditions are 71.6% solvent A. 8.4% solvent B and 20% solvent C at 0.8 mL/min constant flow-rate for 6 mins. 00:06:00

48 The gradient is developed from 6 min to 35 min by increasing solvent B from 8.4% to 27.2% with concomitant decrease in solvent A from 71.6% to 52.8 %. 00:29:00

49 From 35 to 37 min solvent B is increased to 57% and solvent A reduced to 23%, then maintained for 2 min up until 39 min. 00:04:00

50 Between 39 and 41 min the solvent ratios are returned to the starting conditions and then maintained for 1 min. 00:03:00

51 At 42 mins the flow-rate is increased to 1.2 mL/min and maintained for 12 min up to 54 min. 00:12:00

52 Between 54 and 60 min the flow-rate is gradually returned to 0.8 mL/min. 00:06:00






53

Individual GSL species are identified by their glucose unit values (GUs), calculated with the HPLC Empower software using a homopolymer dextran ladder (Fig. 2).

53.1 To help identify GSL peaks in the plasma samples, compare the GUs of the peaks in the plasma samples with those in a GSL standard mix.

53.2 

Prepare a mixture of commercially available authentic GSLs in a screw-cap tube:

 50 µl	Matreya Neutral GSL mix	Cat No 1505 (Glc-Cer, Lac-Cer, Gb3, Gb4)	1.0 mg/ml
 100 µl	Matreya Ganglioside mix	Cat No 1510 (Lac-Cer, GM3, GD3)	0.5 mg/ml
 100 µl	Matreya Ganglioside mix	Cat No 1511 (GA1, GM1a, GD1a, GD1b, GT1b)	0.5 mg/ml
 30 µl	Matreya GM2 Ganglioside	Cat No 1502 (GM2)	0.5 mg/ml
 30 µl	Sigma Asialo-GM2	Cat No G9398 (GA2)	1.0 mg/ml

53.3 Dry down the mixture under nitrogen and then digest, label and clean up as done for the plasma samples in steps 25-41.

53.4 Aliquot and store at  **-20 °C**.

53.5 In order to calculate molar quantities from integrated peaks in the chromatogram, inject a calibration standard containing 2.5 pmol 2AA-labelled chitotriose (Ludger) with each sample set (not shown).

Figures 1- 4: HPLC profile and sugar structures for GSLs in human plasma together with the GSL biosynthetic pathway

54



Figure 1 GSL Biosynthetic pathway in the background together with a typical HPLC trace below. GSLs detected in the profile are highlighted in the scheme.

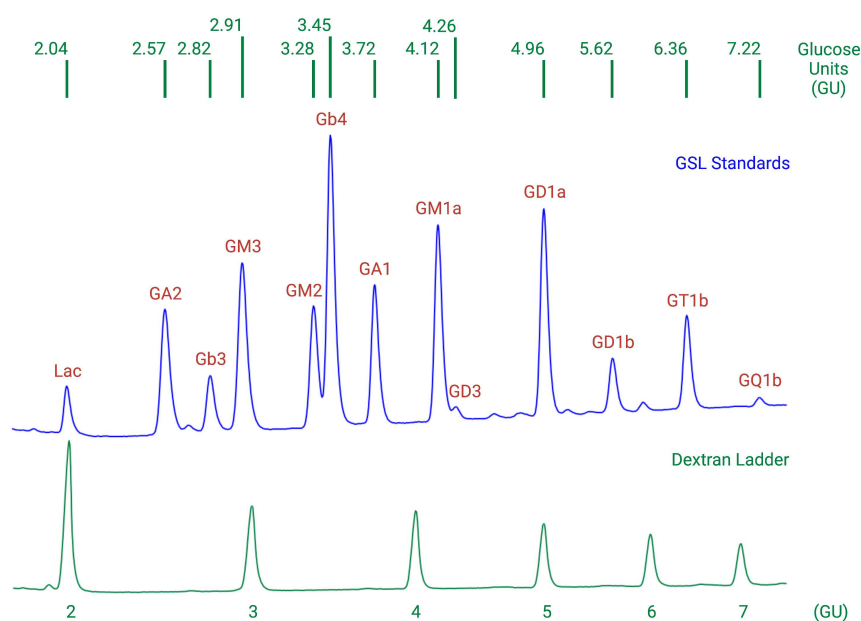


Figure 2 HPLC profiles for GSL standards and 2AA-labelled homopolymer dextran ladder. Glucose Units were calculated using the retention times of a 2AA-labelled homopolymer dextran ladder and are shown for each of the GSL standard peaks.

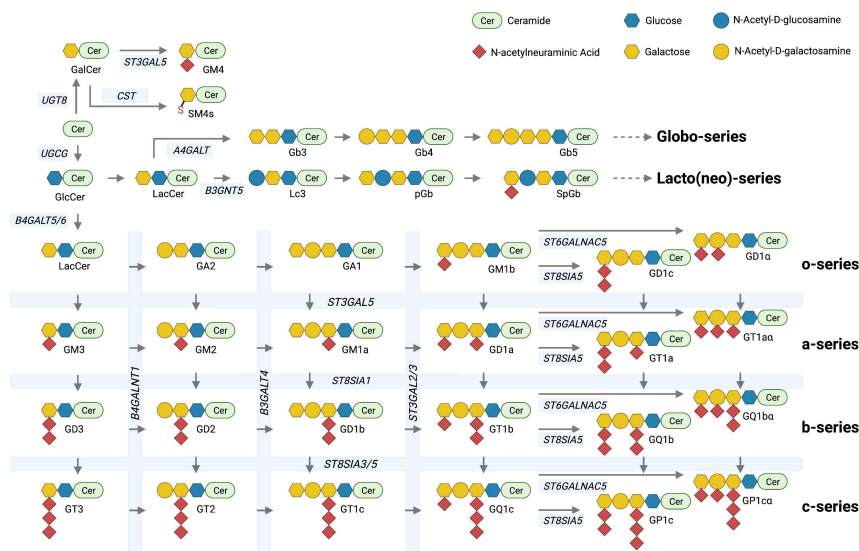


Figure 3 GSL biosynthesis. Biosynthetic enzymes genes are indicated in the blue grid. Ganglioside names are abbreviated according to Svennerholm [1]. [1] L. Svennerholm, Designation and schematic structure of gangliosides and allied glycosphingolipids, Prog Brain Res 101 (1994) XI-XIV.

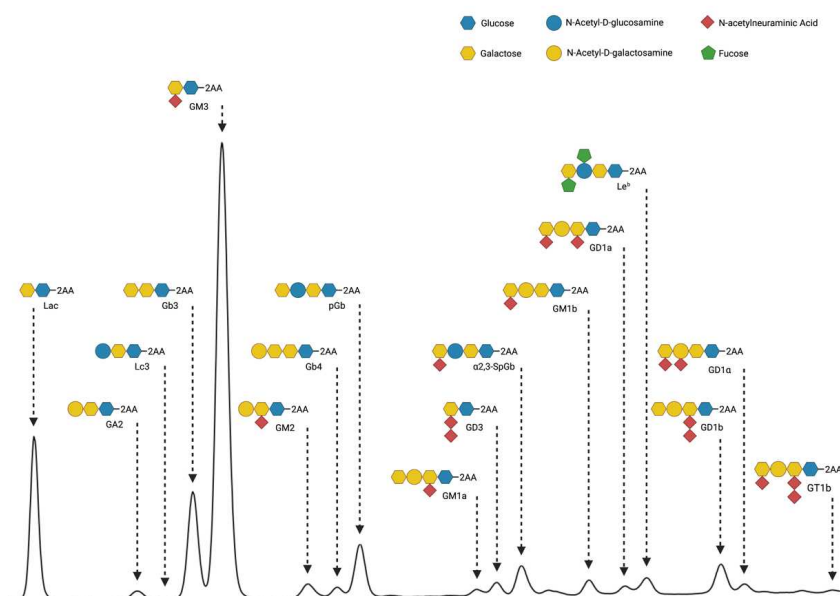


Figure 4 Typical HPLC profile for 2AA-labelled glycans released from human plasma glycosphingolipids with sugar structures indicated for each peak. In order to calculate values for total GSL concentrations, use the sum of the areas for all peaks from LacCer to GT1b. GlcCer, which elutes before LacCer, is not included because it is not fully hydrolysed by EGC'ase I.

Figures created with [BioRender.com](https://www.biorender.com)

