



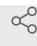
Version 2 ▼

Oct 25, 2022

Western Blotting (Fly Heads) V.2

Mel Feany^{1,2}¹Brigham and Women's Hospital; ²Harvard Medical School

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.8epv5j96jl1b/v2

Daniel's workspace



Daniel El Kodsí

ABSTRACT

This protocol describes how to perform a Western Blotting technique using fly heads.

DOI

dx.doi.org/10.17504/protocols.io.8epv5j96jl1b/v2

PROTOCOL CITATION

Mel Feany 2022. Western Blotting (Fly Heads). **protocols.io**<https://dx.doi.org/10.17504/protocols.io.8epv5j96jl1b/v2>

Version created by Daniel El Kodsí



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Oct 25, 2022

LAST MODIFIED

Oct 25, 2022

PROTOCOL INTEGER ID

71786

- 1 Homogenize desired number of fly heads in 1 X Laemmli sample buffer.

10m

2 Heat samples to 🔥 **100 °C** for ⌚ **00:10:00** , spin briefly before loading.

3 Load premade gel into western blotting apparatus. Fill reservoir with Running Buffer:
Running Buffer:

📏 **6 g** Tris-HCL

📏 **28.9 g** glycine

Fill to 📏 **1 L** with distilled water

Add 📏 **5 mL** 20% SDS

4 Load samples on gel and attach electrodes.

5 Run gel at 120 V until dye front reaches the bottom of the gel, ⌚ **01:00:00** . Run longer for ^{1h} greater separation.

6 Remove gel and transfer using Trans-Blot Turbo.

7 Perform antigen retrieval by microwaving ⌚ **00:09:00** in PBS.

9m

8 Block membrane in 1X PBS with 0.05% Tween-20 and 3% dry milk for ⌚ **01:00:00** .

1h


9 Add primary antibody at correct dilution in PBSTween + milk and incubate with shaking
⌚ **Overnight** at 🔥 **4 °C** .

10 Wash blot 3x in PBSTween, ⌚ **00:05:00** each, with shaking.

5m

Add secondary antibody at the correct dilution in PBSTween + milk, incubate with shaking at ^{3h}

11  **Room temperature** for  **03:00:00** .

12 Wash blot in PBSTween  **00:30:00** with frequent wash changes.

30m

13 Develop with ECL substrate or image fluorescence, as appropriate.