

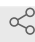


Sep 19, 2022

Creating Gelvatol Solution

Daniela Boassa^{1,2}, Mark Ellisman^{1,2}¹National Center for Microscopy and Imaging Research, University of California, San Diego, La Jolla CA;²Department of Neurosciences, University of California, San Diego, La Jolla CA

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.4r3l2o2rqv1y/v1NCMIR@UCSDNCMIR
University of California, San Diego

ABSTRACT

A gel-like mounting media (with an anti-fade ingredient) for fluorescence slides.

DOI

dx.doi.org/10.17504/protocols.io.4r3l2o2rqv1y/v1

PROTOCOL CITATION

Daniela Boassa, Mark Ellisman 2022. Creating Gelvatol Solution. **protocols.io**
<https://protocols.io/view/creating-gelvatol-solution-b52xq8fn>

FUNDERS ACKNOWLEDGEMENT

NIH/NIGMS
Grant ID: R24 GM137200-1

KEYWORDS

Gelvatol Solution, NCMIR, UCSD, Fluorescence Slides

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Mar 07, 2022

LAST MODIFIED

Sep 19, 2022

PROTOCOL INTEGER ID

59191

MATERIALS TEXT

- DABCO (in small fridge below oven; Sigma #D-2522)
- Polyvinyl Alcohol (in chemical bin#11; Sigma #D-8136)
- Glycerol (in chemical bin #7; Sigma #G-9012)
- 0.2M Tris buffer @ pH 8.5 (12ml = 2.4ml 1M Tris + 9.6 ml double distilled water)
- 1.5 ml microcentrifuge tubes
- 90 ml centrifuge tubes (including tubes for balance)

SAFETY WARNINGS

Be careful when handling DABCO, as it is corrosive to skin (must be treated as hazardous dry or sharps waste).

Label all tubes with "Gelvatol" and date prepared (in case a batch turns bad or is incorrectly prepared).

BEFORE STARTING

- Add a small drop of Gelvatol to specimen.
- Carefully place coverslip on the drop(s) avoiding bubbles.
- The mounting media will set overnight at 4 degrees or at room temperature 3-6 hours.

- 1 Add 2.4 g of [Polyvinyl Alcohol](#) (Mol. Wt. 30,000-70,000) to 6ml of [glycerol](#).
- 2 Stir well to mix.
- 3 Add 6ml of double distilled water and leave for at least 2 hours at room temperature (rotate overnight)
- 4 Add 12ml of 0.2M Tris (pH 8.5).

2.4 ml Tris 1M
9.6 ml H₂O

- 5 Heat to 50° C for 10 minutes with occasional mixing.
- 6 After polyvinyl alcohol is dissolved, clarify by centrifugation (5000g = ~6800 rpm) for 15 minutes at 20° C.
- 7 Use the large tabletop 5810R Eppendorf centrifuge.
- 8 Collect supernatant liquid and measure amount collected.
- 9 Add DABCO to 2.5% as antifade medium:

2.5g DABCO / 100ml of liquid
- 10 pH to 8.6 with pH paper (the paper will disintegrate, so pipette onto paper rather than putting paper into container).
- 11 Aliquot in 0.5ml microcentrifuge tubes and store at -20° C.
- 12 Stocks of Gelvatol are supposed to be stable at room temperature for several weeks after thawing.