

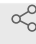



Aug 31, 2022

# Isolation of PBMCs From Whole Blood

Gerardo Ongari<sup>1</sup><sup>1</sup>IRCCS Mondino Foundation, Pavia

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.j8nlkwbm6l5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkwbm6l5r/v1) Gerardo Ongari

## ABSTRACT

This protocol details methods for the isolation of peripheral blood mononuclear cells (PBMCs) from Whole Blood

## DOI

[dx.doi.org/10.17504/protocols.io.j8nlkwbm6l5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlkwbm6l5r/v1)

## PROTOCOL CITATION

Gerardo Ongari 2022. Isolation of PBMCs From Whole Blood. **protocols.io**  
<https://protocols.io/view/isolation-of-pbmcs-from-whole-blood-cfzvtp66>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Aug 31, 2022

## LAST MODIFIED

Aug 31, 2022

## PROTOCOL INTEGER ID

69397

## MATERIALS TEXT

### Materials

- 4 EDTA Lavender top tubes (9mL/each)
- Density gradient medium (Histopaque-1077 Sigma-Aldrich)
- Phosphate Buffered Saline (DPBS without calcium without magnesium)
- Tubes for centrifugations
- 1.8 ml Cryotube vials (Nunc 055004)






#### SAFETY WARNINGS

Please refer to the Safety Data Sheets (SDS) for health and environmental hazards.

#### BEFORE STARTING

**Pre-procedure:** Ensure all reagents are at room temperature (15 - 25°C). The procedure must be carried out under a cell culture hood.








#### Step 1: Plasma Isolation

- 1 Collect  **36 mL** whole blood into EDTA tubes
- 2 Centrifuge the blood samples at  **1000 x g, Room temperature, 00:20:00** 20m
- 3 Carefully remove collection tubes from centrifuge; plasma is the top layer above the rest of the whole blood layer
- 4 Using caution not bring up the remaining whole blood layer, remove plasma and aliquot it into  **1.5 mL** microcentrifuge tubes
- 5 Centrifuge the isolated plasma at  **1600 x g, Room temperature, 00:20:00** , to remove <sup>20m</sup> residual cells and debris
- 6 After centrifugation, aliquot the plasma supernatant into 1mL aliquots by using microcentrifuge tubes labeled with patient ID
- 7 Freeze at  **-80 °C** .

#### Step 2: PBMC isolation

1h

- 8 Dilute the rest of the whole blood sample to a 1:1 volume ratio with the PBS

- 9 Add  **15 mL** density gradient medium to a fresh 50mL conical tube and gently layer the diluted blood on top of the density gradient medium. Take care not to mix the two layers.
- 10 Centrifuge at  **800 x g, Room temperature, 00:20:00 , with brake OFF** 20m
- 11 Carefully harvest the cells by inserting the pipette directly through the upper plasma layer to the mononuclear cells at the interface. Alternatively, you can first remove the upper layer and then collect the cells
- 12 Wash the harvested PBMC in  **12 mL** of PBS
- 13 Centrifuge at  **300 x g, Room temperature, 00:15:00** 15m
- 14 Discard the supernatant and wash the PBMC pellet in 5mL of PBS
- 15 Centrifuge at  **300 x g, Room temperature, 00:15:00** 15m
- 16 Discard the supernatant and resuspend the cell pellet in  **5 mL** of PBS
- 17 Count the cells using Trypan blue staining.
- 18 Resuspend PBMCs to  $5 \times 10^6$  cells/mL in PBS and aliquot them in cryotube vials (1mL/vial) labeled with patient ID
- 19 Centrifuge at  **800 x g, Room temperature, 00:10:00** 10m

20 Discard the supernatant and store the cell pellets at  $-80^{\circ}\text{C}$