



Jun 09, 2021

Cell growth and harvest for E. coli-based Cell-free Protein Synthesis

In 1 collection

Weston Kightlinger¹¹SwiftScale Biologics, Inc

1

Works for me



Share

dx.doi.org/10.17504/protocols.io.bvh4n38w

SwiftScale Biologics

 Emily Pennell

ABSTRACT

Cell-free protein synthesis reactions require cell-free extract derived from rapidly dividing E. coli cells to produce high yields of protein therapeutics, antigens, enzymes, vaccines, and other proteins. This protocol describes the cell growth and harvest steps for E. coli to produce extract.

Read more on cell-free protein synthesis here: <https://www.swiftscalebio.com/blog/cell-free-protein-synthesis>

ATTACHMENTS

[dn6kbkfm.pdf](#)

DOI

dx.doi.org/10.17504/protocols.io.bvh4n38w

PROTOCOL CITATION

Weston Kightlinger 2021. Cell growth and harvest for E. coli-based Cell-free Protein Synthesis. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bvh4n38w>

COLLECTIONS ⓘ



Collection of protocols for cell-free protein synthesis

KEYWORDS

Cell growth, E. coli, Cell-free Protein Synthesis, Harvest

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jun 03, 2021



LAST MODIFIED

Jun 09, 2021

OWNERSHIP HISTORY







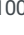






Jun 03, 2021  UrmilasJun 09, 2021  Emily Pennell

Materials:

1.  **2.5 L** Tunair flask with baffles
2. Autoclavable  **250 mL** glass bottles
3. Media components (see below)

A	B
2xYTPG Total Volume (mL) =	1000
2xYTP*** Total Volume (mL) =	750
NaCl (g)	5
Tryptone (g)	16
Yeast Extract (g)	10
Potassium Phosphate, dibasic (K ₂ HPO ₄) (g)	7
Potassium Phosphate, monobasic (KH ₂ PO ₄) (g)	3
Glucose Total Volume (mL) =	250
Glucose (g)	18

***Adjust 2XYTP pH to 7.2 using 5 N KOH.

4. MilliQ water
5.  **2.5 L** Tunair flask with baffles
6. Chosen E. coli strain (often BL21 (DE3) in standard protocols)
7.  **1 L** 2xYTPG ( **750 mL** of 2xYTP media in Tunair flask and autoclaved,  **250 mL** of  **72 g/L** glucose in a  **250 mL** glass bottle and autoclaved)
8. LB Media
9. S30 Buffer Components:
 - a. 100X:  **1 Molarity (M)** Tris-Acetate  **pH8.2** from Trizma Base, adjusted with Glacial Acetic Acid.
 - b. 100X:  **1.4 Molarity (M)** Mg(Acetate)₂ from Mg(OAc)₂ Tetrahydrate
 - c. 100X:  **6.0 Molarity (M)** K(Acetate)  **588.90 g** K(OAc)
 - d. Optionally, DTT can be added at  **2 Milimolar (mM)** using  **1 Molarity (M)** DTT as a 500X stock.

Cell growth and harvest for E. coli-based Cell-free Protein Synthesis

5m

1 

Start a  **50 mL** seed culture of E. coli strain and grow in LB media  **Overnight** at  **37 °C** .

2 

Add **250 mL** glucose to **750 mL** 2xYTP media in Tunair flask to obtain 2xYTPG media.

- 3 Inoculate Tunair flask culture at approximately $OD_{600} = 0.05$.



Grow culture at **37 °C** with shaking at **250 rpm** to $OD_{600} = 0.6$ and induce T7 with IPTG (if necessary). Add to **1 Milimolar (mM)** IPTG.



5m

Harvest culture at $OD_{600} = 3$ by centrifugation at **8000 x g** for **00:05:00** and washing with S30 buffer three times.

- 6 After the last wash, pour off S30 supernatant to obtain the pellet. Resuspend this pellet in S30 and lyse or immediately flash freeze on liquid nitrogen and store at **-80 °C** until ready for cell lysis.