


Nov 03, 2022

# DNA extraction for long-read sequencing of bacteria

Eby Sim<sup>1,2</sup><sup>1</sup>The University of Sydney Infectious Diseases Institute;<sup>2</sup>Centre for Infectious Diseases & Microbiology Public Health

1 Works for me

Share

[dx.doi.org/10.17504/protocols.io.4r3l2ox24v1y/v1](https://dx.doi.org/10.17504/protocols.io.4r3l2ox24v1y/v1) Eby Sim

## ABSTRACT

Utilisation of long-read sequencing can reliably generate complete bacterial genomes. Here, we present a DNA extraction method which introduces minor modifications to the DNeasy® UltraClean® Microbial Kit (Qiagen) to generate DNA suitable for long-read sequencing. When sequenced, the extracted DNA should yield median read lengths greater than 7 kb. In addition, this protocol uses cultures growing on solid media as a starting point which will be useful for laboratories that do not routinely use liquid cultures.

## DOI

[dx.doi.org/10.17504/protocols.io.4r3l2ox24v1y/v1](https://dx.doi.org/10.17504/protocols.io.4r3l2ox24v1y/v1)

## PROTOCOL CITATION

Eby Sim 2022. DNA extraction for long-read sequencing of bacteria .  
**protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.4r3l2ox24v1y/v1>



## KEYWORDS

DNA extraction, Long read sequencing

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Aug 16, 2022

LAST MODIFIED

Nov 03, 2022

PROTOCOL INTEGER ID

68691

## GUIDELINES


Users are reminded to be very gentle and deliberate in pipetting as this could shear DNA. Usage of wide-bore pipette tips (if accessible) will be preferable.

## MATERIALS TEXT

Please refer to the DNeasy UltraClean Microbial Kit Handbook for a list of all required equipment. No additional reagents are required.

## BEFORE STARTING

Users should take note that this protocol uses both mechanical lysis and spin-columns for DNA extraction. This protocol should **not** be used if reads lengths of > 60 kb is desired. The following image shows the typical read length distribution obtained from this methodology.

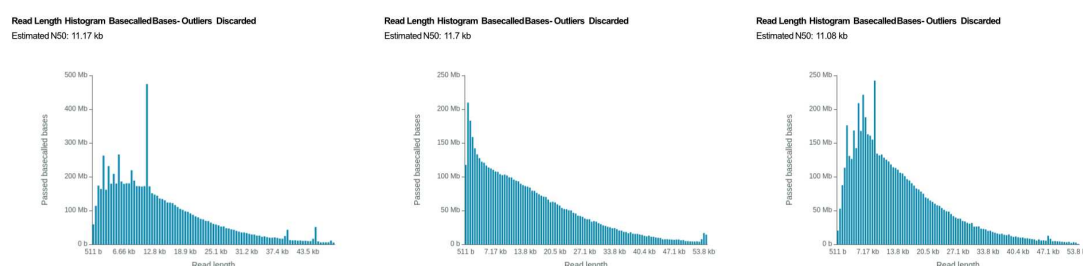























Read length distributions of three different ONT runs on DNA extracted with this methodology. Libraries were prepared using the Rapid Barcoding Kit (SQK-RBK004) following manufacturer's instructions.

## General bacteria culture

- 1 Streak bacteria of interest onto their respective, optimum solid media and incubate plates at optimal growth conditions.
- 2 On the day of DNA extraction, observe the plate to ensure purity. Do not attempt extraction if different colony morphologies are observed.

## DNA extraction

- 3 To a clean 2 mL Powerbead tube, add  **300 µL Powerbead Solution** and  **50 µL Solution SL** .
- 4 Take a 1 µL Inoculation loop and pick up 4 streaks from the first quadrant and dislodge the biomass into the Powerbead tube containing both the PowerBead solution and Solution SL.
- 5 Using a fresh 1 µL Inoculation loop, pick up 4 streaks from the second quadrant and dislodge biomass into the Powerbead tube containing both the PowerBead solution and Solution SL.
- 6 Briefly vortex the Powerbead tube to mix reagents and bacterial biomass.
- 7 Carefully affix the Vortex adaptor onto the Vortex-Genie 2 vortex and place Powerbead tubes<sup>2m</sup> horizontally, with the cap facing inwards onto the Vortex adaptor. Vortex at maximum speed for  **00:02:00** .
- 8 Remove Powerbead tubes from the Vortex adaptor and centrifuge at 1m  
 **10000 x g, Room temperature, 00:01:00** .
- 9 Carefully aspirate  **300 µL** of supernatant without disturbing the pellet or picking up beads. Slowly dispense entire volume into a 2 mL collection tube.
- 10 Add  **100 µL Solution IRS** and gently finger-flick the collection tube to mix. If there is liquid stuck on the underside of the lid after mixing, give the tube a quick wrist-flick to collect the liquid.
- 11 Incubate the tube  **On ice** for  **00:06:00** . 6m
- 12 Centrifuge the tube at  **10000 x g, Room temperature, 00:02:00** . 2m

- 13 Carefully aspirate  **300 µL** of supernatant, without disturbing the pellet, and slowly dispense into a 2 mL collection tube.
- 14 Add  **600 µL Solution SB** and gently invert tube to mix. If there is liquid stuck on the underside of the lid after mixing, pulse centrifuge collect the liquid.
- 15 Slowly aspirate  **700 µL** of the DNA mixture and dispense into a spin column and centrifuge at  **10000 x g, Room temperature, 00:00:30** . Discard the flow through. 30s
- 16 Slowly aspirate the remaining  **200 µL** and dispense into the same spin column (as step 15) and centrifuge at  **10000 x g, Room temperature, 00:00:30** . Discard the flow through and replace the collection tube with a new collection tube. 30s
- 17 Add  **300 µL Solution CB** and centrifuge at 30s  
 **10000 x g, Room temperature, 00:00:30** . Discard the flow through and replace the collection tube with a new collection tube.
- 18 Perform a dry centrifuge at  **10000 x g, Room temperature, 00:01:00** to remove residual ethanol. Replace collection tube with a new 1.5 mL Lo-Bind tube. 1m
- 19 Add  **50 µL Solution EB** to the centre of the membrane. Incubate at 3m  
 **Room temperature** for  **00:03:00** .
- 20 Centrifuge the spin column at  **10000 x g, Room temperature, 00:00:30** . Discard spin column 30s

#### DNA quality control

- 21 Assess quality on a Spectrophotometers (NanoDrop™ (ThermoFisher) or equivalent). Key quality metrics are listed below. Samples that do not fall within this value should not be sent for long-read sequencing. While nucleic acid concentration is also measured by the

spectrophotometer, it is not a key parameter at this stage.

A260/A280: 1.8 - 2.0  
A260/A230: 2.0 - 2.22

- 22 Assess the concentration of dsDNA on a fluorometer (Qubit or equivalent). Ensure that concentration and amount of extracted dsDNA meets the requirement of the sequencing technology.

Due to the higher DNA input requirements for long-read sequencing, the broad range assay would be the most appropriate kit to use on the Qubit.

- 23 Assess the integrity of the extracted dsDNA via electrophoretic separation (either via a 0.6% (w/v) agarose gel electrophoresis or TapeStation) using an appropriately sized ladder.

Ensure that majority of the dsDNA fragments are greater than 20 kb.

- 24 Store DNA at -80°C.

Do not freeze DNA before finishing steps 21-23. Keep DNA in 4°C until all quality checks are done.